

(19) World Intellectual Property  
Organization  
International Bureau



(43) International Publication Date  
22 July 2004 (22.07.2004)

PCT

(10) International Publication Number  
**WO 2004/061410 A2**

(51) International Patent Classification<sup>7</sup>: **G01N**  
(21) International Application Number:  
PCT/US2003/037090  
(22) International Filing Date:  
16 December 2003 (16.12.2003)  
(25) Filing Language: English  
(26) Publication Language: English  
(30) Priority Data:  
60/434,075 18 December 2002 (18.12.2002) US

(71) Applicants (for all designated States except US): **CIPHERGEN BIOSYSTEMS, INC.** [US/US]; 6611 Dumbarton Circle, Fremont, CA 94555 (US). **QUEEN ELIZABETH HOSPITAL** [CN/CN]; 30 Gascoigne Road, Kowloon, Hong Kong SAR (CN).

(72) Inventors; and

(75) Inventors/Applicants (for US only): **YIP, Timothy, Tak, Chun** [CN/CN]; 1, Kapok Path, Westwood, Palm Springs, Yuen Long N.T., Hong Kong SAR (CN). **CHO, Chi, Shing**

[—/CN]; Flat 9E, Block 4, The Tolo Place, Sunshine City, Ma On Shan, N.T., Hong Kong SAR (CN). **AU, Siu, Kie** [—/CN]; c/o Department of Clinical Oncology, Queen Elizabeth Hospital, 30 Gascoigne Road, Kowloon, Hong Kong SAR (CN). **YIP, Tai-Tung** [US/US]; 1532 Aster Court, Cupertino, CA 95014 (US). **YIP, Christine, L.** [US/US]; 1532 Aster Court, Cupertino, CA 95014 (US). **YIP, Victor, F.** [US/US]; 33008 Compton CT, Union City, CA 94587 (US).

(74) Agents: **BENT, Stephen, A.** et al.; Foley & Lardner, Washington Harbour, 3000 K Street, N.W. Suite 500, Washington, DC 20007-5101 (US).

(81) Designated States (national): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BW, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, EG, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NI, NO, NZ, OM, PG, PH, PL, PT, RO, RU, SC, SD, SE, SG, SK, SL, SY, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, YU, ZA, ZM, ZW.

(84) Designated States (regional): ARIPO patent (BW, GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM),

[Continued on next page]

(54) Title: **SERUM BIOMARKERS IN LUNG CANCER**

MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION
IM-1	2011	A	IM-37	3893	A	IM-72	54026	A	IM-109	2882	B
IM-2	2030	A	IM-38	3960	A	IM-73	60170	A	IM-110	2967	B
IM-3	2069	A	IM-39	3972	A	IM-75	74372	A	IM-111	2977	B
IM-4	2128	A	IM-40	3984	A	IM-76	75545	A	IM-112	2994	B
IM-5	2148	A	IM-41	4068	A	IM-77	77543	A	IM-113	3031	B
IM-6	2188	A	IM-42	4178	A	IM-78	79507	A	IM-114	3048	B
IM-7	2232	A	IM-43	4287	A	IM-79	89854	A	IM-115	3148	B
IM-8	2277	A	IM-44	4297	A	IM-80	101831	A	IM-116	3166	B
IM-9	2295	A	IM-45	4309	A	IM-81	104301	A	IM-117	3283	B
IM-10	2318	A	IM-46	4484	A	IM-82	125160	A	IM-118	3308	B
IM-11	2411	A	IM-47	4649	A	IM-83	132976	A	IM-119	3332	B
IM-12	2434	A	IM-48	4798	A	IM-84	149099	A	IM-120	3432	B
IM-13	2467	A	IM-49	5104	A	IM-85	2018	B	IM-121	3450	B
IM-14	2482	A	IM-50	5918	A	IM-86	2029	B	IM-122	3561	B
IM-15	2498	A	IM-51	6122	A	IM-87	2144	B	IM-123	3615	B
IM-16	2565	A	IM-52	6192	A	IM-88	2130	B	IM-124	3714	B
IM-17	2574	A	IM-53	6452	A	IM-89	2168	B	IM-125	3730	B
IM-18	2586	A	IM-54	6680	A	IM-90	2184	B	IM-126	3834	B
IM-19	2605	A	IM-55	7788	A	IM-91	2200	B	IM-127	3899	B
IM-20	2722	A	IM-56	8145	A	IM-92	2284	B	IM-128	3969	B
IM-21	2746	A	IM-57	8954	A	IM-93	2299	B	IM-129	3986	B
IM-22	2788	A	IM-58	9312	A	IM-94	2314	B	IM-130	3997	B
IM-23	2855	A	IM-59	9449	A	IM-95	2414	B	IM-131	4013	B
IM-24	2871	A	IM-60	10272	A	IM-96	2428	B	IM-132	4181	B
IM-25	2984	A	IM-61	11663	A	IM-97	2451	B	IM-133	4297	B
IM-26	3030	A	IM-62	13376	A	IM-98	2468	B	IM-134	4311	B
IM-27	3144	A	IM-63	14698	A	IM-99	2483	B	IM-135	4465	B
IM-28	3243	A	IM-64	15190	A	IM-100	2565	B	IM-136	4484	B
IM-29	3273	A	IM-65	15951	A	IM-101	2583	B	IM-137	4579	B
IM-30	3290	A	IM-66	15712	A	IM-102	2597	B	IM-138	4608	B
IM-31	3368	A	IM-67	15925	A	IM-103	2697	B	IM-139	4669	B
IM-32	3445	A	IM-68	23436	A	IM-104	2715	B	IM-140	4747	B
IM-33	3483	A	IM-69	39794	A	IM-105	2740	B	IM-141	4862	B
IM-34	3678	A	IM-70	44166	A	IM-106	2762	B	IM-142	4891	B
IM-35	3779	A	IM-71	46890	A	IM-107	2767	B	IM-143	5033	B
IM-36	3793	A				IM-108	2805	B	IM-144	5077	B

(57) Abstract: Certain biomarkers and biomarker combinations are useful in a qualifying lung cancer status in a subject. A diagnostic methodology employing these biomarkers and combinations can detect whether a subject has lung cancer.

WO 2004/061410 A2

BEST AVAILABLE COPY



European patent (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IT, LU, MC, NL, PT, RO, SE, SI, SK, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).

*For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.*

**Published:**

- *without international search report and to be republished upon receipt of that report*

## SERUM BIOMARKERS IN LUNG CANCER

### BACKGROUND OF THE INVENTION

[0001] The present invention relates generally to the field of serum biomarkers in lung carcinoma. More particularly, the invention relates to serum biomarkers that can distinguish lung cancer from normal.

[0002] Lung cancer is the leading cause of cancer death worldwide, resulting in 150,000 deaths per year in the United States. The mortality rate from lung cancer is greater than the combined mortality from breast, prostate and colorectal cancers. On the basis of morphology, lung cancer can be broadly classified into four main categories namely, adenocarcinoma, squamous cell carcinoma, large cell undifferentiated carcinoma and small cell carcinoma. In Hong Kong from 1990 to 1996, the proportions for adenocarcinoma, squamous cell carcinoma, large cell undifferentiated carcinoma and small cell carcinoma are 45.5%, 27.5%, 4.7% and 10.3% respectively. Both squamous cell carcinoma and small cell carcinoma are strongly associated with a smoking history.

[0003] Adenocarcinoma, squamous cell carcinoma, and large cell undifferentiated carcinoma are usually referred as "non-small cell carcinoma." They are relatively chemo-resistant, and hence the mainstay of treatment is surgery. By contrast, small cell carcinoma has a higher propensity for distant metastases and is mainly treated by chemotherapy.

[0004] Biopsy can be used to diagnose lung cancer, but it is an invasive procedure and, therefore, less than desirable. Other diagnostic methods for lung cancer include ultrasound and computed tomography (CT) scan.

[0005] It would be highly desirable to have a biomarker or combination of biomarkers capable of distinguishing between lung cancer and normal cells. In addition, a simple test could aid in tracking treatment progress and even identify molecular targets for therapy. The literature on lung cancer diagnosis has not disclosed heretofore such a biomarker or combination of biomarkers, however.

## SUMMARY OF THE INVENTION

[0006] In accordance with the present invention, biomarkers and combinations of biomarkers are used to identify lung cancer. The method successfully distinguishes between lung cancer and normal states, and can be used to identify the particular type of lung cancer. In one embodiment, a method for qualifying lung carcinoma status in a subject (e.g., a patient) comprises analyzing a biological sample from the subject for one or more of the top 50 biomarkers as shown in Figure 2 or Figures 4A and 4B. Thus, to assess overall lung cancer risk versus normal, a biomarker is selected from the group consisting of

(A) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-157, IM-176, IM-445, IM-177, IM-440, IM-468, IM-438, IM-547, IM-359, IM-436, IM-106, IM-455, IM-444, IM-158, IM-265, IM-50, IM-159, IM-156, IM-439, IM-157, IM-508, IM-514, IM-478, IM-473, IM-360, IM-435, IM-150, IM-151, IM-110, IM-51, IM-163, IM-437, IM-546, IM-153, and IM-268, or

(B) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70, WM-369, WM-17, WM-473, WM-47, WM-203, WM-276, WM-279, WM-62, WM-366, WM-456, WM-428, WM-384, WM-287, WM-420, WM-292, WM-431, WM-455, WM-20, WM-340, WM-105, WM-389, WM-63, WM-354, WM-450, WM-466, WM-296, WM-343, WM-341, WM-339, WM-55, WM-66, WM-48, WM-38, WM-138, and WM-310,

[0007] wherein the biomarker is differentially present in samples of a subject with lung cancer and a so-called "normal" subject that is free of lung cancer.

[0008] More preferably, one or more of the top 15 biomarkers as shown in Figure 2 or Figures 4A and 4B is used to qualify lung cancer status. Thus, for assessing overall lung cancer status versus normal, the protein is selected from the group consisting of

(A) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-471, IM-510, IM-544, IM-474, and IM-155, or

(B) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70.

[0009] Still more preferably, one or more of the top 5 biomarkers as shown in Figure 2 or Figures 4A and 4B is used to qualify lung cancer status. In this instance, for overall lung cancer status versus normal, the biomarker is selected from the group consisting of

(A) IM-522, IM-273, IM-520, IM-519, and IM-454, or

(B) WM-61, WM-447, WM-446, WM-133, and WM-119.

[0010] In one embodiment, the method measures a plurality of biomarkers. The plurality of biomarkers can be measured simultaneously.

[0011] Biomarkers that, by themselves, are able to identify lung cancer include the WM-446 and WM-447 protein biomarkers, and these are particularly preferred.

[0012] The present invention also provides a method for qualifying lung cancer status in a subject (e.g., a patient), comprising (A) providing a spectrum generated by subjecting a biological sample from said subject to mass spectroscopic analysis that includes profiling on a chemically-derivatized affinity surface, and (B) putting the spectrum through pattern-recognition analysis that is keyed to at least one peak selected from the top 50 biomarkers as shown in Figure 2 or Figures 4A and 4B.

Thus, for qualifying overall lung cancer status, the biomarker is selected from the group consisting of

(i) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-157, IM-176, IM-445, IM-177, IM-440, IM-468, IM-438, IM-547, IM-359, IM-436, IM-106, IM-455, IM-444, IM-158, IM-265, IM-50, IM-159, IM-156, IM-439, IM-157, IM-508, IM-514, IM-478, IM-473, IM-360, IM-435, IM-150, IM-151, IM-110, IM-51, IM-163, IM-437, IM-546, IM-153, and IM-268 or

(B) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70, WM-369, WM-17, WM-473, WM-47, WM-203, WM-276, WM-279, WM-62, WM-366, WM-456, WM-428, WM-384, WM-287, WM-420, WM-292, WM-431, WM-455, WM-20, WM-340, WM-105, WM-389, WM-63, WM-354, WM-450, WM-466, WM-296, WM-343, WM-341, WM-339, WM-55, WM-66, WM-48, WM-38, WM-138, and WM-310.

[0013] For assessing the overall lung cancer status, the pattern-recognition analysis may, for example, be paired to a pair of peaks selected from the group consisting of

(A) IM-266 and IM-474, IM-266 and IM-38, IM-266 and IM-454, IM-266 and IM-522, IM-266 and IM-544, IM-266 and IM-471, IM-474 and IM-151, IM-474 and IM-156, IM-474 and IM-544, IM-474 and IM-38, IM-522 and IM-507, IM-522 and IM-156, and IM-522 and IM-440;

or

(B) WM-447 and WM-59, WM-447 and WM-19, WM-447 and WM-118, WM-447 and WM-473, WM-19 and WM-59, WM-19 and WM-473, WM-19 and WM-369, WM-61 and WM-154, WM-61 and WM-369, WM-118 and WM-59 and WM-282 and WM-127.

[0014] More preferably, for assessing overall lung cancer status, the pattern-recognition analysis is keyed to a pair of peaks selected from the group consisting of

(A) IM-266 and IM-474, IM-266 and IM-544, and IM-156 and IM-522;

or

(B) WM-447 and WM-59, WM-447 and WM-19, and WM-19 and WM-59.

[0015] Alternatively, the pattern-recognition analysis for assessing overall lung cancer status may be keyed to a triplet of peaks selected from the group consisting of

(A) IM-266, IM-454 and IM-474; and IM-266, IM-474 and IM-544;

or

(B) WM-447, WM-19 and WM-473.

[0016] In other embodiments, the pattern-recognition analysis may be keyed to a combination of more than three peaks, more particularly to a combination of 4, 5 or 6 peaks, where the combination is selected from among the combinations shown in Tables 1 and 2 herein.

[0017] In each case, the biomarker is differentially present in samples of a subject with lung cancer and a normal subject.

[0018] The invention also contemplates a kit for detecting and diagnosing lung cancer, thereby to assess lung cancer status. Kits within the invention comprise, for example, (i) an adsorbent attached to a substrate that retains one or more of the biomarkers shown in Figure 2 or Figures 4A and 4B, and (ii) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent. An inventive kit may further comprise a washing solution and/or instructions for making a washing solution. The kits may include more than type of adsorbent, each present on a different substrate, *e.g.*, on a WCX and IMAC biochip. In addition, the kits may comprise one or more containers with biomarker samples, to be used as standard(s) for calibration. The substrate comprising the adsorbent may be designed to engage a probe interface and, hence, function as a probe in gas phase ion spectrometry, preferably mass spectrometry. Alternatively, the kit may further comprise a second substrate adapted to engage the probe interface, on which the substrate comprising the adsorbent is mounted.

[0019] The method and kit according to the invention produce an article of manufacture in which one or more biomarkers according to the invention are bound to an adsorbent, optionally contacted with a matrix or energy absorbing molecule.

[0020] The present invention also provides software for qualifying lung carcinoma status in a subject, comprising an algorithm for analyzing data extracted from a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, wherein said data relates to one or more biomarkers according to the invention. In one embodiment, the algorithm carries out a pattern-recognition analysis that is keyed to data relating to at least one of the biomarkers. In another embodiment, the algorithm comprises classification tree analysis that is keyed to data relating to at least one of the biomarkers. In yet another embodiment, the algorithm

comprises an artificial neural network analysis that is keyed to data relating to at least one of the biomarkers.

[0021] In certain embodiments, the present invention provides methods and kits that use serum amyloid a protein or a fragment thereof to qualify lung carcinoma status in a subject. In one of these embodiments, the serum amyloid a biomarker has an apparent molecular weight of about 2803, 3168, 3277, 3552, 3897, 4300, 4490, 4655, 5927, 6874, 7776, 7941, 8152, 8952, 9233, 10300, 10866, or 10851 Daltons. In another embodiment, the serum amyloid a biomarker has an apparent molecular weight of about 3168, 3277, 3552, 3897, 4300, 4490, 4655, 7776, 7941, 8152, 8952, or 10851 Daltons. In yet another embodiment, the serum amyloid a biomarker has an apparent molecular weight of about 11.5 to 11.7 kD.

#### BRIEF DESCRIPTION OF THE DRAWINGS

[0022] Figures 1A-1D show all biomarkers identified with a Cu(II) IMAC3 ProteinChip® array format.

[0023] Figure 2 shows the top 50 biomarkers identified with a Cu(II) IMAC3 ProteinChip® array format.

[0024] Figures 3A-3O show all biomarkers identified with a WCX ProteinChip® array format.

[0025] Figures 4A and 4B show the top 50 biomarkers identified with a WCX ProteinChip® array format.

[0026] Figure 5 shows fragments of serum amyloid A (SAA) that are biomarkers according to the present invention.

[0027] Figure 6 shows identification of SAA biomarkers with an anti-SAA antibody.

[0028] Figures 7-16 are spectra from WCX chips in which all of the top 15 WCX marker peaks are labeled, along with various other peaks from among the top 50 WCX peaks. Red shows spectra from lung cancer patients and gray shows normals.

[0029] Figures 17-28 are spectra from IMAC chips in which all of the top 15 WCX marker peaks are labeled, along with various other peaks from among the top 50 IMAC peaks. Blue shows spectra from lung cancer patients and gray shows normals.



## DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENTS

[0030] In accordance with the present invention, a series of biomarkers associated with lung cancer has been discovered. In the present context, a biomarker is an organic biomolecule, particularly a polypeptide or protein, which is differentially present in a sample taken from a subject having lung cancer as compared to a comparable sample taken from a normal subject. A biomarker also may be differentially present in a sample taken from a subject with one type of lung cancer, *e.g.*, small cell carcinoma, as compared to a comparable sample taken from a subject with a different type of lung cancer, *e.g.*, adenocarcinoma or squamous cell carcinoma, or differentially present at different stages of a type of lung cancer. A biomarker is differentially present in samples taken from two groups of subjects if it is present at an elevated level or a decreased level in samples of the first group as compared to samples of the second group. More particularly, a biomarker is a polypeptide that is characterized by an apparent molecular weight, as determined by mass spectrometry, and that is present in samples from lung cancer subjects in an elevated or decreased level, as compared to subjects that do not have lung cancer. A biomarker is differentially present between two sets of samples if the amount of the biomarker in one sample set differs in a statistically significant way ( $p < 0.01$ ) from the amount of biomarker in the other sample set.

[0031] The biomarkers of the invention can be used to assess lung cancer status in a subject. For example, they are capable of identifying lung cancer and successfully distinguishing it from normal subjects, thereby providing a way of diagnosing the presence or absence of lung cancer, including the presence or absence of a particular kind of lung cancer. In addition, the biomarkers are useful in assessing the risk of developing lung cancer, in staging of lung cancer and in assessing the effectiveness of treatment. Thus, "lung cancer status" in the context of the present invention includes, *inter alia*, the presence or absence of disease, the risk of developing disease, the stage of the disease, and the effectiveness of treatment of disease. Based on this status, further procedures may be indicated, including additional diagnostic tests or therapeutic procedures or regimens, such as endoscopy, biopsy, surgery, chemotherapy, immunotherapy, and radiation therapy.

[0032] In some instances, a single biomarker is capable of identifying lung cancer with a sensitivity or specificity of at least 85%, whereas, in other instances, a combination or plurality of biomarkers is used to obtain a sensitivity or specificity of at least 85%. The biomarkers and combinations of biomarkers thus can be used to qualify lung cancer status in a subject or patient.

[0033] The biomarkers according to the invention are present in serum. The biological sample used according to the present invention, however, need not be a serum sample. Thus, a biological sample for qualifying lung cancer status may be a serum, plasma or blood sample, although serum samples are preferred.

[0034] All of the biomarkers are characterized by molecular weight. A list of all the biomarkers obtained with the Cu(II) IMAC3 ProteinChip® array (Ciphergen Biosystems, Inc., Fremont, California, USA) is provided in Figures 1A-1D, and Figure 2 lists the top 50 biomarkers that distinguish between lung cancer and normal subjects that are identified by Cu(II) IMAC3 protocol described herein. Figures 3A-3O comprise a list of all the biomarkers obtained with the WCX2 ProteinChip® array, and Figures 4A and 4B comprise a ranking of the top 50 biomarkers that distinguish between (i) lung cancer and normal subjects, (ii) subjects with each of four types of lung cancer and normal subjects, and (iii) two types of lung cancer, e.g., adenocarcinoma versus squamous cell carcinoma, as identified by WCX2 protocol described herein.

[0035] The top 50 biomarkers were determined by decision tree analysis using Biomarker Patterns™ software from Ciphergen Biosystems, Inc. Biomarkers other than those within the top 50 also are useful in distinguishing between subjects with lung cancer and normal subjects and may, in particular, appear in decision trees with multiple nodes. In preferred embodiments, one or more of the top 15 biomarkers are used, and in even more preferred embodiments, one or more of the top 5 biomarkers are used.

[0036] In each of Figures 1A-1D and 3A-3O, the number in the first column is the biomarker identifier. Thus, the first row in Figures 1A-1D relates to biomarker IM-1, the second row relates to biomarker IM-2, and so forth ("IM-" denoting biomarkers identified with the IMAC chip). Similarly, the first row in Figures 3A-3O relates to

biomarker WM-1 and the second row relates to biomarker WM-2 ("WM-" denoting biomarkers identified with the WCX2 chip). The number in the second column in Figures 1A-1D is the apparent molecular weight of the biomarker in daltons, as determined by mass spectrometry. In Figures 3A-3O, the apparent molecular weights for the biomarkers identified in the first column are reported in columns 3 through 11. The letter in the second column of Figures 1A-1D and the third column of Figures 3A-3O denotes the fraction in which the biomarker elutes in the protocol described herein; that is, biomarkers with an "A" elute in the first fraction, biomarkers with a "B" elute in the second fraction, and so forth. The fraction in which the biomarker elutes correlates with its pI, which biomarkers eluting at higher pH having a higher pI, and biomarkers eluting at lower pH having a lower pI.

[0037] Presenting the mass and affinity characteristics of a given biomarker within the invention, as in this description, characterizes that biomarker so as allow one to obtain and measure it, in accordance with the teachings herein. If desired, any of the biomarkers can be sequenced, in order to obtain an amino acid sequence, but this is not required to practice the present invention.

[0038] For example, a biomarker can be peptide mapped with a number of enzymes, such as trypsin and V8 protease, and the molecular weights of the digestion fragments can be used to search databases for sequences that match the molecular weights of the digestion fragments generated by the various enzymes. Alternatively, if the biomarkers are not proteins included in known databases, degenerate probes can be made based on the N-terminal amino acid sequence of the biomarker, which then are used to screen a genomic or cDNA library created from a sample from which the biomarker was initially detected. The positive clones can be identified, amplified, and their recombinant DNA sequences can be subcloned using techniques which are well known. Finally, protein biomarkers can be sequenced using protein ladder sequencing. Protein ladders can be generated by fragmenting the molecules and subjecting fragments to enzymatic digestion or other methods that sequentially remove a single amino acid from the end of the fragment. The ladder is then analyzed by mass spectrometry. The difference in masses of the ladder fragments identifies the amino acid removed from the end of the molecule.

[0039] Several biomarkers identified in accordance with the teachings of the present invention fit to serum amyloid A (SAA) or to a fragment of SAA. SAA is a well-known acute phase inflammatory marker. A number of the SAA biomarkers are identified in Figure 5 by both molecular mass and amino acid sequence. Most of these markers bound anti-SAA antibodies, as shown in Figure 6. The intact mass of SAA is 11.5 to 11.7 kD, and these biomarkers also have been identified by the present methodology. Fragments preferably have a molecular mass of at least about 200 Daltons, more preferably at least about 500 Daltons. In even more preferred embodiments, fragments have a molecular mass of at least about 800 Daltons, and most preferably at least about 1 Kilodalton.

[0040] In one embodiment, the fragments of SAA include a sequence of amino acids that is recognized by an epitope of an anti-SAA antibody. One way of identifying suitable fragments for use in the present invention is to enzymatically digest SAA and test the resulting fragments for the ability to bind to an anti-SAA antibody. Fragments that bind anti-SAA antibody can be sequenced using techniques well-known in the art, although the sequence of the fragment is not needed to practice the invention. In order to practice the invention with a fragment from the enzymatic digest that is identified as binding anti-SAA antibody, all that is required is to subject to the fragment to mass spectrometry to determine its mass.

[0041] The serum biomarkers according to the present invention were identified by comparing mass spectra of samples derived from sera from two groups of newly-diagnosed subjects, subjects with lung cancer and normal subjects. The subjects were diagnosed according to standard clinical criteria. Lung cancer subjects were histologically confirmed, and subjects without lung cancer were followed for at least 18 months following serum collection for any sign of lung cancer, to exclude subjects with asymptomatic lung cancer.

[0042] Sera from each group of subjects was collected, and fractionated with Q Ceramic HyperDF ion exchange resin (Biosepra SA, France) into six fractions which eluted at different pH. Fraction A comprised the flow through plus pH 9 eluant, Fraction B comprised the pH 7 eluant, Fraction C comprised the pH 5 eluant, Fraction D comprised the pH 4 eluant, Fraction E comprised the pH 3 eluant, and Fraction F

comprised isopropyl alcohol/acetonitrile TFA eluant. Fractions A through F are identified on Figures 7-28 as Fractions 1 through 6, respectively.

[0043] Each fraction was diluted and applied to a ProteinChip® array, either a Cu(II) IMAC3 or WCX2 chip array. Both of these chip arrays are produced by CIPHERGEN Biosystems, Inc. (Fremont, CA).

[0044] The Cu(II) IMAC3 is an "immobilized metal affinity-capture" chip, with a nitrilotriacetic acid surface for high-capacity copper binding and subsequent affinity capture of proteins with metal binding residues. Imidazole may be used in binding and washing solutions to moderate protein binding, including binding of non-specific proteins. Increasing the concentration of imidazole in the washing buffers reduces the binding of the target proteins. It is produced by photopolymerizing 5-methylacrylamido-2-(N,N-bis(carboxymethyl)amino)pentanoic acid (7.5 wt%) and N,N'-methylenebisacrylamide (0.4 wt%) using (-) riboflavin (0.02 wt%) as a photoinitiator. The monomer solution is deposited onto the chip substrate and irradiated to photopolymerize. The chip then is activated with Cu(II).

[0045] The WCX2 is a weak cation exchange array with a carboxylate surface to bind cationic proteins. The negatively charged carboxylate groups on the surface of the WCX2 chip interact with the positive charges exposed on the target proteins. The binding of the target proteins is reduced by increasing the concentration of salt or by increasing the pH of the washing buffers.

[0046] Following application of the eluant fraction, the chips were incubated to allow the polypeptides in the eluant to bind to the sites on the chip by an affinity interaction. After incubation, each chip array was washed to remove polypeptides that bind non-specifically and buffer contaminants. That chip then was dried, and an energy absorbing molecule or matrix was applied to it, to facilitate desorption and ionization in a mass spectrometer.

[0047] In the mass spectrometer, retained polypeptides were desorbed from the chip array by laser desorption and ionization in a ProteinChip® Reader, which is integrated with ProteinChip® Software and a personal computer to analyze proteins captured on chip arrays. The ion optic and laser optic technologies in the ProteinChip® Reader detects proteins ranging from small peptides of less than 1000 Da up to proteins of

300 kilodaltons or more, and calculates the mass based on time-of-flight. Ionized polypeptides were detected and their mass accurately determined by this Time-of-Flight (TOF) Mass Spectrometry.

[0048] The mass spectra obtained for each group were subjected to scatter plot analysis, to eliminate run-to-run variation. Protein clusters on the scatter plot that had the same pattern for both lung cancer and normal subjects, *i.e.*, protein clusters that were either elevated in both groups of subjects or depressed in both groups of subjects, were eliminated as potential biomarkers. The remaining polypeptides were further analyzed for their ability to accurately identify subjects with lung cancer. Because the molecular weights were derived from scatter plot analysis, and because of limits on the ability of mass spectrometry to resolve molecular weights, the "absolute" molecular weight values given in Figures 1A-1D and 3A-3O actually represent approximate molecular weights.

[0049] The biomarkers of this invention are characterized by their mass-to-charge ratio as determined by mass spectrometry. The mass-to-charge ratio of each biomarker is provided in Figures 1A-1D and 3A-3O. For example, IM-1 in Figure 1A has a measured mass-to-charge ratio of 2011. The mass-to-charge ratios were determined from mass spectra generated on a Ciphergen Biosystems, Inc. PBS II mass spectrometer. This instrument has a mass accuracy of about +/- 0.15 percent. Additionally, the instrument has a mass resolution of about 400 to 1000 m/dm, where m is mass and dm is the mass spectral peak width at 0.5 peak height. The mass-to-charge ratio of the biomarkers was determined using Biomarker Wizard™ software (Ciphergen Biosystems). Biomarker Wizard assigns a mass-to-charge ratio to a biomarker by clustering the mass-to-charge ratios of the same peaks from all the spectra analyzed, as determined by the PBSII, taking the maximum and minimum mass-to-charge-ratio in the cluster, and dividing by two. Accordingly, the masses provided reflect these specifications.

[0050] The biomarkers of this invention are further characterized by the shape of their spectral peak in time-of-flight mass spectrometry. Mass spectra showing peaks representing the biomarkers are presented in Figures 7-28. The biomarker identifier numbers from Figures 2 and 4A-4B, respectively, are shown next to the peak, along

with their rank, which is indicated in parentheses below the biomarker identifier number.

[0051] The biomarkers of this invention are further characterized by their binding properties on chromatographic surfaces. Most of the biomarkers bind to IMAC (Cu) or WCX adsorbents (e.g., the CIPHERGEN® IMAC (Cu) or WCX ProteinChip® arrays) after washing as described herein.

[0052] Thus, a given molecular weight for a biomarker herein should be interpreted as the midpoint of a molecular-weight range. The accuracy of the mass spectrometer is  $\pm 0.15\%$ , and the actual molecular weight for a biomarker is therefore the value given,  $\pm 0.15\%$ . For example, the actual molecular weight for biomarker IM-273 is  $11705 \pm 0.15\%$ , or between 11687 and 11722. Often, the range surrounding the "absolute" value given in the figure is no more than  $\pm 5$  daltons (2006 to 2016 for IM-1), generally no more than  $\pm 3$  daltons (2008 to 2014 for IM-1), and often as small as  $\pm 1$  dalton (2010 to 2012 daltons for IM-1).

[0053] CART® (Salford Systems, San Diego, CA), a classification and regression tree software, was used to determine whether a potential biomarker had predictive value in assessing lung cancer. A software macro randomly selected a subset of 15% of the peaks from Figures 1A-1D or Figures 3A-3O. The peaks and peak heights from each sample were provided to the CART® software for analysis. The software performed an iterative analysis until a single decision tree was generated that was capable of distinguishing between cancerous and non-cancerous. Each node in the resulting decision tree sorted based on the peak height of a single biomarker. A tree may contain any number of nodes, but generally contains from 1 to 6 nodes. From a practical standpoint in a commercial diagnostic test, a decision tree with fewer nodes is preferred. A total of 2000 decision trees, each based on a different 15% subset of the peaks from Figures 1A-1D or Figures 3A-3O, were generated.

[0054] The CART® software assigned a score to each biomarker in the subset, based on its relative importance. A score of 100 is very high and a score of 0 is very low. The CART® software also determined the sensitivity and specificity of each decision tree.

[0055] The data generated by the decision tree analysis was subjected to further analysis. The biomarkers were ranked based on their average scores, which were determined by adding up a biomarker's scores for each decision tree in which it appeared, and dividing by the total number of decision trees in which the biomarker appeared. Approximately 500 of the potential biomarkers showed up in at least one tree, and most of the biomarkers showed up in about 150 to 400 of the two thousand trees. The top 50 biomarkers for the IMAC and WCX chip arrays as determined by this method are shown in Figures 2 and 4A-4B, respectively.

[0056] All of the trees having sensitivities and specificities greater than 85% also were identified. All trees capable of distinguishing lung cancer from normal and having from 1 to 6 nodes that meet the 85/85 criterion are shown in Tables 1 and 2.

**TABLE 1. Decision trees with IMAC Biomarkers.**

2 Nodes				
474	151			
474	156			
522	507		2 trees	
522	440		2 trees	
3 Nodes				
266	454	474		
474	156	153		
474	40	156		
520	276	113		
520	265	401		
522	151	474		
522	478	153		
522	156	474		
4 Nodes				
148	521	508	251	



266	544	474	493	
266	157	126	420	
266	544	474	482	
266	471	474	38	
266	544	474	38	
266	514	471	203	
522	58	266	474	
5 Nodes				
266	544	473	151	437
266	454	474	153	264
273	143	544	401	199

TABLE 2. Decision Trees with WCX Biomarkers.

1 Node					
446					
447					
2 Nodes					
282	127				
3 Nodes					
61	16	27			
61	119	154			
61	120	154			
61	369	184			
61	184	129			
61	19	282			
133	282	319			
282	59	218			
282	111	65			

446	19	16			
<b>4 Nodes</b>					
61	369	282	184		
61	48	203	3		
446	369	111	67		
446	466	58	120		
446	19	59	113		
446	282	19	47		
447	118	59	417		
447	118	59	473		
447	65	59	275		
447	19	59	282		
447	369	59	206		
447	19	59	253		
447	19	47	70		
<b>5 Nodes</b>					
61	369	128	184	197	
61	17	425	366	341	
133	139	363	216	273	
282	133	48	19	253	
369	310	19	109	384	
446	282	15	319	66	
447	19	71	473	31	
447	19	17	473	438	
447	47	31	365	59	
<b>6 Nodes</b>					
369	366	192	471	19	439

[0057] Each of the biomarker combinations of Tables 1 and 2 are preferred combinations for distinguishing lung cancer subjects from normal subjects in accordance with the present invention.

[0058] All biomarkers that appeared in at least two of the trees that met the 85/85 criterion were identified. For these biomarkers, Tables 3 and 4 provide the number of times the biomarker occurred in a trees that met the criterion, as well as the ranking of that biomarker on the top 50 lists of Figures 2 and 4A-4B.

**TABLE 3. Correlation of IMAC biomarker decision tree frequencies and ranking.**

Peak		# times		Rank
266		9		9
522		8		1
474		4		14
520		2		3
148		1		8
273		1		2

**TABLE 4. Correlation of WCX biomarker decision tree frequencies and ranking.**

Peak		# times		Rank
447		11		2
61		10		1
446		7		3
282		4		9
369		2		8
133		2		4

[0059] Biomarkers that occurred frequently in the highly discriminatory trees occurred among the top 50 ranked biomarkers, and typically had a top 10 ranking. In addition, certain pairs of biomarkers reappear, *e.g.*, WM-447 and WM-59, WM-447 and WM-19, WM-19 and WM-59, IM-266 and IM-474, IM-266 and IM-38, IM-266 and IM-454, IM-522 and IM-266. There also are repeats among triplets of biomarkers, such as IM-266, IM-266 and IM-38, and WM-447, WM-19 and WM-473. Other repeating pairs and trios of biomarkers can be seen in Tables 3 and 4, and are preferred.

[0060] Biomarkers and combinations of biomarkers identified in accordance with the present description may be used to qualify lung cancer status in a subject. In particular, a biomarker or combination of biomarkers can be used to distinguish lung cancer patients from normal patients with a high degree of specificity or sensitivity, *i.e.*, greater than at least 85%, preferably greater than at least 90%, and more preferably greater than 95%.

[0061] According to one aspect of the invention, therefore, the detection of biomarkers for diagnosis of lung cancer status entails contacting a sample from a subject with a substrate, *e.g.*, a SELDI probe, having an adsorbent thereon, under conditions that allow binding between the biomarker and the adsorbent, and then detecting the biomarker bound to the adsorbent by gas phase ion spectrometry, for example, mass spectrometry. Other detection paradigms that can be employed to this end include optical methods, electrochemical methods (voltametry and amperometry techniques), atomic force microscopy, and radio frequency methods, *e.g.*, multipolar resonance spectroscopy. Illustrative of optical methods, in addition to microscopy, both confocal and non-confocal, are detection of fluorescence, luminescence, chemiluminescence, absorbance, reflectance, transmittance, and birefringence or refractive index (*e.g.*, surface plasmon resonance, ellipsometry, a resonant mirror method, a grating coupler waveguide method or interferometry).

[0062] In one aspect, the markers of this invention are detect by gas phase ion spectrometry, which refers to the use of a gas phase ion spectrometer to detect gas phase ions. A gas phase ion spectrometer is an apparatus that detects gas phase ions. Gas phase ion spectrometers include an ion source that supplies gas phase ions. Gas

phase ion spectrometers include, for example, mass spectrometers, ion mobility spectrometers, and total ion current measuring devices.

[0063] "Mass spectrometer" refers to a gas phase ion spectrometer that measures a parameter which can be translated into mass-to-charge ratios of gas phase ions. Mass spectrometers generally include an ion source and a mass analyzer. Examples of mass spectrometers are time-of-flight, magnetic sector, quadrupole filter, ion trap, ion cyclotron resonance, electrostatic sector analyzer and hybrids of these. "Mass spectrometry" refers to the use of a mass spectrometer to detect gas phase ions. "Laser desorption mass spectrometer" refers to a mass spectrometer which uses laser as a means to desorb, volatilize, and ionize an analyte.

[0064] "Mass analyzer" refers to a sub-assembly of a mass spectrometer that comprises means for measuring a parameter which can be translated into mass-to-charge ratios of gas phase ions. In a time-of flight mass spectrometer the mass analyzer comprises an ion optic assembly, a flight tube and an ion detector.

[0065] "Ion source" refers to a sub-assembly of a gas phase ion spectrometer that provides gas phase ions. In one embodiment, the ion source provides ions through a desorption/ionization process. Such embodiments generally comprise a probe interface that positionally engages a probe in an interrogatable relationship to a source of ionizing energy (e.g., a laser desorption/ionization source) and in concurrent communication at atmospheric or subatmospheric pressure with a detector of a gas phase ion spectrometer.

[0066] Forms of ionizing energy for desorbing/ionizing an analyte from a solid phase include, for example: (1) laser energy; (2) fast atoms (used in fast atom bombardment); (3) high energy particles generated via beta decay of radionuclides (used in plasma desorption); and (4) primary ions generating secondary ions (used in secondary ion mass spectrometry). The preferred form of ionizing energy for solid phase analytes is a laser (used in laser desorption/ionization), in particular, nitrogen lasers, Nd-Yag lasers and other pulsed laser sources. "Fluence" refers to the laser energy delivered per unit area of interrogated image. Typically, a sample is placed on the surface of a probe, the probe is engaged with the probe interface and the probe

surface is struck with the ionizing energy. The energy desorbs analyte molecules from the surface into the gas phase and ionizes them.

[0067] Other forms of ionizing energy for analytes include, for example: (1) electrons which ionize gas phase neutrals; (2) strong electric field to induce ionization from gas phase, solid phase, or liquid phase neutrals; and (3) a source that applies a combination of ionization particles or electric fields with neutral chemicals to induce chemical ionization of solid phase, gas phase, and liquid phase neutrals.

[0068] A preferred mass spectrometric technique for use in the invention is Surface Enhanced Laser Desorption and Ionization (SELDI), as described, for example, in U.S. patents No. 5,719,060 and No. 6,225,047, both to Hutchens and Yip, in which the surface of a probe that presents the analyte (here, one or more of the biomarkers) to the energy source plays an active role in desorption/ionization of analyte molecules. In this context, "probe" refers to a device adapted to engage a probe interface and to present an analyte to ionizing energy for ionization and introduction into a gas phase ion spectrometer, such as a mass spectrometer. A probe typically includes a solid substrate, either flexible or rigid, that has a sample-presenting surface, on which an analyte is presented to the source of ionizing energy.

[0069] One version of SELDI, called Surface-Enhanced Affinity Capture" or "SEAC," involves the use of probes comprised of a chemically selective surface ("SELDI probe"). A "chemically selective surface" is one to which is bound either the adsorbent, also called a "binding moiety" or "capture reagent," or a reactive moiety that is capable of binding a capture reagent, *e.g.*, through a reaction forming a covalent or coordinate covalent bond.

[0070] The phrase "reactive moiety" here denotes a chemical moiety that is capable of binding a capture reagent. Epoxide and carbodiimidazole are useful reactive moieties to covalently bind polypeptide capture reagents such as antibodies or cellular receptors. Nitriloacetic acid and iminodiacetic acid are useful reactive moieties that function as chelating agents to bind metal ions that interact non-covalently with histidine containing peptides. A "reactive surface" is a surface to which a reactive moiety is bound. An "adsorbent" or "capture reagent" can be any material capable of

binding a biomarker of the invention. Suitable adsorbents for use in SELDI, according to the invention, are described in U.S. patent No. 6,225,047, *supra*.

[0071] One type of adsorbent is a "chromatographic adsorbent," which is a material typically used in chromatography. Chromatographic adsorbents include, for example, ion exchange materials, metal chelators, immobilized metal chelates, hydrophobic interaction adsorbents, hydrophilic interaction adsorbents, dyes, simple biomolecules (e.g., nucleotides, amino acids, simple sugars and fatty acids), mixed mode adsorbents (e.g., hydrophobic attraction/electrostatic repulsion adsorbents).

"Biospecific adsorbent" is another category, for adsorbents that contain a biomolecule, e.g., a nucleotide, a nucleic acid molecule, an amino acid, a polypeptide, a polysaccharide, a lipid, a steroid or a conjugate of these (e.g., a glycoprotein, a lipoprotein, a glycolipid). In certain instances the biospecific adsorbent can be a macromolecular structure such as a multiprotein complex, a biological membrane or a virus. Illustrative biospecific adsorbents are antibodies, receptor proteins, and nucleic acids. A biospecific adsorbent typically has higher specificity for a target analyte than a chromatographic adsorbent.

[0072] Another version of SELDI is Surface-Enhanced Neat Desorption (SEND), which involves the use of probes comprising energy absorbing molecules that are chemically bound to the probe surface ("SEND probe"). The phrase "Energy absorbing molecules" (EAM) denotes molecules that are capable of absorbing energy from a laser desorption ionization source and, thereafter, contributing to desorption and ionization of analyte molecules in contact therewith. The EAM category includes molecules used in MALDI, frequently referred to as "matrix," and is exemplified by cinnamic acid derivatives, sinapinic acid (SPA), cyano-hydroxy-cinnamic acid (CHCA) and dihydroxybenzoic acid, ferulic acid, and hydroxyaceto-phenone derivatives. The category also includes EAMs used in SELDI, as enumerated, for example, by U.S. 5,719,060 and U.S. 60/351,971 (Kitagawa), filed January 25, 2002.

[0073] Another version of SELDI, called Surface-Enhanced Photolabile Attachment and Release (SEPAR), involves the use of probes having moieties attached to the surface that can covalently bind an analyte, and then release the analyte through breaking a photolabile bond in the moiety after exposure to light, e.g., to laser light.

For instance, see U.S. 5,719,060. SEPAR and other forms of SELDI are readily adapted to detecting a biomarker or biomarker profile, pursuant to the present invention.

[0074] The detection of the biomarkers according to the invention can be enhanced by using certain selectivity conditions, *e.g.*, adsorbents or washing solutions. The phrase “wash solution” refers to an agent, typically a solution, which is used to affect or modify adsorption of an analyte to an adsorbent surface and/or to remove unbound materials from the surface. The elution characteristics of a wash solution can depend, for example, on pH, ionic strength, hydrophobicity, degree of chaotropism, detergent strength, and temperature.

[0075] Pursuant to one aspect of the present invention, a sample is analyzed by means of a “biochip,” a term that denotes a solid substrate, having a generally planar surface, to which a capture reagent (adsorbent) is attached. Frequently, the surface of a biochip comprises a plurality of addressable locations, each of which has the capture reagent bound there. A biochip can be adapted to engage a probe interface and, hence, function as a probe in gas phase ion spectrometry preferably mass spectrometry. Alternatively, a biochip of the invention can be mounted onto another substrate to form a probe that can be inserted into the spectrometer.

[0076] A variety of biochips is available for the capture of biomarkers, in accordance with the present invention, from commercial sources such as CIPHERGEN Biosystems (Fremont, CA), PERKIN ELMER (Packard BioScience Company (Meriden CT), ZYOMYX (Hayward, CA), and PHYLOS (Lexington, MA). Exemplary of these biochips are those described in U.S. patents No. 6,225,047, *supra*, and No. 6,329,209 (Wagner *et al.*), and in PCT publications WO 99/51773 (Kuimelis and Wagner) and WO 00/56934 (Englert *et al.*).

[0077] More specifically, biochips produced by CIPHERGEN Biosystems have surfaces, presented on an aluminum substrate in strip form, to which are attached, at addressable locations, chromatographic or biospecific adsorbents. The surface of the strip is coated with silicon dioxide.

[0078] Illustrative of CIPHERGEN ProteinChip® arrays are biochips H4, SAX-2, WCX-2, and IMAC-3, which include a functionalized, cross-linked polymer in the



form of a hydrogel, physically attached to the surface of the biochip or covalently attached through a silane to the surface of the biochip. The H4 biochip has isopropyl functionalities for hydrophobic binding. The SAX-2 biochip has quaternary ammonium functionalities for anion exchange. The WCX-2 biochip has carboxylate functionalities for cation exchange. The IMAC-3 biochip has nitriloacetic acid functionalities that adsorb transition metal ions, such as  $\text{Cu}^{++}$  and  $\text{Ni}^{++}$ , by chelation. These immobilized metal ions, in turn, allow for adsorption of biomarkers by coordinate covalent bonding. Thus, Ciphergen's IMAC ProteinChip® arrays are sold with reactive moieties that become adsorbent upon the addition by the user of a metal solution.

[0079] In keeping with the above-described principles, a substrate with an adsorbent is contacted with the sample, containing serum, for a period of time sufficient to allow biomarker that may be present to bind to the adsorbent. In one embodiment of the invention, more than one type of substrate with adsorbent thereon is contacted with the biological sample. For example, a sample may be applied to both a WCX and an IMAC chip. This technique can allow for even more definitive assessment of cancer status. After the incubation period, the substrate is washed to remove unbound material. Any suitable washing solutions can be used; preferably, aqueous solutions are employed.

[0080] An energy absorbing molecule then is applied to the substrate with the bound biomarkers. As noted, an energy absorbing molecule is a molecule that absorbs energy from an energy source such as a laser, thereby assisting in desorption of biomarkers from the substrate. Exemplary energy absorbing molecules include, as noted above, cinnamic acid derivatives, sinapinic acid and dihydroxybenzoic acid. Preferably sinapinic acid is used.

[0081] The biomarkers bound to the substrates are detected in a gas phase ion spectrometer such as a time-of-flight mass spectrometer. The biomarkers are ionized by an ionization source such as a laser, the generated ions are collected by an ion optic assembly, and then a mass analyzer disperses and analyzes the passing ions. The detector then translates information of the detected ions into mass-to-charge

ratios. Detection of a biomarker typically will involve detection of signal intensity. Thus, both the quantity and mass of the biomarker can be determined.

[0082] Data generated by desorption and detection of biomarkers can be analyzed with the use of a programmable digital computer. The computer program analyzes the data to indicate the number of markers detected, and optionally the strength of the signal and the determined molecular mass for each biomarker detected. Data analysis can include steps of determining signal strength of a biomarker and removing data deviating from a predetermined statistical distribution. For example, the observed peaks can be normalized, by calculating the height of each peak relative to some reference. The reference can be background noise generated by the instrument and chemicals such as the energy absorbing molecule which is set as zero in the scale.

[0083] The computer can transform the resulting data into various formats for display. The standard spectrum can be displayed, but in one useful format only the peak height and mass information are retained from the spectrum view, yielding a cleaner image and enabling biomarkers with nearly identical molecular weights to be more easily seen. In another useful format, two or more spectra are compared, conveniently highlighting unique biomarkers and biomarkers that are up- or down-regulated between samples. Using any of these formats, one can readily determine whether a particular biomarker is present in a sample.

[0084] Software used to analyze the data can include code that applies an algorithm to the analysis of the signal to determine whether the signal represents a peak in a signal that corresponds to a biomarker according to the present invention. The software also can subject the data regarding observed biomarker peaks to classification tree or ANN analysis, to determine whether a biomarker peak or combination of biomarker peaks is present that indicates lung cancer status. Analysis of the data may be "keyed" to a variety of parameters that are obtained either directly or indirectly from the mass spectrometric analysis of the sample. These parameters include, but are not limited to, the presence or absence of one or more peaks, the height of one or more peaks, the log of the height of one or more peaks, and other arithmetic manipulations of peak height data.

[0085] In another aspect, the present invention provides kits for aiding in the diagnosis of lung cancer status, which kits are used to detect biomarkers according to the invention. The kits screen for the presence of biomarkers and combinations of biomarkers that are differentially present in samples from normal subjects and subjects with lung cancer.

[0086] In one embodiment, the kit comprises a substrate having an adsorbent thereon, wherein the adsorbent is suitable for binding a biomarker according to the invention, and a washing solution or instructions for making a washing solution, in which the combination of the adsorbent and the washing solution allows detection of the biomarker using gas phase ion spectrometry, e.g., mass spectrometry. The kit may include more than type of adsorbent, each present on a different substrate.

[0087] In another embodiment, a kit of the invention may include a first substrate, comprising an adsorbent thereon, and a second substrate onto which the first substrate is positioned to form a probe, which can be inserted into a gas phase ion spectrometer, e.g., a mass spectrometer. In another embodiment, an inventive kit may comprise a single substrate that can be inserted into the spectrometer.

[0088] In a further embodiment, such a kit can comprise instructions for suitable operational parameters in the form of a label or separate insert. For example, the instructions may inform a consumer how to collect the sample or how to wash the probe. In yet another embodiment the kit can comprise one or more containers with biomarker samples, to be used as standard(s) for calibration.

[0089] In a preferred embodiment, the detection of biomarkers for diagnosis of lung cancer in a subject entails contacting a sample from a subject or patient, preferably a serum sample, with a substrate having an adsorbent thereon under conditions that allow binding between the biomarker and the adsorbent, and then detecting the biomarker bound to the adsorbent by gas phase ion spectrometry, preferably by Surface Enhanced Laser Desorption/Ionization (SELDI) mass spectrometry. The biomarkers are ionized by an ionization source such as a laser. The generated ions are collected by an ion optic assembly and accelerated toward an ion detector. Ions that strike the detector generate an electric potential that is digitized by a high speed time-array recording device that digitally captures the analog signal. Ciphergen's

ProteinChip® system employs an analog-to-digital converter (ADC) to accomplish this. The ADC integrates detector output at regularly spaced time intervals into time-dependent bins. The time intervals typically are one to four nanoseconds long. Furthermore, the time-of-flight spectrum ultimately analyzed typically does not represent the signal from a single pulse of ionizing energy against a sample, but rather the sum of signals from a number of pulses. This reduces noise and increases dynamic range. This time-of-flight data is then subject to data processing. In Ciphergen's ProteinChip® software, data processing typically includes TOF-to-M/Z transformation, baseline subtraction, high frequency noise filtering. Thus, both the quantity and mass of the biomarker can be determined.

[0090] The detection of the biomarkers can be enhanced by using certain selectivity conditions, *e.g.*, adsorbents or washing solutions. In one embodiment, the same or similar selectivity conditions that were used to discover the biomarkers are used in the method of detecting the biomarker in the sample. For example, immobilized metal affinity capture chips such as the Cu(II) IMAC3 and weak cationic exchange chips such as the WCX2 chips are preferred as the adsorbents for biomarker detection. However, other adsorbents can be used, as long as they have the binding characteristics suitable for binding the biomarkers.

[0091] More particularly, armed with the information regarding the biomarkers identified herein, various methods can be used to recognize patterns of doublets, triplets, and higher combinations of biomarkers according to the invention. These methods take raw data regarding which peaks are present and their intensity and provide a differential diagnosis of lung cancer versus normal for a sample.

[0092] Thus, the process can be divided into the learning phase and the classification phase. In the learning phase, a learning algorithm is applied to a data set that includes members of the different classes that are meant to be classified, for example, data from a plurality of samples diagnosed as cancer and data from a plurality of samples assigned a negative diagnosis. The methods used to analyze the data include, but are not limited to, artificial neural network, support vector machines, genetic algorithm and self-organizing maps and classification and regression tree analysis. These methods are described, for example, in WO 01/31579, May 3, 2001

(Barnhill *et al.*); WO 02/06829, January 24, 2002 (Hitt *et al.*) and WO 02/42733, May 30, 2002 (Paulse *et al.*). The learning algorithm produces a classifying algorithm. The classifier is keyed to elements of the data, such as particular markers and particular intensities of markers, usually in combination, that can classify an unknown sample into one of the two classes. The classifier is ultimately used for diagnostic testing.

[0093] Software, both freeware and proprietary software, is readily available to analyze such patterns in data, and to devise additional patterns with any predetermined criteria for success. Those biomarkers which by themselves are predictive of a differential diagnosis of lung cancer versus normal do not require pattern recognition software to analyze the data.

[0094] The following examples are offered by way of illustration, and are not limiting.

#### **Example I. Fractionation of serum**

Buffers:

1. U9 (9M urea, 2% CHAPS, 50mM Tris-HCl pH9)
2. U1 (1M urea, 0.22% CHAPS, 50mM Tris-HCl pH9)
3. wash buffer 1: 50mM Tris-HCl with 0.1% n-octyl  $\beta$ -D-Glucopyranoside (OGP) pH9
4. wash buffer 2: 100mM sodium phosphate with 0.1% OGP pH7
5. wash buffer 3: 100mM sodium acetate with 0.1% OGP pH5
6. wash buffer 4: 100mM sodium acetate with 0.1% OGP pH4
7. wash buffer 5: 50mM sodium citrate with 0.1% OGP pH3
8. wash buffer 6: 33.3% isopropanol / 16.7% acetonitrile / 0.1% trifluoroacetic acid in water.

[0095] Thirty microliters of U9 buffer were added to 20 $\mu$ L of serum in a tube and were mixed at 4°C for 20 minutes. Ion exchange resin (Q Ceramic HyperDF ion exchange resin, Biosepra SA, France) was washed 3 times with 5 bed volumes of 50mM Tris-HCl pH9 and stored in 50% suspension. To each well of a 96-well filter plate (96-well Silent Screen filter plate, Lonrodvne membrane, 0.45 micron pore,

Nalge Nunc International, USA), 125  $\mu$ L of ion exchange resin (50% suspension) was added on a Biomek 2000 Automation Workstation (Beckman Coulter, Fullerton, CA), washed 3 times with 150 $\mu$ L U1 buffer, and vacuum dried. Urea-treated serum was transferred to each well of ion exchange resin. The serum tube was rinsed with 50 $\mu$ L of U1 buffer, which was also transferred to the corresponding well in filter plate. The filter plate was mixed on a platform shaker at 4°C for 30 minutes. Flow-through fraction was collected in a 96-well plate by vacuum suction (Fraction 1). Then, 100 $\mu$ L of wash buffer 1 was added to each well of filter plate and mixed for 10 minutes at room temperature. Eluant was collected into the same 96-well plate (Fraction 1). Resins in the filter plate were subsequently washed two times each with 100 $\mu$ L wash buffers 2, 3, 4, 5 and 6. Each eluant (total volume of 200 $\mu$ L) was collected in a 96-well plate (Fractions 2,3,4,5 and 6).

#### **Example 2. SELDI analysis of fractionated serum**

[0096] ProteinChip® Arrays were set up in 96-well bioprocessors. Buffer delivery and sample incubation were performed on a Biomek 2000 Automation Workstation. Each serum fraction was analyzed on IMAC3 (loaded with copper) and WCX2 ProteinChip® Arrays in duplicates. IMAC3 copper and WCX2 arrays (CIPHERGEN Biosystems Inc, Fremont, CA) were equilibrated two times with 150 $\mu$ L of binding buffer (100mM sodium phosphate + 0.5M NaCl pH7 for IMAC3, 100mM sodium acetate pH4 for WCX2). Each serum fraction was diluted in the corresponding binding buffer (1/5 dilution for IMAC3 and 1/10 dilution for WCX2) and 100 $\mu$ L was applied to each ProteinChip® array. Incubation was performed on a platform shaker at room temperature for 30 minutes. Each array was washed three times with 150 $\mu$ L of corresponding binding buffer and rinsed two times with water. ProteinChip® arrays were air-dried. Sinapinic acid matrix (prepared in 50% acetonitrile, 0.5% trifluoroacetic acid) was applied to each array. ProteinChip® arrays were read on a ProteinChip® PBSII Reader (CIPHERGEN Biosystems Inc.) A total of 253 laser shots were averaged for each array.

[0097] All publications and patent documents cited in this application are incorporated by reference in their entirety for all purposes to the same extent as if

each individual publication or patent document were so individually denoted. By their citation of various references in this document Applicants do not admit that any particular reference is "prior art" to their invention.

What we claim is:

1. A method for qualifying lung carcinoma status in a subject, comprised of analyzing a biological sample from said subject for a diagnostic level of a protein selected from either a first group consisting of

(i) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-157, IM-176, IM-445, IM-177, IM-440, IM-468, IM-438, IM-547, IM-359, IM-436, IM-106, IM-455, IM-444, IM-158, IM-265, IM-50, IM-159, IM-156, IM-439, IM-157, IM-508, IM-514, IM-478, IM-473, IM-360, IM-435, IM-150, IM-151, IM-110, IM-51, IM-163, IM-437, IM-546, IM-153, and IM-268,

or from a second group consisting of

(ii) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70, WM-369, WM-17, WM-473, WM-47, WM-203, WM-276, WM-279, WM-62, WM-366, WM-456, WM-428, WM-384, WM-287, WM-420, WM-292, WM-431, WM-455, WM-20, WM-340, WM-105, WM-389, WM-63, WM-354, WM-450, WM-466, WM-296, WM-343, WM-341, WM-339, WM-55, WM-66, WM-48, WM-38, WM-138, and WM-310,

wherein the biomarker is differentially present in samples of a subject with lung cancer and a normal subject that is free of lung cancer.

2. The method according to claim 1, wherein the protein is selected from either a first group consisting of

(i) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, and IM-155,

or from a second group consisting of

(ii) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, and WM-70.

3. The method according to claim 1, wherein the protein is selected from either a first group consisting of

(i) IM-522, IM-273, IM-520, IM-519, and IM-454,

or from a second group consisting



(ii) WM-61, WM-447, WM-446, WM-133, and WM-119.

4. The method according to claim 1, which uses a single biomarker selected from the group consisting of the WM-446 and WM-447.

5. A method for qualifying lung carcinoma risk in a subject, comprising  
(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from either a first group consisting of

(i) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-157, IM-176, IM-445, IM-177, IM-440, IM-468, IM-438, IM-547, IM-359, IM-436, IM-106, IM-455, IM-444, IM-158, IM-265, IM-50, IM-159, IM-156, IM-439, IM-157, IM-508, IM-514, IM-478, IM-473, IM-360, IM-435, IM-150, IM-151, IM-110, IM-51, IM-163, IM-437, IM-546, IM-153, and IM-268,

or from a second group consisting of

(ii) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70, WM-369, WM-17, WM-473, WM-47, WM-203, WM-276, WM-279, WM-62, WM-366, WM-456, WM-428, WM-384, WM-287, WM-420, WM-292, WM-431, WM-455, WM-20, WM-340, WM-105, WM-389, WM-63, WM-354, WM-450, WM-466, WM-296, WM-343, WM-341, WM-339, WM-55, WM-66, WM-48, WM-38, WM-138, and WM-310.

6. The method according to claim 5, wherein the pattern-recognition analysis is keyed to a pair of peaks selected either from a first group consisting of

(i) IM-266 and IM-474, IM-266 and IM-38, IM-266 and IM-454, IM-266 and IM-522, IM-266 and IM-544, IM-266 and IM-471, IM-474 and IM-151, IM-474 and IM-156, IM-474 and IM-544, IM-474 and IM-38, IM-522 and IM-507, IM-522 and IM-156, and IM-522 and IM-440;

or from a second group consisting of

(ii) WM-447 and WM-59, WM-447 and WM-19, WM-447 and WM-118, WM-447 and WM-473, WM-19 and WM-59, WM-19 and WM-473, WM-19 and WM-369, WM-61 and WM-154, WM-61 and WM-369, WM-118 and WM-59 and WM-282 and WM-127.

7. The method according to claim 5, wherein the pattern-recognition analysis is keyed to a pair of peaks selected from either a first group consisting of

(i) IM-266 and IM-474, IM-266 and IM-544, and IM-156 and IM-522;

or from a second group consisting of

(ii) WM-447 and WM-59, WM-447 and WM-19, and WM-19 and WM-59.

8. The method according to claim 5, wherein the pattern-recognition analysis is keyed to a triplet of peaks selected from

(i) IM-266, IM-454 and IM-474; and IM-266, IM-474 and IM-544;

or wherein the analysis is keyed to

(ii) WM-447, WM-19 and WM-473.

9. A kit for detecting and diagnosing lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers selected from either a first group consisting of

(i) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-157, IM-176, IM-445, IM-177, IM-440, IM-468, IM-438, IM-547, IM-359, IM-436, IM-106, IM-455, IM-444, IM-158, IM-265, IM-50, IM-159, IM-156, IM-439, IM-157, IM-508, IM-514, IM-478, IM-473, IM-360, IM-435, IM-150, IM-151, IM-110, IM-51, IM-163, IM-437, IM-546, IM-153, and IM-268,

or from a second group consisting of

(ii) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70, WM-369, WM-17, WM-473, WM-47, WM-203, WM-276, WM-279, WM-62, WM-366, WM-456, WM-428, WM-384, WM-287, WM-420, WM-292, WM-431, WM-455, WM-20, WM-340, WM-105, WM-389, WM-63, WM-354, WM-450, WM-466, WM-296, WM-343, WM-341, WM-339, WM-55, WM-66, WM-48, WM-38, WM-138, and WM-310, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

10. A kit according to claim 9, further comprising a washing solution or instructions for making a washing solution.

11. A kit according to claim 9, wherein the substrate is a SELDI probe that comprises either (i) functionalities that adsorb transition metal ions by chelation or (ii) functionalities that allow for cation exchange.

12. A method for qualifying lung adenocarcinoma status in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-447, WM-652, WM-61, WM-446, WM-290, WM-363, WM-133, WM-341, WM-285, WM-366, WM-282, WM-362, WM-310, WM-292, WM-120, WM-134, WM-276, WM-428, WM-277, WM-20, WM-119, WM-340, WM-48, WM-389, WM-450, WM-47, WM-343, WM-17, WM-583, WM-70, WM-706, WM-346, WM-466, WM-646, WM-384, WM-336, WM-294, WM-339, WM-473, WM-369, WM-38, WM-283, WM-685, WM-66, WM-55, WM-650, WM-307, WM-278, WM-342, and WM-429.

13. The method according to claim 12, wherein the protein is selected from the group consisting of WM-447, WM-652, WM-61, WM-446, WM-290, WM-363, WM-133, WM-341, WM-285, WM-366, WM-282, WM-362, WM-310, WM-292, and WM-120.

14. The method according to claim 12, wherein the protein is selected from the group consisting of WM-447, WM-652, WM-61, WM-446, WM-290.

15. A method for qualifying status of lung adenocarcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from either a first group consisting of WM-447, WM-652, WM-61, WM-446, WM-290, WM-363, WM-133, WM-341, WM-285, WM-366, WM-282, WM-362, WM-310, WM-292, WM-120, WM-134, WM-276, WM-428, WM-277, WM-20, WM-119, WM-340, WM-48, WM-

389, WM-450, WM-47, WM-343, WM-17, WM-583, WM-70, WM-706, WM-346, WM-466, WM-646, WM-384, WM-336, WM-294, WM-339, WM-473, WM-369, WM-38, WM-283, WM-685, WM-66, WM-55, WM-650, WM-307, WM-278, WM-342, and WM-429.

16. The method according to claim 15, wherein the protein is selected from the group consisting of WM-447, WM-652, WM-61, WM-446, WM-290, WM-363, WM-133, WM-341, WM-285, WM-366, WM-282, WM-362, WM-310, WM-292, and WM-120.

17. The method according to claim 15, wherein the protein is selected from the group consisting of WM-447, WM-652, WM-61, WM-446, WM-290.

18. A kit for detecting and diagnosing lung adenocarcinoma, comprising  
(A) an adsorbent attached to a substrate that retains one or more of biomarkers selected from the group consisting of WM-447, WM-652, WM-61, WM-446, WM-290, WM-363, WM-133, WM-341, WM-285, WM-366, WM-282, WM-362, WM-310, WM-292, WM-120, WM-134, WM-276, WM-428, WM-277, WM-20, WM-119, WM-340, WM-48, WM-389, WM-450, WM-47, WM-343, WM-17, WM-583, WM-70, WM-706, WM-346, WM-466, WM-646, WM-384, WM-336, WM-294, WM-339, WM-473, WM-369, WM-38, WM-283, WM-685, WM-66, WM-55, WM-650, WM-307, WM-278, WM-342, and WM-429, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

19. A kit according to claim 18, further comprising a washing solution or instructions for making a washing solution.

20. A kit according to claim 18, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

21. A method for qualifying squamous cell lung carcinoma status in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-447, WM-61, WM-277, WM-446, WM-133, WM-134, WM-363, WM-362, WM-276, WM-706, WM-203, WM-466, WM-366, WM-65, WM-70, WM-341, WM-429, WM-347, WM-17, WM-47, WM-431, WM-62, WM-473, WM-384, WM-438, WM-652, WM-282, WM-389, WM-290,

WM-278, WM-456, WM-673, WM-340, WM-55, WM-455, WM-645, WM-138, WM-420, WM-450, WM-369, WM-279, WM-342, WM-471, WM-674, WM-120, WM-20, WM-287, WM-83, WM-154, and WM-128.

22. The method according to claim 21, wherein the protein is selected from the group consisting of WM-447, WM-61, WM-277, WM-446, WM-133, WM-134, WM-363, WM-362, WM-276, WM-706, WM-203, WM-466, WM-366, WM-65, and WM-70.

23. The method according to claim 21, wherein the protein is selected from the group consisting of WM-447, WM-61, WM-277, WM-446, and WM-133.

24. A method for qualifying status of squamous cell lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from either a first group consisting of WM-447, WM-61, WM-277, WM-446, WM-133, WM-134, WM-363, WM-362, WM-276, WM-706, WM-203, WM-466, WM-366, WM-65, WM-70, WM-341, WM-429, WM-347, WM-17, WM-47, WM-431, WM-62, WM-473, WM-384, WM-438, WM-652, WM-282, WM-389, WM-290, WM-278, WM-456, WM-673, WM-340, WM-55, WM-455, WM-645, WM-138, WM-420, WM-450, WM-369, WM-279, WM-342, WM-471, WM-674, WM-120, WM-20, WM-287, WM-83, WM-154, and WM-128.

25. The method according to claim 24, wherein the protein is selected from the group consisting of WM-447, WM-61, WM-277, WM-446, WM-133, WM-134, WM-363, WM-362, WM-276, WM-706, WM-203, WM-466, WM-366, WM-65, and WM-70.

26. The method according to claim 24, wherein the protein is selected from the group consisting of WM-447, WM-61, WM-277, WM-446, and WM-133.

27. A kit for detecting and diagnosing squamous cell lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers selected from the group consisting of WM-447, WM-61, WM-277, WM-446, WM-133, WM-134, WM-363, WM-362, WM-276, WM-706, WM-203, WM-466, WM-366, WM-65, WM-70, WM-341, WM-429, WM-347, WM-17, WM-47, WM-431, WM-62, WM-473, WM-384, WM-438, WM-652, WM-282, WM-389, WM-290, WM-278, WM-456, WM-673, WM-340, WM-55, WM-455, WM-645, WM-138, WM-420, WM-450, WM-369, WM-279, WM-342, WM-471, WM-674, WM-120, WM-20, WM-287, WM-83, WM-154, and WM-128, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

28. A kit according to claim 27, further comprising a washing solution or instructions for making a washing solution.

29. A kit according to claim 27, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

30. A method for qualifying small cell lung carcinoma status in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-70, WM-706, WM-369, WM-447, WM-61, WM-652, WM-282, WM-446, WM-456, WM-134, WM-203, WM-646, WM-455, WM-65, WM-685, WM-473, WM-343, WM-466, WM-341, WM-340, WM-363, WM-339, WM-457, WM-86, WM-506, WM-72, WM-287, WM-82, WM-528, WM-85, WM-73, WM-138, WM-384, WM-83, WM-450, WM-310, WM-277, WM-79, WM-207, WM-278, WM-290, WM-366, WM-472, WM-420, WM-147, WM-55, WM-669, WM-357, WM-429, and WM-279.

31. The method according to claim 30, wherein the protein is selected from the group consisting of WM-70, WM-706, WM-369, WM-447, WM-61, WM-652, WM-282, WM-446, WM-456, WM-134, WM-203, WM-646, WM-455, WM-65, and WM-685.

32. The method according to claim 30, wherein the protein is selected from the group consisting of WM-70, WM-706, WM-369, WM-447, and WM-61.

33. A method for qualifying status of small cell lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from either a first group consisting of WM-70, WM-706, WM-369, WM-447, WM-61, WM-652, WM-282, WM-446, WM-456, WM-134, WM-203, WM-646, WM-455, WM-65, WM-685, WM-473, WM-343, WM-466, WM-341, WM-340, WM-363, WM-339, WM-457, WM-86, WM-506, WM-72, WM-287, WM-82, WM-528, WM-85, WM-73, WM-138, WM-384, WM-83, WM-450, WM-310, WM-277, WM-79, WM-207, WM-278, WM-290, WM-366, WM-472, WM-420, WM-147, WM-55, WM-669, WM-357, WM-429, and WM-279.

34. The method according to claim 33, wherein the protein is selected from the group consisting of WM-70, WM-706, WM-369, WM-447, WM-61, WM-652, WM-282, WM-446, WM-456, WM-134, WM-203, WM-646, WM-455, WM-65, and WM-685.

35. The method according to claim 33, wherein the protein is selected from the group consisting of WM-70, WM-706, WM-369, WM-447, and WM-61.

36. A kit for detecting and diagnosing small cell lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers selected from the group consisting of WM-70, WM-706, WM-369, WM-447, WM-61, WM-652, WM-282, WM-446, WM-456, WM-134, WM-203, WM-646, WM-455, WM-65, WM-685, WM-473, WM-343, WM-466, WM-341, WM-340, WM-363, WM-339, WM-457, WM-86, WM-506, WM-72, WM-287, WM-82, WM-528, WM-85, WM-73, WM-138, WM-384, WM-83, WM-450, WM-310, WM-277, WM-79, WM-207, WM-278, WM-290, WM-366, WM-472, WM-420, WM-147, WM-55, WM-669, WM-357, WM-429, and WM-279, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

37. A kit according to claim 36, further comprising a washing solution or instructions for making a washing solution.

38. A kit according to claim 36, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

39. A method for qualifying non-small cell lung carcinoma status in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-341, WM-342, WM-343, WM-48, WM-340, WM-346, WM-47, WM-339, WM-389, WM-669, WM-447, WM-652, WM-154, WM-587, WM-456, WM-450, WM-283, WM-207, WM-436, WM-384, WM-61, WM-167, WM-382, WM-285, WM-650, WM-203, WM-119, WM-282, WM-686, WM-383, WM-429, WM-11, WM-208, WM-451, WM-473, WM-220, WM-685, WM-338, WM-71, WM-266, WM-70, WM-545, WM-675, WM-446, WM-120, WM-267, WM-466, WM-347, WM-153, and WM-38.

40. The method according to claim 39, wherein the protein is selected from the group consisting of WM-341, WM-342, WM-343, WM-48, WM-340, WM-346, WM-47, WM-339, WM-389, WM-669, WM-447, WM-652, WM-154, WM-587, and WM-456.

41. The method according to claim 39, wherein the protein is selected from the group consisting of WM-341, WM-342, WM-343, WM-48, and WM-340.

42. A method for qualifying status of non-small cell lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from the group consisting of WM-341, WM-342, WM-343, WM-48, WM-340, WM-346, WM-47, WM-339, WM-389, WM-669, WM-447, WM-652, WM-154, WM-587, WM-456, WM-450, WM-283, WM-207, WM-436, WM-384, WM-61, WM-167, WM-382, WM-285, WM-650, WM-203, WM-119, WM-282, WM-686, WM-383, WM-429, WM-11, WM-208, WM-451, WM-473, WM-220, WM-685, WM-338, WM-71, WM-266, WM-70, WM-545, WM-675, WM-446, WM-120, WM-267, WM-466, WM-347, WM-153, and WM-38.



43. The method according to claim 42, wherein the protein is selected from the group consisting of WM-341, WM-342, WM-343, WM-48, WM-340, WM-346, WM-47, WM-339, WM-389, WM-669, WM-447, WM-652, WM-154, WM-587, and WM-456.

44. The method according to claim 42, wherein the protein is selected from the group consisting of WM-341, WM-342, WM-343, WM-48, and WM-340.

45. A kit for detecting and diagnosing non-small cell lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers WM-341, WM-342, WM-343, WM-48, WM-340, WM-346, WM-47, WM-339, WM-389, WM-669, WM-447, WM-652, WM-154, WM-587, WM-456, WM-450, WM-283, WM-207, WM-436, WM-384, WM-61, WM-167, WM-382, WM-285, WM-650, WM-203, WM-119, WM-282, WM-686, WM-383, WM-429, WM-11, WM-208, WM-451, WM-473, WM-220, WM-685, WM-338, WM-71, WM-266, WM-70, WM-545, WM-675, WM-446, WM-120, WM-267, WM-466, WM-347, WM-153, and WM-38, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

46. A kit according to claim 45, further comprising a washing solution or instructions for making a washing solution.

47. A kit according to claim 45, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

48. A method for qualifying large cell lung carcinoma status in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-16, WM-26, WM-499, WM-134, WM-647, WM-277, WM-310, WM-363, WM-446, WM-221, WM-648, WM-657, WM-290, WM-328, WM-447, WM-684, WM-183, WM-190, WM-686, WM-397, WM-466, WM-20, WM-17, WM-545, WM-47, WM-191, WM-147, WM-480, WM-590, WM-218, WM-285, WM-652, WM-651, WM-366, WM-403, WM-418, WM-430, WM-456, WM-714, WM-646, WM-109, WM-302, WM-587, WM-375, WM-131, WM-706, WM-398, WM-309, WM-55, and WM-488.

49. The method according to claim 48, wherein the protein is selected from the group consisting of WM-16, WM-26, WM-499, WM-134, WM-647, WM-277, WM-310, WM-363, WM-446, WM-221, WM-648, WM-657, WM-290, WM-328, and WM-447.

50. The method according to claim 48, wherein the protein is selected from the group consisting of WM-16, WM-26, WM-499, WM-134, and WM-647.

51. A method for qualifying status of large cell lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from the group consisting of WM-16, WM-26, WM-499, WM-134, WM-647, WM-277, WM-310, WM-363, WM-446, WM-221, WM-648, WM-657, WM-290, WM-328, WM-447, WM-684, WM-183, WM-190, WM-686, WM-397, WM-466, WM-20, WM-17, WM-545, WM-47, WM-191, WM-147, WM-480, WM-590, WM-218, WM-285, WM-652, WM-651, WM-366, WM-403, WM-418, WM-430, WM-456, WM-714, WM-646, WM-109, WM-302, WM-587, WM-375, WM-131, WM-706, WM-398, WM-309, WM-55, and WM-488.

52. The method according to claim 51, wherein the protein is selected from the group consisting of WM-16, WM-26, WM-499, WM-134, WM-647, WM-277, WM-310, WM-363, WM-446, WM-221, WM-648, WM-657, WM-290, WM-328, and WM-447.

53. The method according to claim 51, wherein the protein is selected from the group consisting of WM-16, WM-26, WM-499, WM-134, and WM-647.

54. A kit for detecting and diagnosing large cell lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers WM-16, WM-26, WM-499, WM-134, WM-647, WM-277, WM-310, WM-363, WM-446, WM-221, WM-648, WM-657, WM-290, WM-328, WM-447, WM-684, WM-183, WM-190, WM-686, WM-397, WM-466, WM-20, WM-17, WM-

545, WM-47, WM-191, WM-147, WM-480, WM-590, WM-218, WM-285, WM-652, WM-651, WM-366, WM-403, WM-418, WM-430, WM-456, WM-714, WM-646, WM-109, WM-302, WM-587, WM-375, WM-131, WM-706, WM-398, WM-309, WM-55, and WM-488, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

55. A kit according to claim 50, further comprising a washing solution or instructions for making a washing solution.

56. A kit according to claim 50, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

57. A method for distinguishing lung adenocarcinoma from squamous lung carcinoma in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-62, WM-415, WM-152, WM-385, WM-347, WM-134, WM-36, WM-108, WM-99, WM-151, WM-289, WM-363, WM-61, WM-117, WM-211, WM-362, WM-133, WM-414, WM-277, WM-141, WM-64, WM-135, WM-447, WM-383, WM-338, WM-63, WM-142, WM-446, WM-186, WM-111, WM-445, WM-455, WM-276, WM-444, WM-181, WM-35, WM-285, WM-456, WM-39, WM-82, WM-17, WM-203, WM-83, WM-412, WM-96, WM-74, WM-457, WM-431, WM-340, and WM-49.

58. The method according to claim 57, wherein the protein is selected from the group consisting of WM-62, WM-415, WM-152, WM-385, WM-347, WM-134, WM-36, WM-108, WM-99, WM-151, WM-289, WM-363, WM-61, WM-117, and WM-211.

59. The method according to claim 57, wherein the protein is selected from the group consisting of WM-62, WM-415, WM-152, WM-385, and WM-347.

60. A method for distinguishing lung adenocarcinoma from squamous lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from the group

consisting of WM-62, WM-415, WM-152, WM-385, WM-347, WM-134, WM-36, WM-108, WM-99, WM-151, WM-289, WM-363, WM-61, WM-117, WM-211, WM-362, WM-133, WM-414, WM-277, WM-141, WM-64, WM-135, WM-447, WM-383, WM-338, WM-63, WM-142, WM-446, WM-186, WM-111, WM-445, WM-455, WM-276, WM-444, WM-181, WM-35, WM-285, WM-456, WM-39, WM-82, WM-17, WM-203, WM-83, WM-412, WM-96, WM-74, WM-457, WM-431, WM-340, and WM-49.

61. The method according to claim 60, wherein the protein is selected from the group consisting of WM-62, WM-415, WM-152, WM-385, WM-347, WM-134, WM-36, WM-108, WM-99, WM-151, WM-289, WM-363, WM-61, WM-117, and WM-211.

62. The method according to claim 60, wherein the protein is selected from the group consisting of WM-62, WM-415, WM-152, WM-385, and WM-347.

63. A kit for distinguishing lung adenocarcinoma from squamous lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers WM-16, WM-26, WM-499, WM-134, WM-647, WM-277, WM-310, WM-363, WM-446, WM-221, WM-648, WM-657, WM-290, WM-328, WM-447, WM-684, WM-183, WM-190, WM-686, WM-397, WM-466, WM-20, WM-17, WM-545, WM-47, WM-191, WM-147, WM-480, WM-590, WM-218, WM-285, WM-652, WM-651, WM-366, WM-403, WM-418, WM-430, WM-456, WM-714, WM-646, WM-109, WM-302, WM-587, WM-375, WM-131, WM-706, WM-398, WM-309, WM-55, and WM-488, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

64. A kit according to claim 63, further comprising a washing solution or instructions for making a washing solution.

65. A kit according to claim 63, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

66. A method for distinguishing lung adenocarcinoma from small cell lung carcinoma in a subject, comprised of analyzing a biological sample from said subject

for a level of a protein selected from the group consisting of WM-457, WM-72, WM-369, WM-78, WM-79, WM-73, WM-64, WM-320, WM-419, WM-85, WM-82, WM-53, WM-412, WM-440, WM-455, WM-313, WM-456, WM-86, WM-70, WM-246, WM-360, WM-190, WM-418, WM-83, WM-257, WM-138, WM-47, WM-252, WM-282, WM-60, WM-68, WM-325, WM-402, WM-411, WM-405, WM-75, WM-417, WM-387, WM-26, WM-410, WM-420, WM-164, WM-67, WM-66, WM-391, WM-340, WM-428, WM-198, WM-312, and WM-152.

67. The method according to claim 66, wherein the protein is selected from the group consisting of WM-457, WM-72, WM-369, WM-78, WM-79, WM-73, WM-64, WM-320, WM-419, WM-85, WM-82, WM-53, WM-412, WM-440, and WM-455.

68. The method according to claim 66, wherein the protein is selected from the group consisting of WM-457, WM-72, WM-369, WM-78, and WM-79.

69. A method for distinguishing lung adenocarcinoma from small cell lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from either a first group consisting of WM-457, WM-72, WM-369, WM-78, WM-79, WM-73, WM-64, WM-320, WM-419, WM-85, WM-82, WM-53, WM-412, WM-440, WM-455, WM-313, WM-456, WM-86, WM-70, WM-246, WM-360, WM-190, WM-418, WM-83, WM-257, WM-138, WM-47, WM-252, WM-282, WM-60, WM-68, WM-325, WM-402, WM-411, WM-405, WM-75, WM-417, WM-387, WM-26, WM-410, WM-420, WM-164, WM-67, WM-66, WM-391, WM-340, WM-428, WM-198, WM-312, and WM-152.

70. The method according to claim 69, wherein the protein is selected from the group consisting of WM-457, WM-72, WM-369, WM-78, WM-79, WM-73, WM-64, WM-320, WM-419, WM-85, WM-82, WM-53, WM-412, WM-440, and WM-455.

71. The method according to claim 69, wherein the protein is selected from the group consisting of WM-457, WM-72, WM-369, WM-78, WM-79.

72. A kit for distinguishing lung adenocarcinoma from small cell lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers WM-276, WM-277, WM-362, WM-257, WM-363, WM-347, WM-53, WM-254, WM-17, WM-252, WM-431, WM-513, WM-446, WM-355, WM-447, WM-133, WM-245, WM-52, WM-96, WM-238, WM-243, WM-138, WM-62, WM-580, WM-134, WM-240, WM-256, WM-203, WM-111, WM-95, WM-247, WM-157, WM-242, WM-556, WM-63, WM-239, WM-234, WM-274, WM-370, WM-301, WM-449, WM-74, WM-261, WM-467, WM-237, WM-262, WM-295, WM-288, WM-384, and WM-37, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

73. A kit according to claim 72, further comprising a washing solution or instructions for making a washing solution.

74. A kit according to claim 72, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

75. A method for distinguishing squamous cell lung carcinoma from small cell lung carcinoma in a subject, comprised of analyzing a biological sample from said subject for a level of a protein selected from the group consisting of WM-276, WM-277, WM-362, WM-257, WM-363, WM-347, WM-53, WM-254, WM-17, WM-252, WM-431, WM-513, WM-446, WM-355, WM-447, WM-133, WM-245, WM-52, WM-96, WM-238, WM-243, WM-138, WM-62, WM-580, WM-134, WM-240, WM-256, WM-203, WM-111, WM-95, WM-247, WM-157, WM-242, WM-556, WM-63, WM-239, WM-234, WM-274, WM-370, WM-301, WM-449, WM-74, WM-261, WM-467, WM-237, WM-262, WM-295, WM-288, WM-384, and WM-37.

76. The method according to claim 75, wherein the protein is selected from the group consisting of WM-276, WM-277, WM-362, WM-257, WM-363, WM-347, WM-53, WM-254, WM-17, WM-252, WM-431, WM-513, WM-446, WM-355, and WM-447.

77. The method according to claim 75, wherein the protein is selected from the group consisting of WM-276, WM-277, WM-362, WM-257, and WM-363.

78. A method for distinguishing squamous cell lung carcinoma from small cell lung carcinoma in a subject, comprising

(A) providing a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, and

(B) extracting data from the spectrum and subjecting the data to pattern-recognition analysis that is keyed to at least one peak selected from either a first group consisting of WM-276, WM-277, WM-362, WM-257, WM-363, WM-347, WM-53, WM-254, WM-17, WM-252, WM-431, WM-513, WM-446, WM-355, WM-447, WM-133, WM-245, WM-52, WM-96, WM-238, WM-243, WM-138, WM-62, WM-580, WM-134, WM-240, WM-256, WM-203, WM-111, WM-95, WM-247, WM-157, WM-242, WM-556, WM-63, WM-239, WM-234, WM-274, WM-370, WM-301, WM-449, WM-74, WM-261, WM-467, WM-237, WM-262, WM-295, WM-288, WM-384, and WM-37.

79. The method according to claim 78, wherein the protein is selected from the group consisting of WM-276, WM-277, WM-362, WM-257, WM-363, WM-347, WM-53, WM-254, WM-17, WM-252, WM-431, WM-513, WM-446, WM-355, and WM-447.

80. The method according to claim 78, wherein the protein is selected from the group consisting of WM-276, WM-277, WM-362, WM-257, and WM-363.

81. A kit for distinguishing squamous cell lung carcinoma from small cell lung carcinoma, comprising

(A) an adsorbent attached to a substrate that retains one or more of the biomarkers WM-276, WM-277, WM-362, WM-257, WM-363, WM-347, WM-53, WM-254, WM-17, WM-252, WM-431, WM-513, WM-446, WM-355, WM-447, WM-133, WM-245, WM-52, WM-96, WM-238, WM-243, WM-138, WM-62, WM-580, WM-134, WM-240, WM-256, WM-203, WM-111, WM-95, WM-247, WM-157, WM-242, WM-556, WM-63, WM-239, WM-234, WM-274, WM-370, WM-301, WM-449, WM-74, WM-261, WM-467, WM-237, WM-262, WM-295, WM-288, WM-384, and WM-37, and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

82. A kit according to claim 81, further comprising a washing solution or instructions for making a washing solution.

83. A kit according to claim 81, wherein the substrate is a SELDI probe that comprises functionalities that allow for cation exchange.

84. Software for qualifying lung carcinoma status in a subject, comprising an algorithm for analyzing data extracted from a spectrum generated by mass spectroscopic analysis of a biological sample taken from the subject, wherein said data relates to one or more biomarkers selected from either a first group consisting of

(i) IM-522, IM-273, IM-520, IM-519, IM-454, IM-507, IM-521, IM-148, IM-266, IM-537, IM-471, IM-510, IM-544, IM-474, IM-155, IM-157, IM-176, IM-445, IM-177, IM-440, IM-468, IM-438, IM-547, IM-359, IM-436, IM-106, IM-455, IM-444, IM-158, IM-265, IM-50, IM-159, IM-156, IM-439, IM-157, IM-508, IM-514, IM-478, IM-473, IM-360, IM-435, IM-150, IM-151, IM-110, IM-51, IM-163, IM-437, IM-546, IM-153, and IM-268,

or from a second group consisting of

(ii) WM-61, WM-447, WM-446, WM-133, WM-119, WM-278, WM-134, WM-363, WM-282, WM-362, WM-120, WM-290, WM-65, WM-277, WM-70, WM-369, WM-17, WM-473, WM-47, WM-203, WM-276, WM-279, WM-62, WM-366, WM-456, WM-428, WM-384, WM-287, WM-420, WM-292, WM-431, WM-455, WM-20, WM-340, WM-105, WM-389, WM-63, WM-354, WM-450, WM-466, WM-296, WM-343, WM-341, WM-339, WM-55, WM-66, WM-48, WM-38, WM-138, and WM-310.

85. Software according to claim 84, wherein said algorithm carries out a pattern-recognition analysis that is keyed to data relating to at least one of the biomarkers.

86. Software according to claim 85, wherein said algorithm comprises classification tree analysis that is keyed to data relating to at least one of the biomarkers.



87. Software according to claim 85, wherein said algorithm comprises artificial neural network analysis that is keyed to data relating to at least one of the biomarkers.

88. A method for qualifying lung carcinoma status in a subject, comprised of analyzing a biological sample from said subject for a diagnostic level of a biomarker that is serum amyloid A protein or a fragment thereof.

89. A method according to claim 88, wherein said serum biomarker has an apparent molecular weight of about 2803, 3168, 3277, 3552, 3897, 4300, 4490, 4655, 5927, 6874, 7776, 7941, 8152, 8952, 9233, 10300, 10866, or 10851 Daltons.

90. A method according to claim 89, wherein said serum biomarker has an apparent molecular weight of about 3168, 3277, 3552, 3897, 4300, 4490, 4655, 7776, 7941, 8152, 8952, or 10851 Daltons.

91. A method according to claim 88, wherein said serum biomarker has an apparent molecular weight of about 11.5 to 11.7 kD.

92. A method according to claim 88, for qualifying risk of lung adenocarcinoma.

93. A method according to claim 88, for qualifying risk of squamous cell lung carcinoma.

94. A method according to claim 88, for qualifying risk of small cell lung carcinoma.

95. A method according to claim 88, for qualifying risk of non-small cell lung carcinoma.

96. A method according to claim 88, for qualifying risk of large cell lung carcinoma.

97. A kit for detecting and diagnosing lung carcinoma, comprising  
(A) an adsorbent attached to a substrate that retains one or more of the biomarkers that are serum amyloid A protein or a fragment thereof.

and

(B) instructions to detect the biomarker(s) by contacting a sample with the adsorbent and detecting the biomarker(s) retained by the adsorbent.

98. A kit according to claim 97, wherein said serum biomarker has an apparent molecular weight of about 2803, 3168, 3277, 3552, 3897, 4300, 4490, 4655, 5927, 6874, 7776, 7941, 8152, 8952, 9233, 10300, 10866, or 10851 Daltons.

99. A kit according to claim 98, wherein said serum biomarker has an apparent molecular weight of about 3168, 3277, 3552, 3897, 4300, 4490, 4655, 7776, 7941, 8152, 8952, or 10851 Daltons.

100. A kit according to claim 97, wherein said serum biomarker has an apparent molecular weight of about 11.5 to 11.7 kD.

101. A kit according to claim 97, further comprising a washing solution or instructions for making a washing solution.

102. A kit according to claim 97, wherein the substrate is a SELDI probe.

FIGURE 1A

MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION
IM-1	2011	A	IM-37	3893	A	IM-72	54026	A	IM-109	2882	B	IM-110	2967	B
IM-2	2030	A	IM-38	3960	A	IM-73	60170	A	IM-111	2977	B	IM-112	2994	B
IM-3	2069	A	IM-39	3972	A	IM-75	74372	A	IM-113	3031	B	IM-114	3048	B
IM-4	2128	A	IM-40	3984	A	IM-76	75545	A	IM-115	3148	B	IM-116	3166	B
IM-5	2146	A	IM-41	4066	A	IM-77	77543	A	IM-117	3283	B	IM-118	3308	B
IM-6	2186	A	IM-42	4178	A	IM-78	79507	A	IM-119	3332	B	IM-120	3432	B
IM-7	2232	A	IM-43	4287	A	IM-79	89854	A	IM-121	3450	B	IM-122	3561	B
IM-8	2277	A	IM-44	4297	A	IM-80	101831	A	IM-123	3615	B	IM-124	3714	B
IM-9	2295	A	IM-45	4309	A	IM-81	104301	A	IM-125	3730	B	IM-126	3834	B
IM-10	2318	A	IM-46	4484	A	IM-82	125160	A	IM-127	3899	B	IM-128	3969	B
IM-11	2411	A	IM-47	4649	A	IM-83	132976	A	IM-129	3986	B	IM-130	3997	B
IM-12	2434	A	IM-48	4798	A	IM-84	149099	A	IM-131	4013	B	IM-132	4181	B
IM-13	2467	A	IM-49	5104	A	IM-85	2016	B	IM-133	4297	B	IM-134	4311	B
IM-14	2482	A	IM-50	5918	A	IM-86	2029	B	IM-135	4465	B	IM-136	4484	B
IM-15	2498	A	IM-51	6122	A	IM-87	2144	B	IM-137	4579	B	IM-138	4608	B
IM-16	2565	A	IM-52	6192	A	IM-88	2130	B	IM-139	4669	B	IM-140	4747	B
IM-17	2574	A	IM-53	6452	A	IM-89	2168	B	IM-141	4862	B	IM-142	4891	B
IM-18	2586	A	IM-54	6660	A	IM-90	2184	B	IM-143	5033	B	IM-144	5077	B
IM-19	2605	A	IM-55	7766	A	IM-91	2200	B						
IM-20	2722	A	IM-56	8145	A	IM-92	2284	B						
IM-21	2746	A	IM-57	8954	A	IM-93	2299	B						
IM-22	2788	A	IM-58	9312	A	IM-94	2314	B						
IM-23	2855	A	IM-59	9449	A	IM-95	2414	B						
IM-24	2871	A	IM-60	10272	A	IM-96	2428	B						
IM-25	2984	A	IM-61	11663	A	IM-97	2451	B						
IM-26	3030	A	IM-62	13376	A	IM-98	2466	B						
IM-27	3144	A	IM-63	14698	A	IM-99	2483	B						
IM-28	3243	A	IM-64	15190	A	IM-100	2565	B						
IM-29	3273	A	IM-65	66758	A	IM-101	2583	B						
IM-30	3290	A	IM-66	15951	A	IM-102	2597	B						
IM-31	3369	A	IM-67	15172	A	IM-103	2697	B						
IM-32	3445	A	IM-68	15925	A	IM-104	2715	B						
IM-33	3483	A	IM-69	23436	A	IM-105	2740	B						
IM-34	3676	A	IM-70	39794	A	IM-106	2752	B						
IM-35	3779	A	IM-71	44166	A	IM-107	2767	B						
IM-36	3793	A		46890	A	IM-108	2865	B						

FIGURE 1B

MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION
IM-145	5099	B	IM-181	16018	B	IM-217	2130	C	IM-253	3733	C	IM-288	83848	C
IM-146	5143	B	IM-182	17346	B	IM-218	2145	C	IM-254	3833	C	IM-287	78892	C
IM-147	5158	B	IM-183	18311	B	IM-219	2167	C	IM-255	3900	C	IM-286	66294	C
IM-148	5272	B	IM-184	22586	B	IM-220	2182	C	IM-256	4010	C	IM-285	60882	C
IM-149	5306	B	IM-185	23422	B	IM-221	2199	C	IM-257	4145	C	IM-284	55898	C
IM-150	5349	B	IM-186	27969	B	IM-222	2211	C	IM-258	4187	C	IM-283	50625	C
IM-151	5364	B	IM-187	32783	B	IM-223	2230	C	IM-259	4299	C	IM-282	47307	C
IM-152	5421	B	IM-188	39808	B	IM-224	2250	C	IM-260	4466	C	IM-281	44460	C
IM-153	5554	B	IM-189	43110	B	IM-225	2280	C	IM-261	4582	C	IM-280	39770	C
IM-154	5711	B	IM-190	44454	B	IM-226	2297	C	IM-262	4813	C	IM-279	33296	C
IM-155	5876	B	IM-191	47215	B	IM-227	2317	C	IM-263	4876	C	IM-278	28001	C
IM-156	5916	B	IM-192	53784	B	IM-228	2412	C	IM-264	5032	C	IM-277	22321	C
IM-157	5931	B	IM-193	55952	B	IM-229	2428	C	IM-265	5365	C	IM-276	15939	C
IM-158	5988	B	IM-194	60573	B	IM-230	2468	C	IM-266	5932	C	IM-275	15114	C
IM-159	6137	B	IM-195	66346	B	IM-231	2481	C	IM-267	7767	C	IM-274	14041	C
IM-160	6200	B	IM-196	73387	B	IM-232	2498	C	IM-268	7973	C	IM-273	11705	C
IM-161	6443	B	IM-197	79203	B	IM-233	2567	C	IM-269	8143	C	IM-272	9293	C
IM-162	6644	B	IM-198	89302	B	IM-234	2585	C	IM-270	9187	C	IM-271	9187	C
IM-163	6958	B	IM-199	94226	B	IM-235	2599	C	IM-271	9187	C	IM-272	9293	C
IM-164	7481	B	IM-200	99358	B	IM-236	2698	C	IM-272	9293	C	IM-273	11705	C
IM-165	7568	B	IM-201	102096	B	IM-237	2715	C	IM-273	11705	C	IM-274	14041	C
IM-166	7765	B	IM-202	107199	B	IM-238	2745	C	IM-274	14041	C	IM-275	15114	C
IM-167	7955	B	IM-203	116936	B	IM-239	2766	C	IM-275	15114	C	IM-276	15939	C
IM-168	8144	B	IM-204	119487	B	IM-240	2867	C	IM-276	15939	C	IM-277	22321	C
IM-169	8698	B	IM-205	122103	B	IM-241	2885	C	IM-277	22321	C	IM-278	28001	C
IM-170	8821	B	IM-206	125431	B	IM-242	2998	C	IM-278	28001	C	IM-279	33296	C
IM-171	8944	B	IM-207	132052	B	IM-243	3052	C	IM-279	33296	C	IM-280	39770	C
IM-172	9138	B	IM-208	138518	B	IM-244	3096	C	IM-280	39770	C	IM-281	44460	C
IM-173	9298	B	IM-209	145147	B	IM-245	3151	C	IM-281	44460	C	IM-282	47307	C
IM-174	9390	B	IM-210	157502	B	IM-246	3167	C	IM-282	47307	C	IM-283	50625	C
IM-175	9516	B	IM-211	168579	B	IM-247	3286	C	IM-283	50625	C	IM-284	55898	C
IM-176	11711	B	IM-212	173391	B	IM-248	3303	C	IM-284	55898	C	IM-285	60882	C
IM-177	11914	B	IM-213	2011	C	IM-249	3335	C	IM-285	60882	C	IM-286	66294	C
IM-178	14033	B	IM-214	2030	C	IM-250	3448	C	IM-286	66294	C	IM-287	78892	C
IM-179	15110	B	IM-215	2050	C	IM-251	3619	C	IM-287	78892	C	IM-288	83848	C
IM-180	15838	B	IM-216	2096	C	IM-252	3709	C	IM-288	83848	C			

MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION
IM-289	89081	C	IM-325	2565	D	IM-361	13857	D	IM-397	2082	E	IM-401	2187	E
IM-290	94147	C	IM-326	2582	D	IM-362	14056	D	IM-398	2128	E	IM-402	2206	E
IM-291	99324	C	IM-327	2597	D	IM-363	15108	D	IM-399	2148	E	IM-403	2232	E
IM-292	107163	C	IM-328	2716	D	IM-364	15844	D	IM-400	2170	E	IM-404	2250	E
IM-293	110350	C	IM-329	2747	D	IM-365	22243	D	IM-401	2187	E	IM-405	2279	E
IM-294	113339	C	IM-330	2767	D	IM-366	25465	D	IM-402	2206	E	IM-406	2296	E
IM-295	116291	C	IM-331	2866	D	IM-367	28022	D	IM-403	2232	E	IM-407	2314	E
IM-296	122769	C	IM-332	2882	D	IM-368	33272	D	IM-404	2250	E	IM-408	2354	E
IM-297	131908	C	IM-333	2994	D	IM-369	40149	D	IM-405	2279	E	IM-409	2394	E
IM-298	145248	C	IM-334	3032	D	IM-370	43113	D	IM-406	2296	E	IM-410	2413	E
IM-299	159252	C	IM-335	3050	D	IM-371	44219	D	IM-407	2314	E	IM-411	2436	E
IM-300	165164	C	IM-336	3148	D	IM-372	47196	D	IM-408	2354	E	IM-412	2457	E
IM-301	174928	C	IM-337	3164	D	IM-373	51062	D	IM-409	2394	E	IM-413	2466	E
IM-302	196003	C	IM-338	3278	D	IM-374	56082	D	IM-410	2413	E	IM-414	2499	E
IM-303	2007	D	IM-339	3334	D	IM-375	58239	D	IM-411	2436	E	IM-415	2566	E
IM-304	2016	D	IM-340	3385	D	IM-376	60285	D	IM-412	2457	E	IM-416	2583	E
IM-305	2030	D	IM-341	3432	D	IM-377	66148	D	IM-413	2466	E	IM-417	2612	E
IM-306	2052	D	IM-342	3451	D	IM-378	73688	D	IM-414	2499	E	IM-418	2662	E
IM-307	2099	D	IM-343	3617	D	IM-379	77433	D	IM-415	2566	E	IM-419	2723	E
IM-308	2130	D	IM-344	3701	D	IM-380	79986	D	IM-416	2583	E	IM-420	2738	E
IM-309	2144	D	IM-345	3725	D	IM-381	80844	D	IM-417	2612	E	IM-421	2750	E
IM-310	2154	D	IM-346	3833	D	IM-382	88962	D	IM-418	2662	E	IM-422	2849	E
IM-311	2166	D	IM-347	3899	D	IM-383	94399	D	IM-419	2723	E	IM-423	2867	E
IM-312	2184	D	IM-348	4008	D	IM-384	99419	D	IM-420	2738	E	IM-424	3036	E
IM-313	2204	D	IM-349	4157	D	IM-385	108395	D	IM-421	2750	E	IM-425	3147	E
IM-314	2231	D	IM-350	4297	D	IM-386	116433	D	IM-422	2849	E	IM-426	3281	E
IM-315	2252	D	IM-351	4580	D	IM-387	123337	D	IM-423	2867	E	IM-427	3319	E
IM-316	2275	D	IM-352	4805	D	IM-388	131977	D	IM-424	3036	E	IM-428	3445	E
IM-317	2299	D	IM-353	6946	D	IM-389	145658	D	IM-425	3147	E	IM-429	3693	E
IM-318	2316	D	IM-354	7053	D	IM-390	152603	D	IM-426	3281	E	IM-430	3731	E
IM-319	2412	D	IM-355	7767	D	IM-391	159524	D	IM-427	3319	E	IM-431	3818	E
IM-320	2435	D	IM-356	7954	D	IM-392	196072	D	IM-428	3445	E	IM-432	3885	E
IM-321	2466	D	IM-357	8139	D	IM-393	2010	E	IM-429	3693	E			
IM-322	2480	D	IM-358	9292	D	IM-394	2029	E	IM-430	3731	E			
IM-323	2499	D	IM-359	11671	D	IM-395	2050	E	IM-431	3818	E			
IM-324	2518	D	IM-360	13727	D	IM-396	2068	E	IM-432	3885	E			

**FIGURE 1D**

MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION	MARKER ID	MW	FRACTION
IM-433	4136	E	IM-469	86211	E	IM-505	4174	F	IM-541	95033	F	IM-541	95033	F	IM-541	95033	F
IM-434	4169	E	IM-470	89407	E	IM-506	4362	F	IM-542	100310	F	IM-542	100310	F	IM-542	100310	F
IM-435	4257	E	IM-471	100270	E	IM-507	4473	F	IM-543	116889	F	IM-543	116889	F	IM-543	116889	F
IM-436	4277	E	IM-472	109638	E	IM-508	4631	F	IM-544	132711	F	IM-544	132711	F	IM-544	132711	F
IM-437	4355	E	IM-473	117132	E	IM-509	4822	F	IM-545	147276	F	IM-545	147276	F	IM-545	147276	F
IM-438	4369	E	IM-474	132843	E	IM-510	5862	F	IM-546	160768	F	IM-546	160768	F	IM-546	160768	F
IM-439	4470	E	IM-475	147160	E	IM-511	6192	F									
IM-440	4486	E	IM-476	152199	E	IM-512	6941	F									
IM-441	4541	E	IM-477	166461	E	IM-513	7626	F									
IM-442	4634	E	IM-478	176635	E	IM-514	7772	F									
IM-443	4841	E	IM-479	2011	F	IM-515	7957	F									
IM-444	5862	E	IM-480	2030	F	IM-516	8150	F									
IM-445	5911	E	IM-481	2128	F	IM-517	8954	F									
IM-446	6649	E	IM-482	2149	F	IM-518	9300	F									
IM-447	6952	E	IM-483	2186	F	IM-519	11545	F									
IM-448	7769	E	IM-484	2207	F	IM-520	11717	F									
IM-449	8148	E	IM-485	2279	F	IM-521	13887	F									
IM-450	8260	E	IM-486	2299	F	IM-522	14073	F									
IM-451	8785	E	IM-487	2319	F	IM-523	15196	F									
IM-452	9301	E	IM-488	2412	F	IM-524	15903	F									
IM-453	10071	E	IM-489	2434	F	IM-525	22460	F									
IM-454	11721	E	IM-490	2467	F	IM-526	23135	F									
IM-455	13910	E	IM-491	2485	F	IM-527	28135	F									
IM-456	15919	E	IM-492	2582	F	IM-528	33577	F									
IM-457	22422	E	IM-493	2605	F	IM-529	39813	F									
IM-458	28233	E	IM-494	2697	F	IM-530	42344	F									
IM-459	33490	E	IM-495	2751	F	IM-531	43274	F									
IM-460	43121	E	IM-496	2865	F	IM-532	44345	F									
IM-461	44558	E	IM-497	3036	F	IM-533	51007	F									
IM-462	4694	E	IM-498	3151	F	IM-534	56318	F									
IM-463	50954	E	IM-499	3372	F	IM-535	60079	F									
IM-464	54478	E	IM-500	3440	F	IM-536	66690	F									
IM-465	60041	E	IM-501	3488	F	IM-537	75122	F									
IM-466	66652	E	IM-502	3717	F	IM-538	78429	F									
IM-467	75580	E	IM-503	3890	F	IM-539	81249	F									
IM-468	79463	E	IM-504	4155	F	IM-540	89384	F									

FIGURE 2

RANK	MW	MARKER ID	RANK	MW	MARKER ID
1	14073	IM-522	39	117132	IM-473
2	11705	IM-273	40	13727	IM-360
3	11717	IM-520	41	4257	IM-435
4	11545	IM-519	42	5349	IM-150
5	11721	IM-454	43	5364	IM-151
6	4473	IM-507	44	2967	IM-110
7	13887	IM-521	45	6122	IM-51
8	5272	IM-148	46	6958	IM-163
9	5365	IM-266	47	4355	IM-437
10	75122	IM-537	48	160768	IM-546
11	100270	IM-471	49	5554	IM-153
12	5862	IM-510	50	7767	IM-268
13	132711	IM-544			
14	132843	IM-474			
15	5876	IM-155			
16	5932	IM-157			
17	11711	IM-176			
18	5911	IM-445			
19	11914	IM-177			
20	4486	IM-440			
21	79463	IM-468			
22	4369	IM-438			
23	100310	IM-542			
24	11671	IM-359			
25	4277	IM-436			
26	2752	IM-106			
27	13910	IM-455			
28	5862	IM-444			
29	5988	IM-158			
30	5347	IM-265			
31	5918	IM-50			
32	6137	IM-159			
33	5916	IM-156			
34	4470	IM-439			
35	5931	IM-157			
36	4631	IM-508			
37	7772	IM-514			
38	176635	IM-478			

Figure 3A

[illegible]



Figure 38

WM-51	A	6827	7561	7770	7576
WM-52	A	7561	7561	7770	7576
WM-53	A	7768	7768	7768	7770
WM-54	A	8139	8139	8139	8143
WM-55	A	8611	8611	8615	8615
WM-56	A	8941	8941	8945	8946
WM-57	A	9142	9142	9144	9142
WM-58	A	9303	9303	9301	9301
WM-59	A	10074	10074	10078	10269
WM-60	A	10279	10279	10278	10269
WM-61	A	10657	10657	10652	10659
WM-62	A	11460	11460	11460	11492
WM-63	A	11699	11699	11697	11691
WM-64	A	12677	12677	12680	12686
WM-65	A	13355	13355	13350	13358
WM-66	A	14701	14701	14702	14702
WM-67	A	15147	15147	15134	15142
WM-68	A	15881	15881	15880	15889
WM-69	A	23472	23472	23462	23472
WM-70	A	26102	26102	26100	26129
WM-71	A	28767	28767	28764	28768
WM-72	A	33328	33328	33341	33399
WM-73	A	39871	39871	39851	39811
WM-74	A	44327	44327	44363	44280
WM-75	A	50167	50167	50098	50037
WM-76	A	51314	51314	51279	51175
WM-77	A	53760	53760	53803	53824
WM-78	A	60059	60059	60103	60398
WM-79	A	62287	62287	62287	62237
WM-80	A	63019	63019	63049	63102
WM-81	A	66655	66655	66660	66571
WM-82	A	68432	68432	68577	68533
WM-83	A	74529	74529	74641	74611
WM-84	A	78052	78052	78059	78078
WM-85	A	88300	88300	88372	89224
WM-86	A	100865	100865	100872	89554
WM-87	A	112640	112640	100884	101195
WM-88	A	124396	124396	112674	112012
WM-89	A	134770	134770	124241	124350
WM-90	A	147895	147895	134748	134581
WM-91	A	171308	171308	147603	147686
WM-92	A	2011	2011	171410	148177
WM-93	B	2087	2087	2012	171699
WM-94	B	2267	2267	2079	2012
WM-95	B	2278	2278	2079	2068
WM-96	B	2276	2276	2228	2081
WM-97	B	2283	2283	2278	2228
WM-98	B	2428	2428	2285	2278
WM-99	B	2501	2501	2294	2285
WM-100	B	2590	2590	2428	2428
WM-101	B	2590	2590	2488	2488

Figure 3c

[illegible]







Figure 3G

WMA-306	D	164336	160431	153873	153876	160722	160347	160310
WMA-307	D	160515	165737	160442	160441	160722	160347	160310
WMA-308	D	165704	165737	165848	165847	165853	165853	165853
WMA-309	D	175390	175478	174282	175447	175360	175446	175360
WMA-310	D	182482	182482	182482	182482	182482	182482	182482
WMA-311	E	2010	2010	2010	2010	2010	2010	2010
WMA-312	E	2029	2029	2029	2029	2029	2029	2029
WMA-313	E	2062	2062	2062	2062	2062	2062	2062
WMA-314	E	2067	2067	2067	2067	2067	2067	2067
WMA-315	E	2081	2081	2081	2081	2081	2081	2081
WMA-316	E	2128	2128	2128	2128	2128	2128	2128
WMA-317	E	2168	2167	2168	2165	2168	2167	2165
WMA-318	E	2185	2186	2185	2185	2186	2181	2181
WMA-319	E	2212	2212	2211	2215	2212	2212	2212
WMA-320	E	2231	2231	2231	2231	2232	2231	2231
WMA-321	E	2278	2278	2278	2278	2278	2278	2278
WMA-322	E	2285	2285	2285	2285	2285	2285	2285
WMA-323	E	2311	2311	2311	2311	2311	2311	2311
WMA-324	E	2389	2389	2389	2389	2389	2389	2389
WMA-325	E	2481	2474	2481	2481	2481	2481	2481
WMA-326	E	2489	2489	2489	2489	2489	2489	2489
WMA-327	E	2567	2567	2567	2567	2567	2567	2567
WMA-328	E	2581	2581	2581	2581	2581	2581	2581
WMA-329	E	2737	2737	2737	2737	2737	2737	2737
WMA-330	E	2862	2862	2862	2862	2862	2862	2862
WMA-331	E	3167	3167	3167	3167	3167	3167	3167
WMA-332	E	3321	3319	3318	3318	3315	3318	3318
WMA-333	E	3443	3442	3439	3439	3439	3439	3439
WMA-334	E	3818	3817	3818	3818	3818	3817	3818
WMA-335	E	3891	3891	3891	3891	3891	3891	3891
WMA-336	E	4084	4085	4085	4085	4085	4085	4085
WMA-337	E	4146	4146	4146	4146	4146	4146	4146
WMA-338	E	4211	4211	4210	4210	4212	4208	4208
WMA-339	E	4259	4259	4259	4259	4257	4258	4258
WMA-340	E	4355	4355	4355	4355	4357	4353	4353
WMA-341	E	4474	4469	4469	4472	4469	4468	4468
WMA-342	E	4541	4541	4540	4543	4540	4537	4537
WMA-343	E	4628	4628	4628	4631	4628	4625	4625
WMA-344	E	4829	4829	4829	4829	4829	4829	4829
WMA-345	E	4859	4859	4859	4859	4859	4854	4854
WMA-346	E	5073	5070	5069	5069	5069	5066	5066
WMA-347	E	5946	5946	5946	5946	5946	5940	5940
WMA-348	E	6462	6462	6461	6461	6461	6461	6461
WMA-349	E	6643	6643	6644	6642	6644	6643	6643
WMA-350	E	6879	6879	6880	6880	6880	6878	6878
WMA-351	E	6885	6885	6885	6885	6885	6884	6884
WMA-352	E	6945	6945	6945	6945	6945	6945	6945
WMA-353	E	7193	7193	7193	7193	7193	7193	7193
WMA-354	E	7770	7770	7769	7776	7772	7769	7769
WMA-355	E	8145	8145	8145	8147	8147	8147	8147
WMA-356	E	8701	8701	8702	8704	8700	8701	8701

Figure 3H

WNA-357 E	8928	8934	8928	8927	8928	8928	8931	8930
WNA-358 E	8949	8955	8955	8952	8955	8955	8945	8949
WNA-359 E	9176	9181	9176	9176	9176	9176	9177	9194
WNA-360 E	9207	9207	9209	9207	9209	9209	9201	9208
WNA-361 E	9470	9478	9454	9470	9454	9454	9470	9477
WNA-362 E	11532	11515	11520	11520	11520	11520	11533	11627
WNA-363 E	11714	11705	11699	11687	11699	11699	11714	11716
WNA-364 E	12482	12450	12452	12452	12452	12452	12481	12459
WNA-365 E	12814	12804	12813	12822	12813	12813	12817	12818
WNA-366 E	13897	13876	13888	13896	13888	13888	13901	13907
WNA-367 E	15148	15160	15154	15154	15154	15154	15155	15154
WNA-368 E	15991	15991	15991	15999	15991	15991	15999	15998
WNA-369 E	17354	17370	17320	17381	17302	17382	17363	17365
WNA-370 E	22359	22345	22348	22322	22348	22314	22364	22332
WNA-371 E	28169	28161	28162	28165	28161	28161	28154	28172
WNA-372 E	33437	33451	33431	33426	33431	33391	33425	33413
WNA-373 E	39455	39235	39235	44591	44578	44578	44551	39238
WNA-374 E	44603	44598	44627	44591	44578	44578	44551	44575
WNA-375 E	51369	51376	51357	51429	51308	51308	51400	51409
WNA-376 E	60921	60942	60826	60847	60744	60744	60747	59799
WNA-377 E	66557	66538	66572	66568	66527	66561	66561	66564
WNA-378 E	75353	75469	75368	75201	75168	75168	75459	75413
WNA-379 E	79373	79334	79338	79336	79767	79767	79378	76327
WNA-380 E	89514	89286	89480	89370	89337	89421	89542	89391
WNA-381 E	100004	89891	89878	89852	100091	100149	100026	100001
WNA-382 E	110157	109954	109994	109816	109772	110474	109921	109921
WNA-383 E	117446	117603	117746	117739	117918	117969	117455	117455
WNA-384 E	132840	132621	132605	132627	132676	132657	132642	132642
WNA-385 E	146553	146507	146553	146206	146470	145912	146010	146853
WNA-386 E	165171	165233	165192	165727	165654	165549	153495	152534
WNA-387 E	182678	182591	182289	182269	182217	182093	165920	165575
WNA-388 E	197393	197348	197208	197254	197265	183083	182701	183356
WNA-389 E	2010	2010	2010	2012	2011	2011	197249	197301
WNA-390 F	2028	2028	2028	2030	2029	2029	2029	2029
WNA-391 F	2033	2033	2033	2033	2034	2034	2032	2032
WNA-392 F	2065	2068	2065	2065	2072	2072	2052	2052
WNA-393 F	2068	2068	2068	2068	2071	2071	2070	2067
WNA-394 F	2065	2065	2063	2064	2064	2064	2064	2060
WNA-395 F	2126	2126	2126	2126	2126	2126	2126	2126
WNA-396 F	2164	2167	2167	2167	2165	2165	2166	2166
WNA-397 F	2214	2213	2214	2214	2212	2213	2213	2213
WNA-398 F	2234	2234	2233	2234	2233	2233	2234	2234
WNA-399 F	2292	2297	2297	2297	2278	2278	2279	2279
WNA-400 F	2290	2296	2297	2297	2297	2297	2296	2296
WNA-401 F	2380	2386	2380	2389	2389	2389	2389	2389
WNA-402 F	2413	2413	2413	2412	2413	2413	2412	2412
WNA-403 F	2462	2462	2476	2474	2482	2482	2462	2462
WNA-404 F	2501	2500	2500	2501	2500	2501	2501	2500
WNA-405 F	2568	2569	2570	2569	2568	2568	2568	2568
WNA-406 F	2581	2582	2581	2581	2582	2581	2581	2582
WNA-407 F	2597	2597	2700	2595	2609	2609	2595	2592

Figure 31

[illegible]



Figure 3J

WM-458 F	33504	33487	33512	33520
WM-460 F	43434	43457	43488	43415
WM-461 F	44066	44078	44094	44082
WM-462 F	51132	51189	51226	51052
WM-463 F	59467	59500	59510	59517
WM-464 F	59469	59728	59452	59470
WM-465 F	59575	59728	59587	59546
WM-466 F	68644	68632	68674	68546
WM-467 F	75418	75313	75455	75438
WM-468 F	75421	75368	75464	75404
WM-469 F	83331	83304	83307	83132
WM-470 F	89519	89546	89576	89588
WM-471 F	89779	89751	89740	89688
WM-472 F	100418	100426	100524	100523
WM-473 F	117035	117130	117127	116803
WM-474 F	132765	132840	132729	132707
WM-475 F	148470	148960	148862	148181
WM-476 A	160334	161492	160783	160570
WM-477 A		2184	2126	
WM-478 A		2414	2139	
WM-479 A		2433	2413	2433
WM-480 A		2568	2434	
WM-481 A			2451	
WM-482 A	2567		2468	
WM-483 A			2598	
WM-484 A			2598	
WM-485 A			2626	
WM-486 A			2637	
WM-487 A			2662	
WM-488 A			2680	
WM-489 A			2684	
WM-490 A			2684	
WM-491 A			2679	
WM-492 A	2881		2684	
WM-493 A			3015	
WM-494 A			3044	
WM-495 A			3061	
WM-496 A			3388	
WM-497 A			3469	
WM-498 A			3480	
WM-499 A			3611	
WM-500 A			3720	
WM-501 A			4161	
WM-502 A			4191	
WM-503 A			4209	
WM-504 A			4525	
WM-505 A			4525	
WM-506 A			4525	
WM-507 A			4525	
WM-508 A			4525	
WM-509 A			4728	
WM-510 A				
WM-511 A				
WM-512 A				
WM-513 A				
WM-514 A				
WM-515 A				
WM-516 A				
WM-517 A				
WM-518 A				
WM-519 A				
WM-520 A				
WM-521 A				
WM-522 A				
WM-523 A				
WM-524 A				
WM-525 A				
WM-526 A				
WM-527 A				
WM-528 A				
WM-529 A				
WM-530 A				
WM-531 A				
WM-532 A				
WM-533 A				
WM-534 A				
WM-535 A				
WM-536 A				
WM-537 A				
WM-538 A				
WM-539 A				
WM-540 A				
WM-541 A				
WM-542 A				
WM-543 A				
WM-544 A				
WM-545 A				
WM-546 A				
WM-547 A				
WM-548 A				
WM-549 A				
WM-550 A				
WM-551 A				
WM-552 A				
WM-553 A				
WM-554 A				
WM-555 A				
WM-556 A				
WM-557 A				
WM-558 A				
WM-559 A				
WM-560 A				
WM-561 A				
WM-562 A				
WM-563 A				
WM-564 A				
WM-565 A				
WM-566 A				
WM-567 A				
WM-568 A				
WM-569 A				
WM-570 A				
WM-571 A				
WM-572 A				
WM-573 A				
WM-574 A				
WM-575 A				
WM-576 A				
WM-577 A				
WM-578 A				
WM-579 A				
WM-580 A				
WM-581 A				
WM-582 A				
WM-583 A				
WM-584 A				
WM-585 A				
WM-586 A				
WM-587 A				
WM-588 A				
WM-589 A				
WM-590 A				
WM-591 A				
WM-592 A				
WM-593 A				
WM-594 A				
WM-595 A				
WM-596 A				
WM-597 A				
WM-598 A				
WM-599 A				
WM-600 A				
WM-601 A				
WM-602 A				
WM-603 A				
WM-604 A				
WM-605 A				
WM-606 A				
WM-607 A				
WM-608 A				
WM-609 A				
WM-610 A				
WM-611 A				
WM-612 A				
WM-613 A				
WM-614 A				
WM-615 A				
WM-616 A				
WM-617 A				
WM-618 A				
WM-619 A				
WM-620 A				
WM-621 A				
WM-622 A				
WM-623 A				
WM-624 A				
WM-625 A				
WM-626 A				
WM-627 A				
WM-628 A				
WM-629 A				
WM-630 A				
WM-631 A				
WM-632 A				
WM-633 A				
WM-634 A				
WM-635 A				
WM-636 A				
WM-637 A				
WM-638 A				
WM-639 A				
WM-640 A				
WM-641 A				
WM-642 A				
WM-643 A				
WM-644 A				
WM-645 A				
WM-646 A				
WM-647 A				
WM-648 A				
WM-649 A				
WM-650 A				
WM-651 A				
WM-652 A				
WM-653 A				
WM-654 A				
WM-655 A				
WM-656 A				
WM-657 A				
WM-658 A				
WM-659 A				
WM-660 A				
WM-661 A				
WM-662 A				
WM-663 A				
WM-664 A				
WM-665 A				
WM-666 A				
WM-667 A				
WM-668 A				
WM-669 A				
WM-670 A				
WM-671 A				
WM-672 A				
WM-673 A				
WM-674 A				
WM-675 A				
WM-676 A				
WM-677 A				
WM-678 A				
WM-679 A				
WM-680 A				
WM-681 A				
WM-682 A				
WM-683 A				
WM-684 A				
WM-685 A				
WM-686 A				
WM-687 A				
WM-688 A				
WM-689 A				
WM-690 A				
WM-691 A				
WM-692 A				
WM-693 A				
WM-694 A				
WM-695 A				
WM-696 A				
WM-697 A				
WM-698 A				
WM-699 A				
WM-700 A				
WM-701 A				
WM-702 A				
WM-703 A				
WM-704 A				
WM-705 A				
WM-706 A				
WM-707 A				
WM-708 A				
WM-709 A				
WM-710 A				
WM-711 A				
WM-712 A				
WM-713 A				
WM-714 A				
WM-715 A				
WM-716 A				
WM-717 A				
WM-718 A				
WM-719 A				
WM-720 A				
WM-721 A				
WM-722 A				
WM-723 A				
WM-724 A				
WM-725 A				
WM-726 A				
WM-727 A				
WM-728 A				
WM-729 A				
WM-730 A				
WM-731 A				
WM-732 A				
WM-733 A				
WM-734 A				
WM-735 A				
WM-736 A				
WM-737 A				
WM-738 A				
WM-739 A				
WM-740 A				
WM-741 A				
WM-742 A				
WM-743 A				
WM-744 A				
WM-745 A				
WM-746 A				
WM-747 A				
WM-748 A				
WM-749 A				
WM-750 A				
WM-751 A				
WM-752 A				
WM-753 A				
WM-754 A				
WM-755 A				
WM-756 A				
WM-757 A				
WM-758 A				
WM-759 A				
WM-760 A				
WM-761 A				
WM-762 A				
WM-763 A				
WM-764 A				
WM-765 A				
WM-766 A				
WM-767 A				
WM-768 A				
WM-769 A				
WM-770 A				
WM-771 A				
WM-772 A				
WM-773 A				
WM-774 A				
WM-775 A				
WM-776 A				
WM-777 A				
WM-778 A				
WM-779 A				
WM-780 A				
WM-781 A				
WM-782 A				
WM-783 A				
WM-784 A				
WM-785 A				
WM-786 A				
WM-787 A				
WM-788 A				
WM-789 A				
WM-790 A				
WM-791 A				
WM-792 A				
WM-793 A				
WM-794 A				
WM-795 A				
WM-796 A				
WM-797 A				
WM-798 A				
WM-799 A				
WM-800 A				
WM-801 A				
WM-802 A				
WM-803 A				
WM-804 A				
WM-805 A				
WM-806 A				
WM-807 A				
WM-808 A				
WM-809 A				
WM-810 A				
WM-811 A				
WM-812 A				
WM-813 A				
WM-814 A				
WM-815 A				
WM-816 A				
WM-817 A				
WM-818 A				
WM-819 A				
WM-820 A				
WM-821 A				
WM-822 A				
WM-823 A				
WM-824 A				
WM-825 A				
WM-826 A				
WM-827 A				
WM-828 A				
WM-829 A				
WM-830 A				
WM-831 A		</		







Figure 3 N

WM-663	E	2599	2598	2488	2487	2723	2723	2723	2723
WM-664	E	2730	2730	2730	2730	2751	2751	2751	2751
WM-665	E	2751	2751	2751	2751	3085	3085	3085	3085
WM-666	E		3375	3375	3375	4128	4128	4128	4128
WM-667	E		4131	4131	4131				
WM-668	E		4162	4162	4162				
WM-669	E								
WM-670	E								
WM-671	E								
WM-672	E								
WM-673	E								
WM-674	E								
WM-675	E								
WM-676	E								
WM-677	E								
WM-678	E								
WM-679	E								
WM-680	E								
WM-681	E								
WM-682	E								
WM-683	E								
WM-684	E								
WM-685	E								
WM-686	E								
WM-687	E								
WM-688	E								
WM-689	E								
WM-690	E								
WM-691	E								
WM-692	E								
WM-693	E								
WM-694	F								
WM-695	F								
WM-696	F								
WM-697	F								
WM-698	F								
WM-699	F								
WM-700	F								
WM-701	F								
WM-702	F								
WM-703	F								
WM-704	F								
WM-705	F								
WM-706	F								
WM-707	F								
WM-708	F								
WM-709	F								
WM-710	F								
WM-711	F								
WM-712	F								
WM-713	F								

Figure 30.



84125  
88712  
162354

161441

84472

84133



WM-714 F  
WM-715 F  
WM-716 F  
WM-717 F

Figure 4A

Rank	Normal vs Cancer	Adeno vs Normal	Squamous vs Normal	Small Cell vs Normal	Non-small Cell vs Normal	Large Cell vs Normal	Adeno vs Squamous	Adeno vs Small Cell	Squamous vs Small Cell
1	WM-61	WM-447	WM-447	WM-70	WM-341	WM-16	WM-62	WM-457	WM-276
2	WM-447	WM-692	WM-61	WM-706	WM-342	WM-26	WM-415	WM-72	WM-277
3	WM-446	WM-61	WM-277	WM-369	WM-343	WM-489	WM-152	WM-369	WM-362
4	WM-133	WM-446	WM-446	WM-447	WM-48	WM-134	WM-385	WM-78	WM-257
5	WM-119	WM-290	WM-133	WM-61	WM-340	WM-847	WM-347	WM-79	WM-363
6	WM-278	WM-363	WM-134	WM-652	WM-346	WM-277	WM-134	WM-73	WM-347
7	WM-134	WM-133	WM-363	WM-282	WM-47	WM-310	WM-38	WM-64	WM-53
8	WM-363	WM-341	WM-362	WM-446	WM-339	WM-363	WM-108	WM-320	WM-254
9	WM-282	WM-285	WM-276	WM-456	WM-389	WM-446	WM-99	WM-419	WM-17
10	WM-362	WM-366	WM-706	WM-134	WM-221	WM-151	WM-151	WM-85	WM-252
11	WM-120	WM-282	WM-203	WM-203	WM-447	WM-648	WM-289	WM-82	WM-431
12	WM-280	WM-382	WM-466	WM-646	WM-652	WM-657	WM-363	WM-53	WM-513
13	WM-65	WM-310	WM-386	WM-455	WM-154	WM-290	WM-61	WM-412	WM-446
14	WM-277	WM-282	WM-65	WM-65	WM-567	WM-328	WM-117	WM-440	WM-355
15	WM-70	WM-120	WM-70	WM-885	WM-456	WM-447	WM-211	WM-455	WM-447
16	WM-369	WM-134	WM-341	WM-473	WM-450	WM-684	WM-362	WM-458	WM-245
17	WM-17	WM-276	WM-429	WM-343	WM-283	WM-183	WM-133	WM-313	WM-133
18	WM-473	WM-428	WM-347	WM-466	WM-207	WM-190	WM-414	WM-56	WM-52
19	WM-47	WM-277	WM-17	WM-341	WM-436	WM-886	WM-277	WM-70	WM-86
20	WM-203	WM-20	WM-47	WM-340	WM-384	WM-397	WM-141	WM-246	WM-238
21	WM-276	WM-119	WM-431	WM-363	WM-61	WM-466	WM-64	WM-360	WM-243
22	WM-279	WM-340	WM-473	WM-339	WM-167	WM-20	WM-135	WM-180	WM-138
23	WM-62	WM-48	WM-62	WM-457	WM-382	WM-17	WM-447	WM-418	WM-62
24	WM-566	WM-389	WM-384	WM-86	WM-285	WM-545	WM-383	WM-83	WM-580
25	WM-456	WM-450	WM-438	WM-506	WM-650	WM-47	WM-338	WM-257	WM-134
26	WM-428	WM-47	WM-652	WM-72	WM-203	WM-191	WM-53	WM-138	WM-240
27	WM-384	WM-343	WM-282	WM-287	WM-119	WM-147	WM-142	WM-47	WM-268
28	WM-287	WM-17	WM-389	WM-82	WM-282	WM-480	WM-446	WM-252	WM-203
29	WM-420	WM-583	WM-290	WM-828	WM-686	WM-590	WM-186	WM-282	WM-111
30	WM-292	WM-70	WM-278	WM-85	WM-383	WM-218	WM-111	WM-50	WM-95
31	WM-431	WM-706	WM-456	WM-73	WM-429	WM-285	WM-445	WM-68	WM-247
32	WM-455	WM-346	WM-673	WM-138	WM-11	WM-652	WM-455	WM-325	WM-157
33	WM-20	WM-466	WM-340	WM-384	WM-208	WM-651	WM-276	WM-402	WM-242
34	WM-340	WM-646	WM-55	WM-83	WM-451	WM-368	WM-444	WM-411	WM-556
35	WM-19	WM-384	WM-455	WM-450	WM-473	WM-403	WM-181	WM-405	WM-53
36	WM-389	WM-336	WM-645	WM-310	WM-220	WM-410	WM-35	WM-75	WM-239
37	WM-53	WM-294	WM-138	WM-277	WM-685	WM-430	WM-285	WM-417	WM-234
38	WM-438	WM-339	WM-420	WM-79	WM-338	WM-456	WM-456	WM-387	WM-274

Figure 4B


39	WM-450	WM-473	WM-450	WM-71	WM-714	WM-39	WM-28	WM-370
40	WM-468	WM-369	WM-369	WM-268	WM-546	WM-82	WM-410	WM-301
41	WM-296	WM-38	WM-279	WM-70	WM-109	WM-17	WM-420	WM-449
42	WM-343	WM-283	WM-342	WM-545	WM-302	WM-203	WM-164	WM-74
43	WM-341	WM-685	WM-471	WM-575	WM-597	WM-83	WM-87	WM-261
44	WM-339	WM-66	WM-674	WM-448	WM-375	WM-412	WM-66	WM-467
45	WM-55	WM-55	WM-120	WM-120	WM-131	WM-96	WM-391	WM-237
46	WM-55	WM-650	WM-20	WM-147	WM-706	WM-74	WM-340	WM-262
47	WM-307	WM-287	WM-466	WM-257	WM-398	WM-457	WM-428	WM-295
48	WM-48	WM-278	WM-347	WM-456	WM-309	WM-431	WM-198	WM-288
49	WM-38	WM-83	WM-154	WM-153	WM-55	WM-340	WM-312	WM-384
50	WM-138	WM-342	WM-128	WM-38	WM-488	WM-49	WM-152	WM-37
	WM-310	WM-429						




[illegible]

## Figure 5


Figure 6

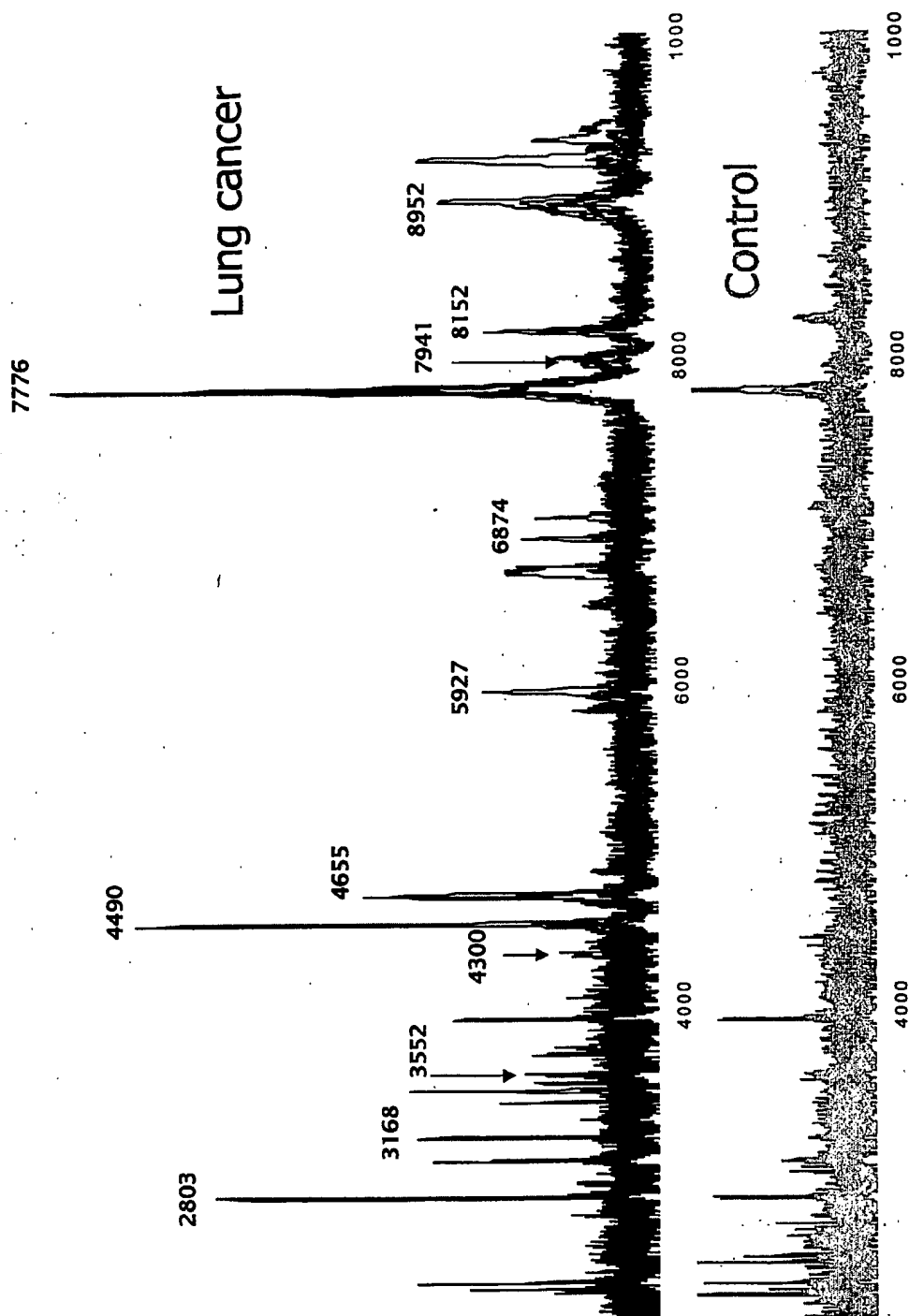


**FIGURE 7**  
**Protein Profile of Selected Samples Q Fraction 1 WCX2**

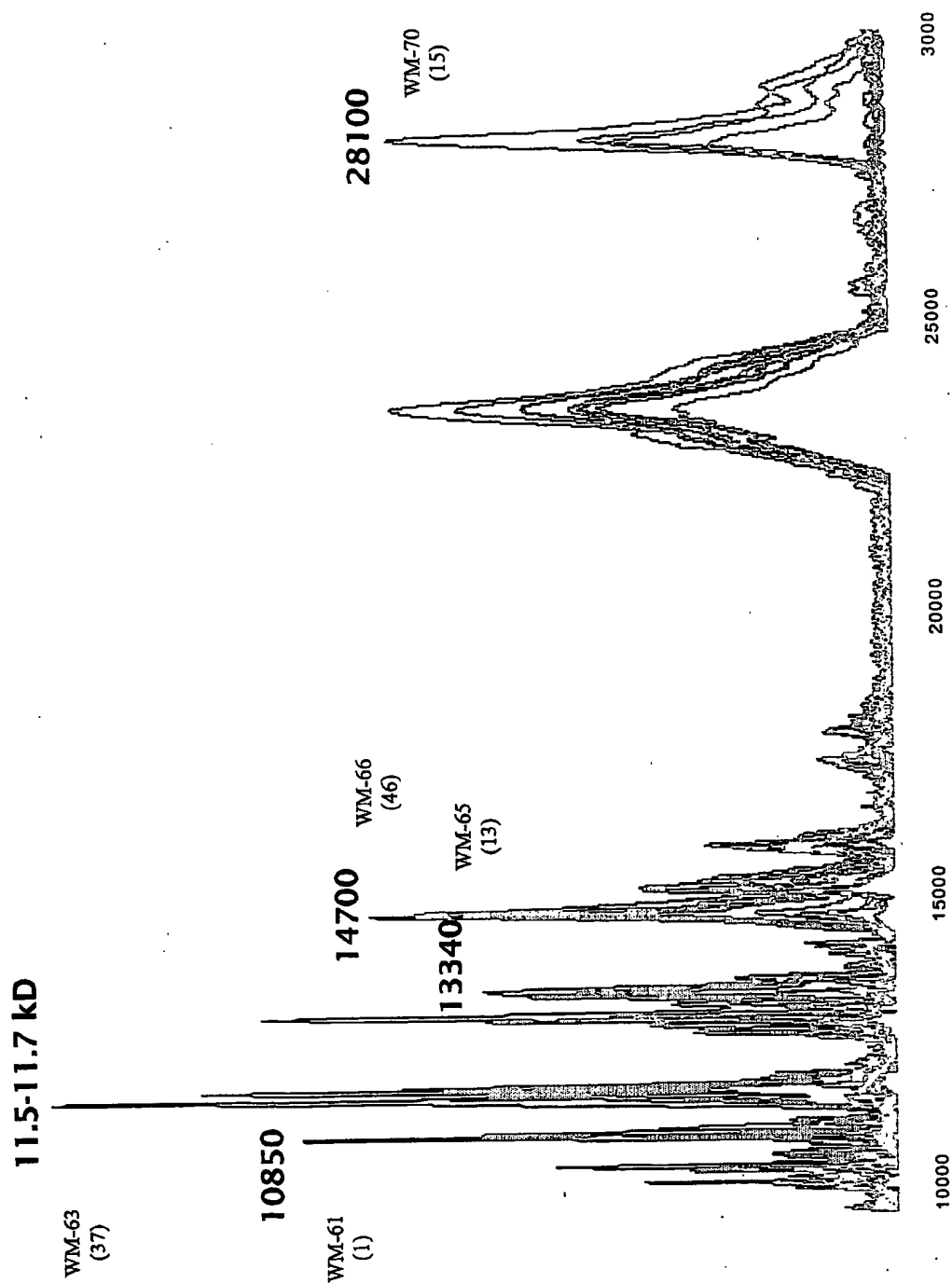


**Figure 8**  
**Protein Profile of Selected Samples**  
**Q Fraction 1 WCX2**

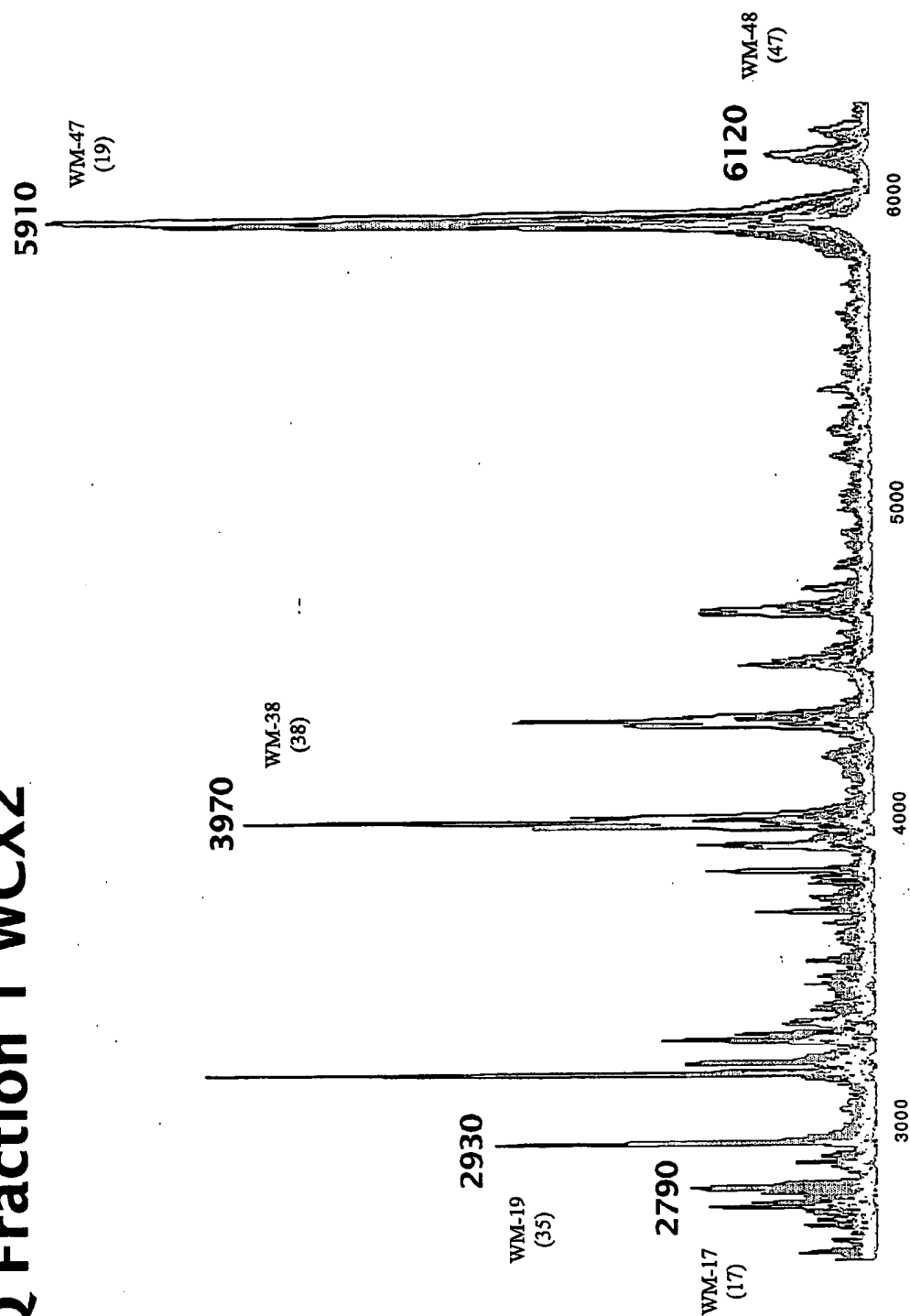


**Figure 9**  
**Protein Profile of Selected Samples**  
**Q Fraction 2 WCX2**

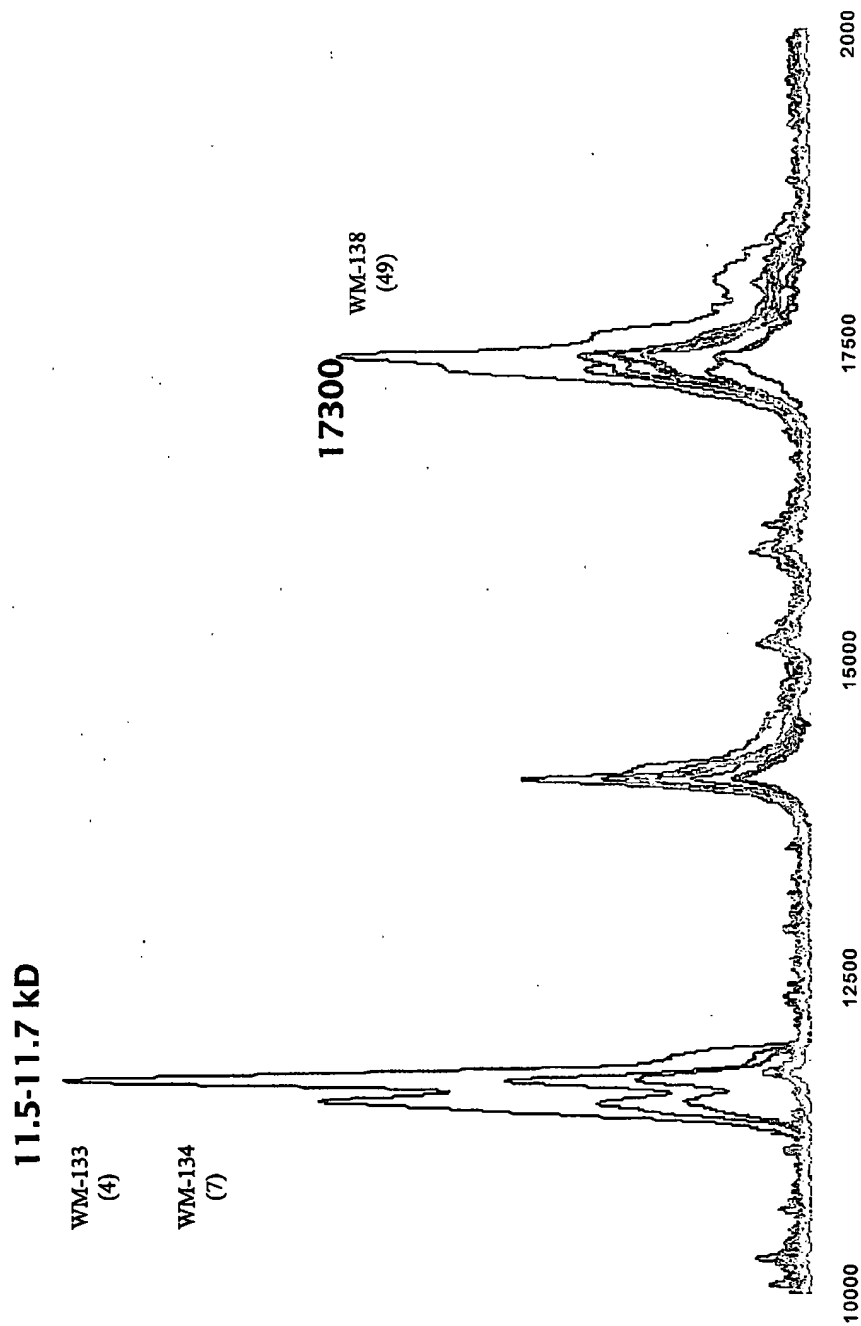


**Figure 10**  
**Protein Profile of Selected Samples**  
**Q Fraction 2 WCX2**

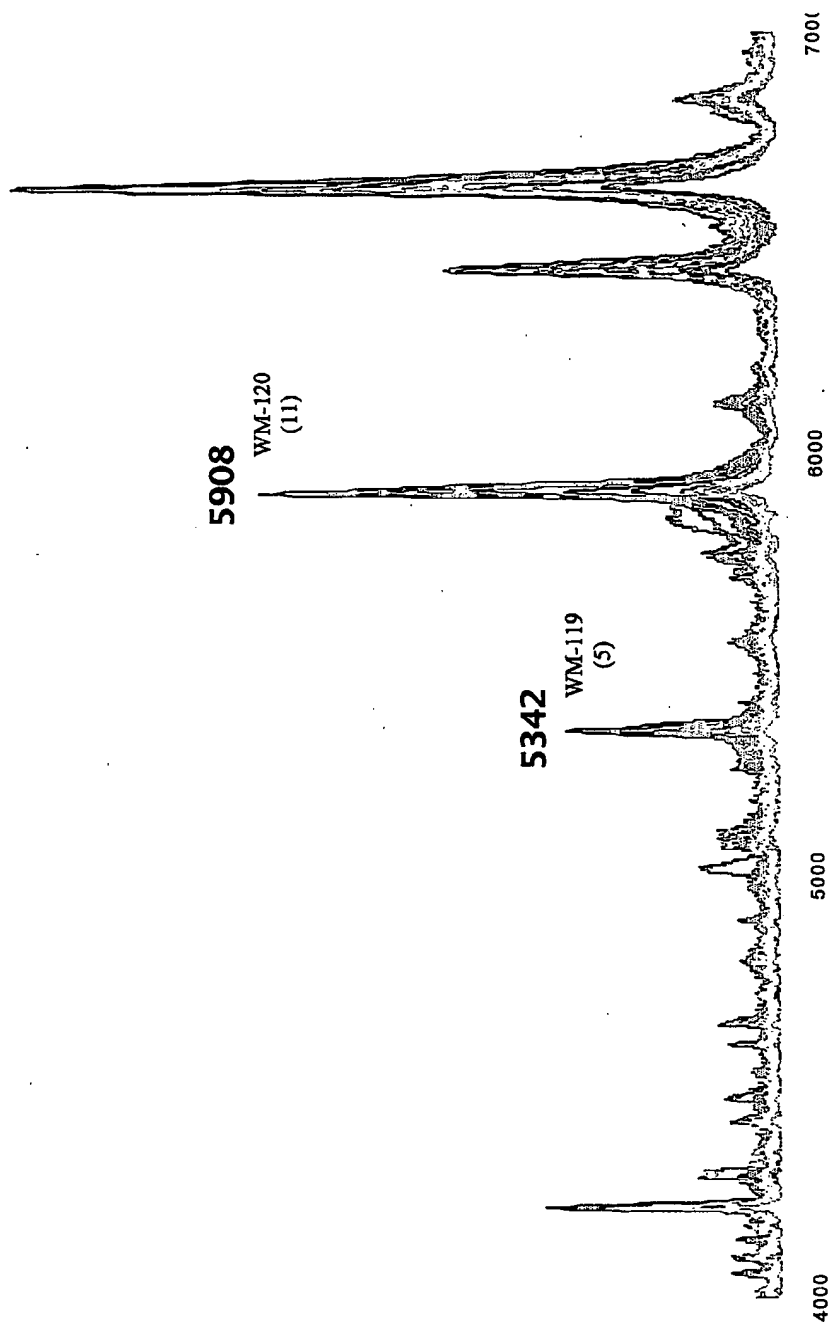

**Figure 11**  
**Protein Profile of Selected Samples**  
**Q Fraction 4 WCX2**

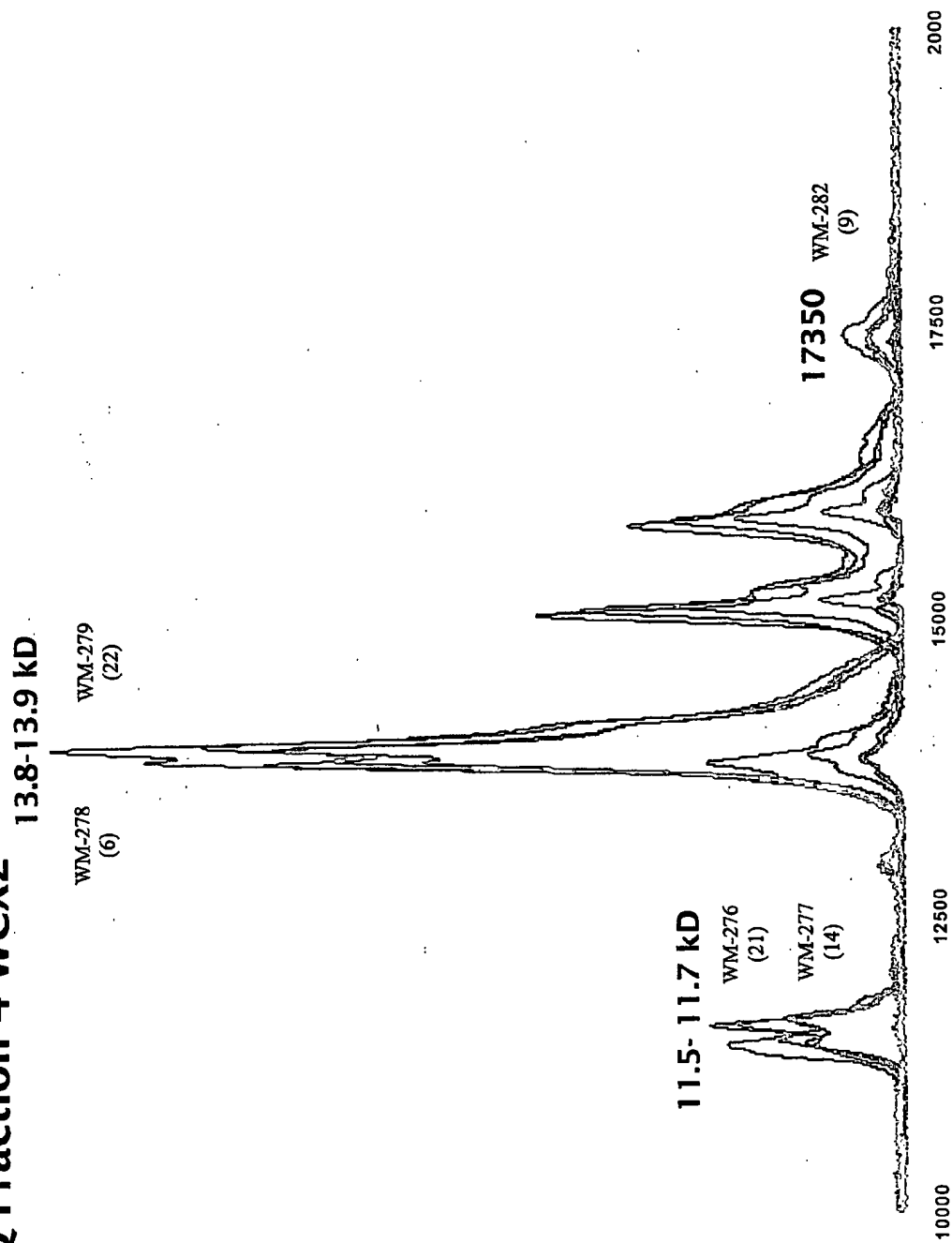


**Figure 12**  
**Protein Profile of Selected Samples**  
**Q Fraction 4 WCX2**

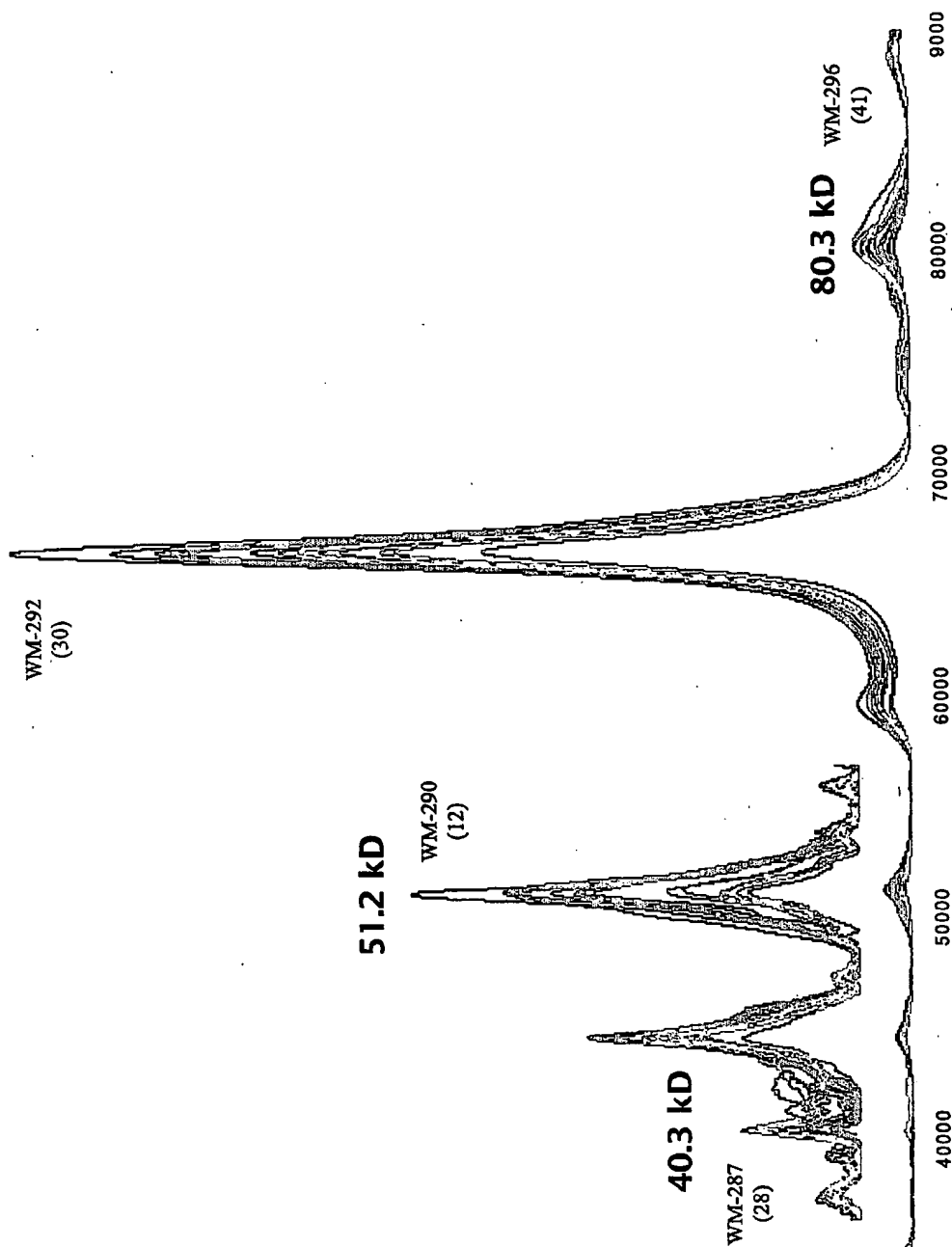





**Figure 13**  
**Protein Profile of Selected Samples**  
**Q Fraction 5 WCX2**

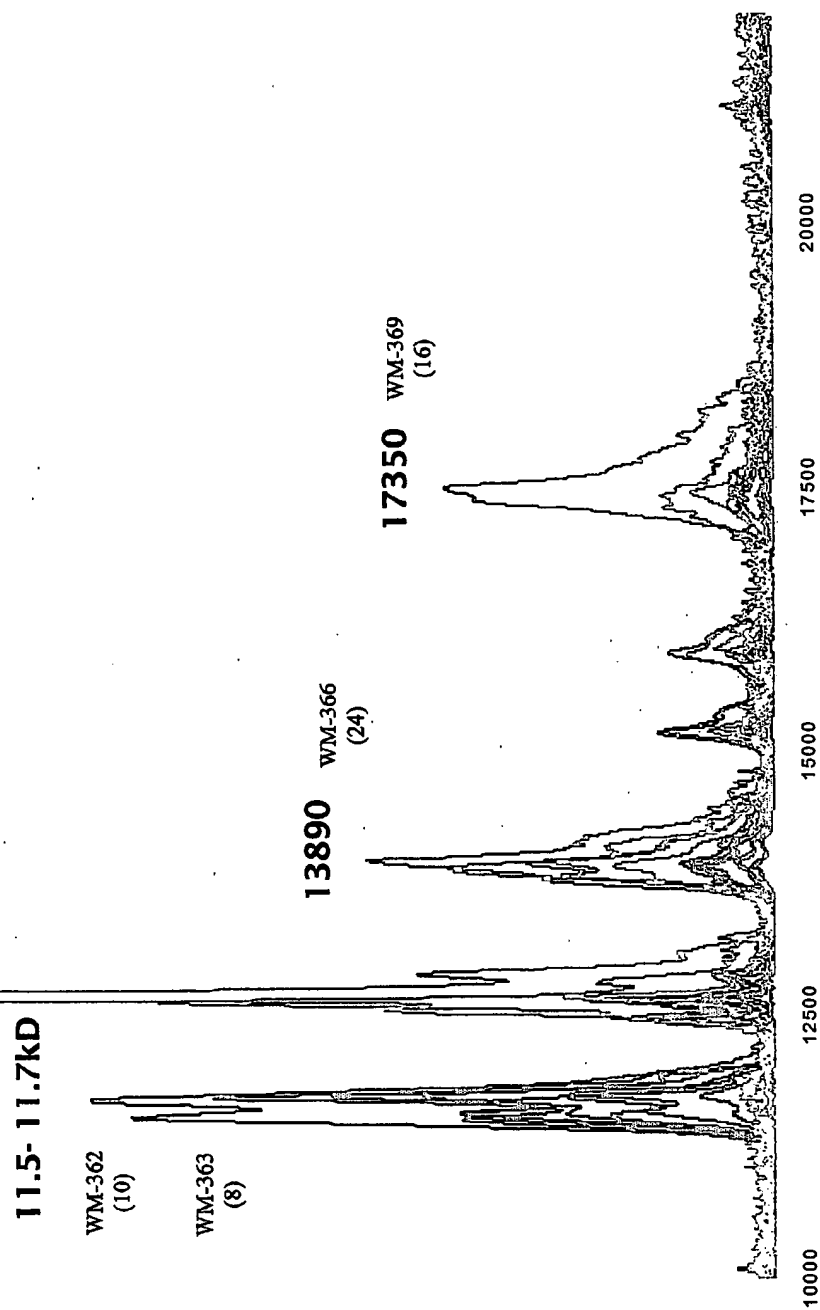


**Figure 14**  
**Protein Profile of Selected Samples**  
**Q Fraction 5 WCX2**

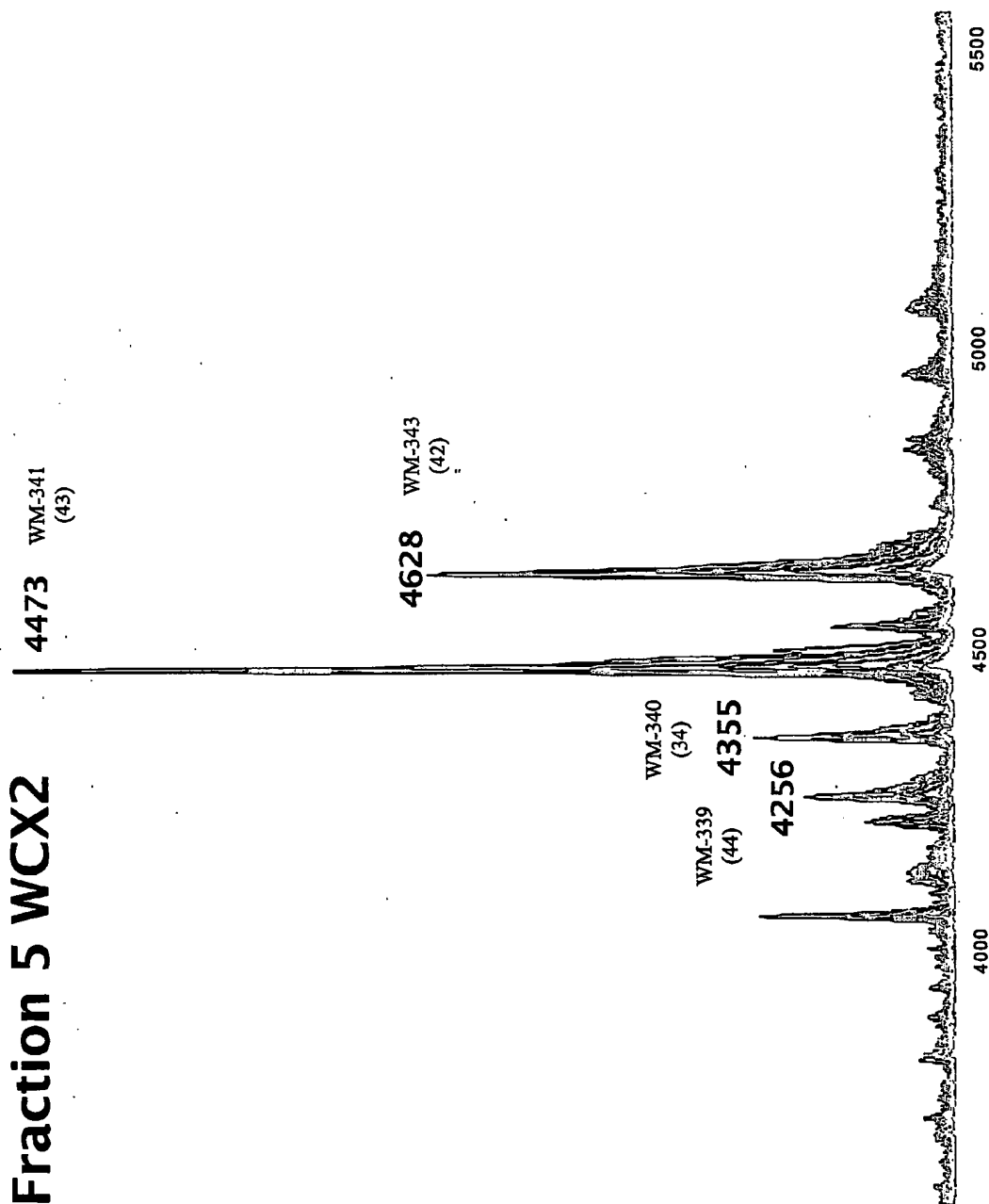


**Figure 15**  
**Protein Profile of Selected Samples**  
**Q Fraction 6 WCX2**

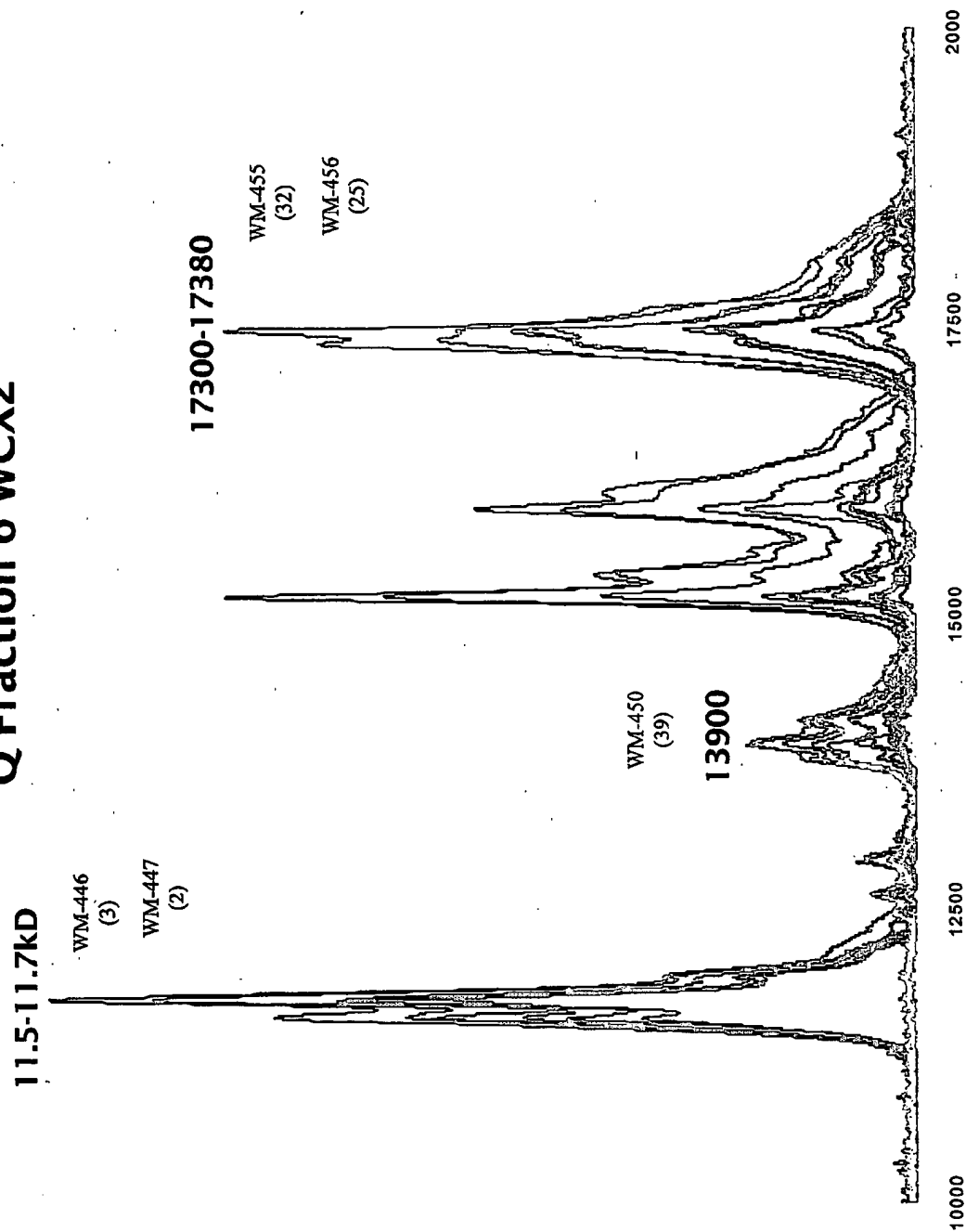


**Figure 16**  
**Protein Profile of Selected Samples**  
**Q Fraction 6 WCX2**

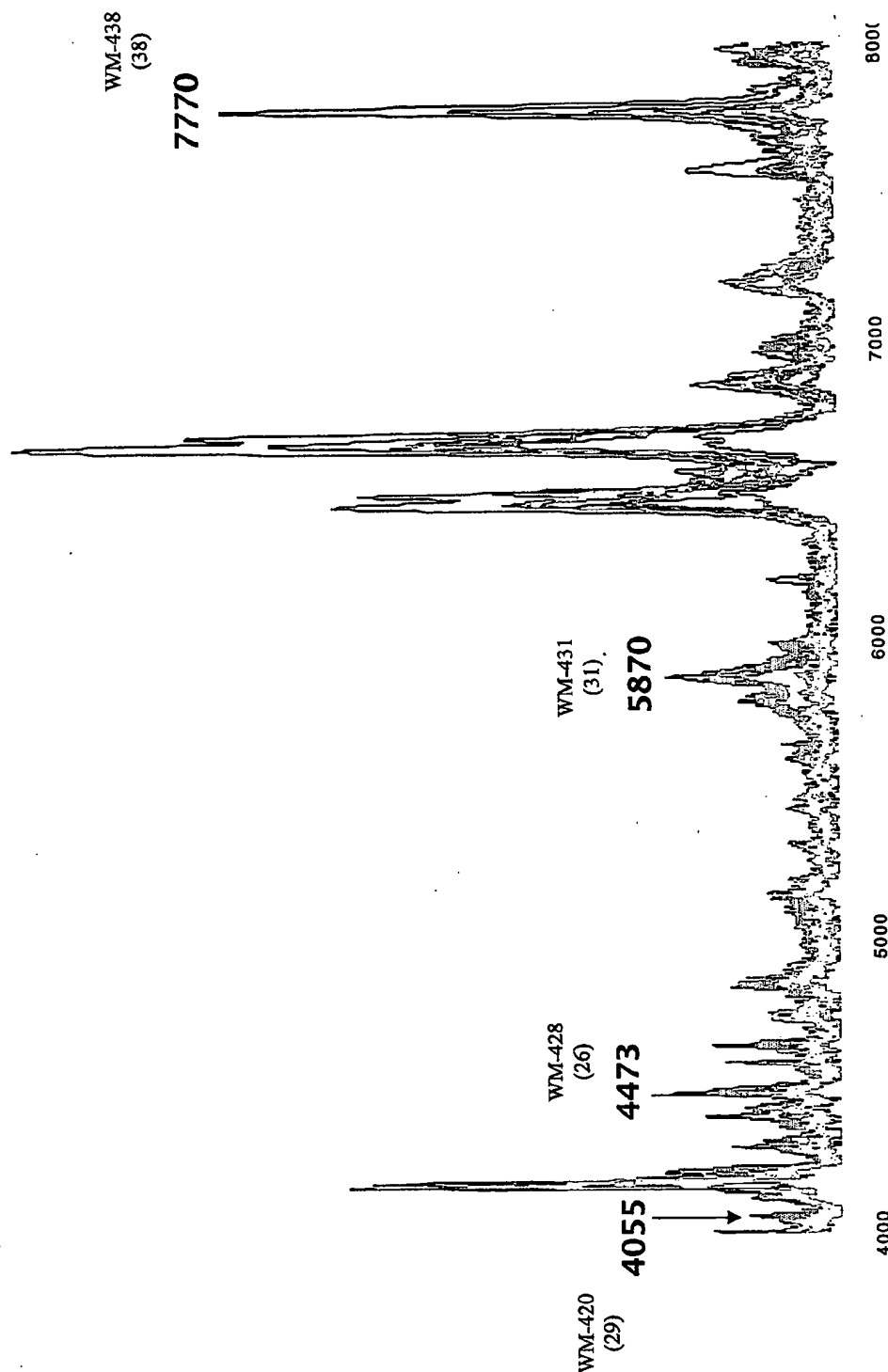


**Figure 17**  
**Protein Profile of Selected Samples**  
**Q Fraction 2 IMAC-Cu(II)**

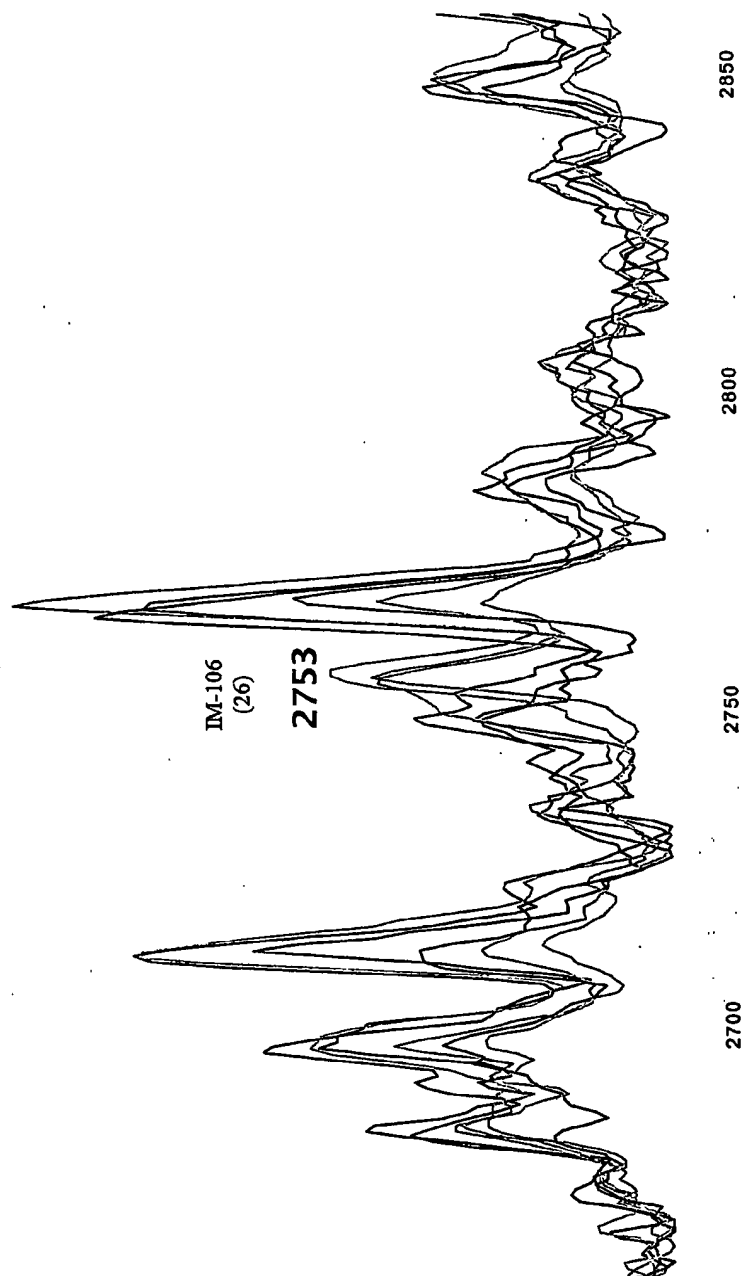


**Figure 18**  
**Protein Profile of Selected Samples**  
**Q Fraction 2 IMAC-Cu(II)**

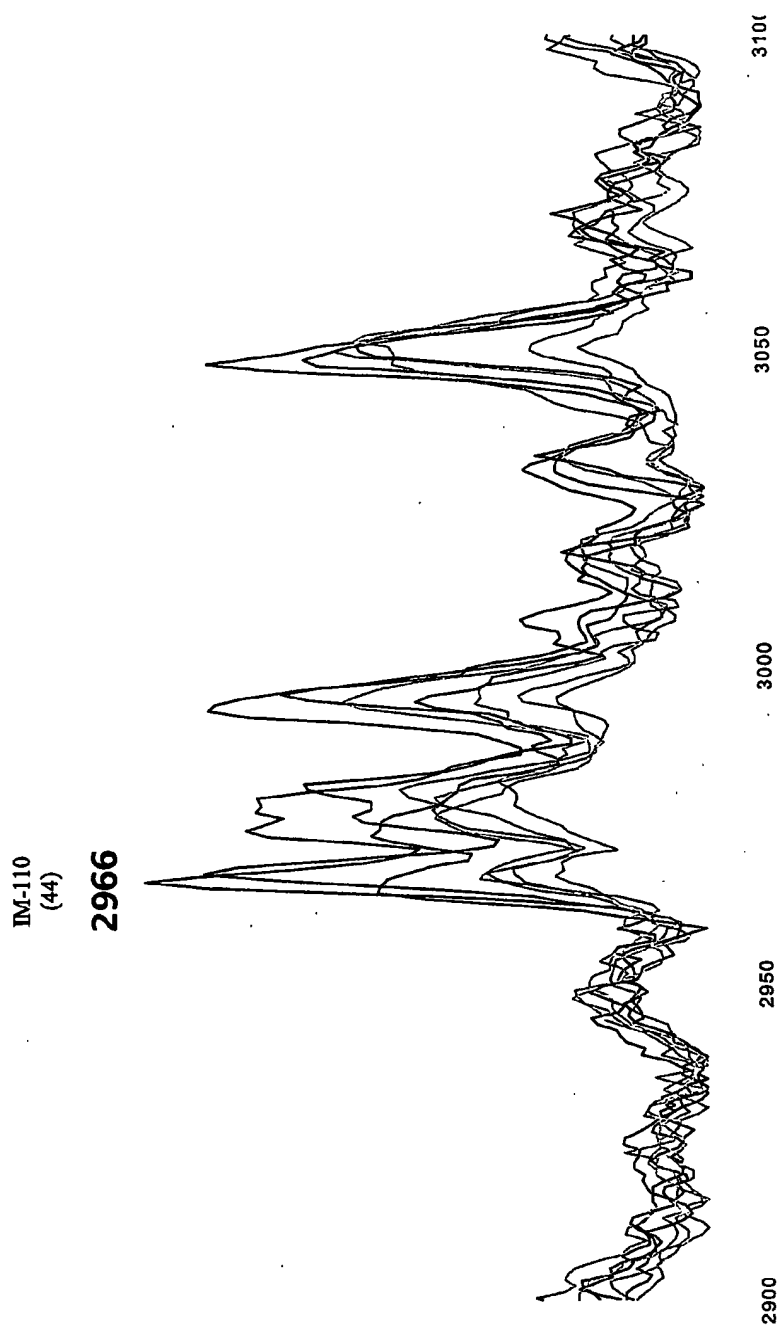

**Figure 19**  
**Protein Profile of Selected Samples**  
**Q Fraction 2 IMAC-Cu(II)**

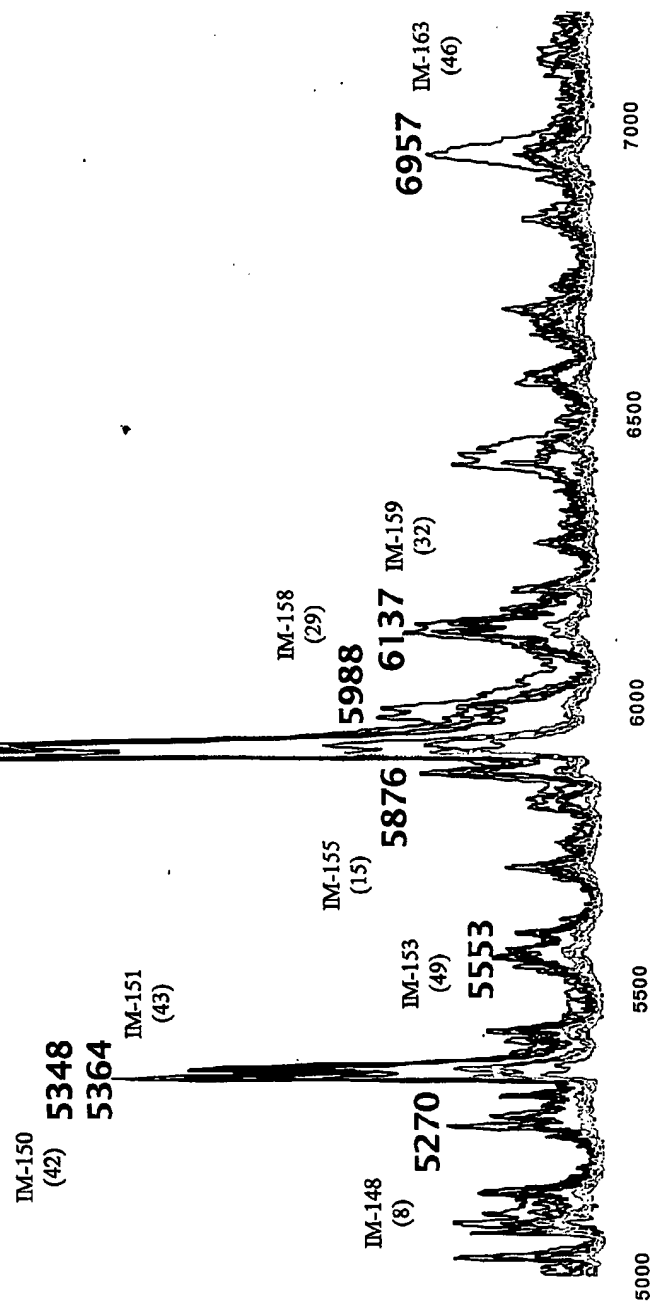


**Figure 20**  
**Protein Profile of Selected Samples**  
**Q Fraction 2 IMAC-Cu(II)**

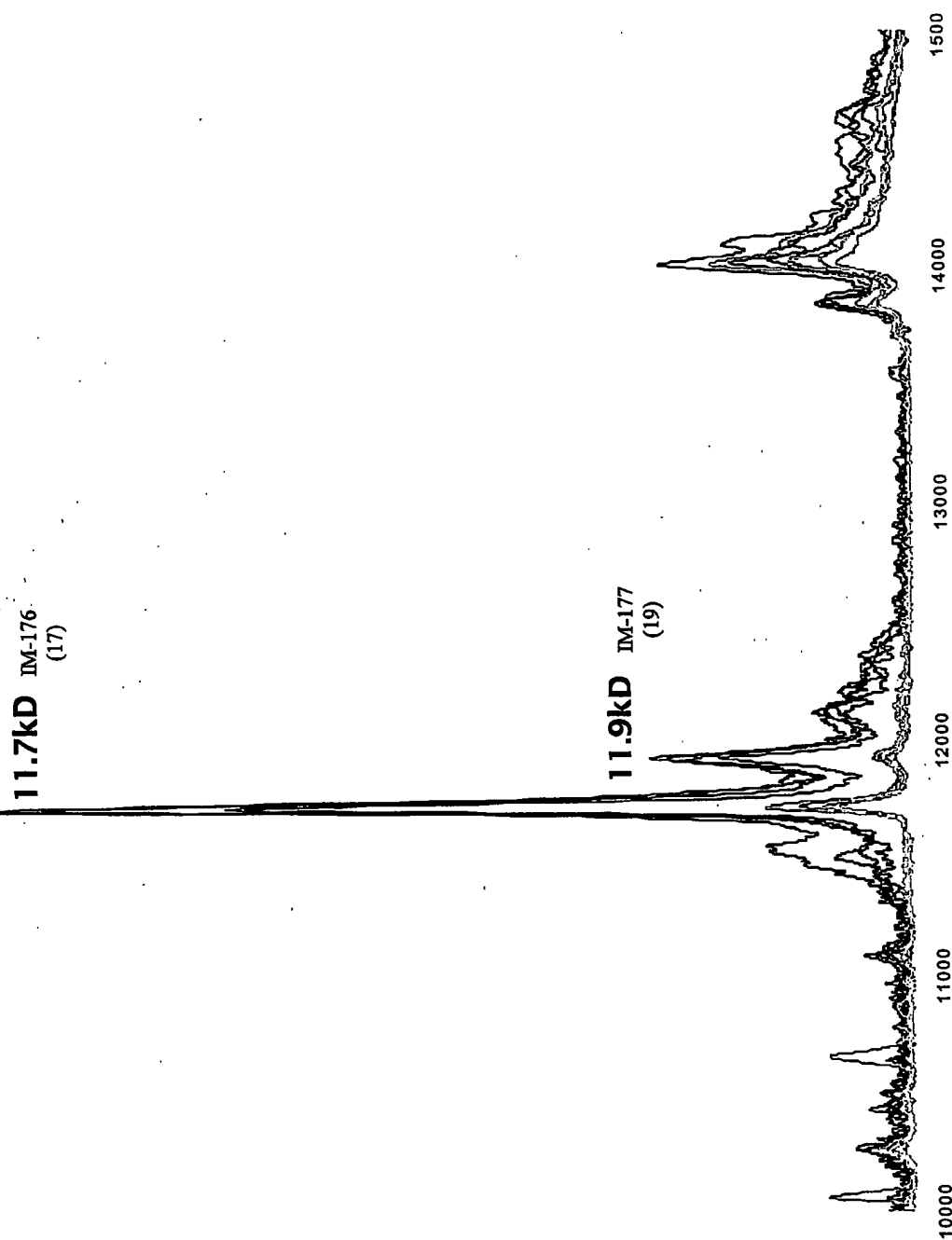




**Figure 21**  
**Protein Profile of Selected Samples**  
**Q Fraction 3 IMAC-Cu(II)**

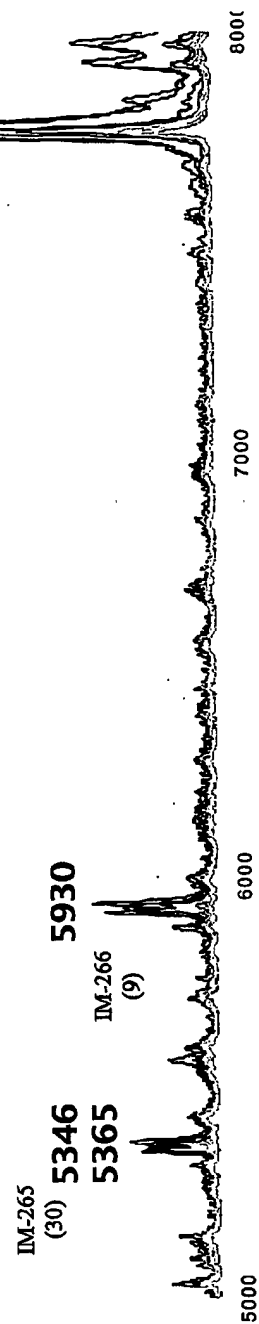


**Figure 22**  
**Protein Profile of Selected Samples**  
**Q Fraction 3 IMAC-Cu(II)**

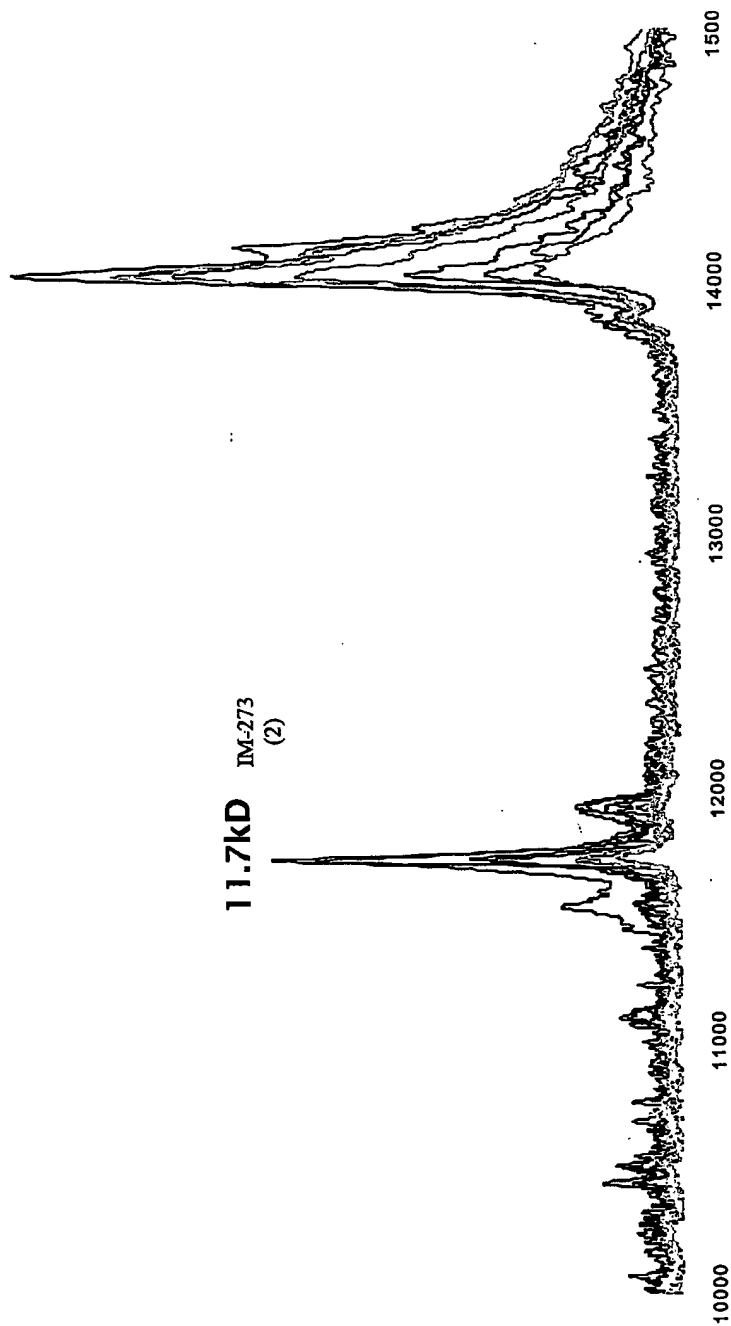


**Figure 23**  
**Protein Profile of Selected Samples**


**Q Fraction 5 IMAC-Cu(II)**

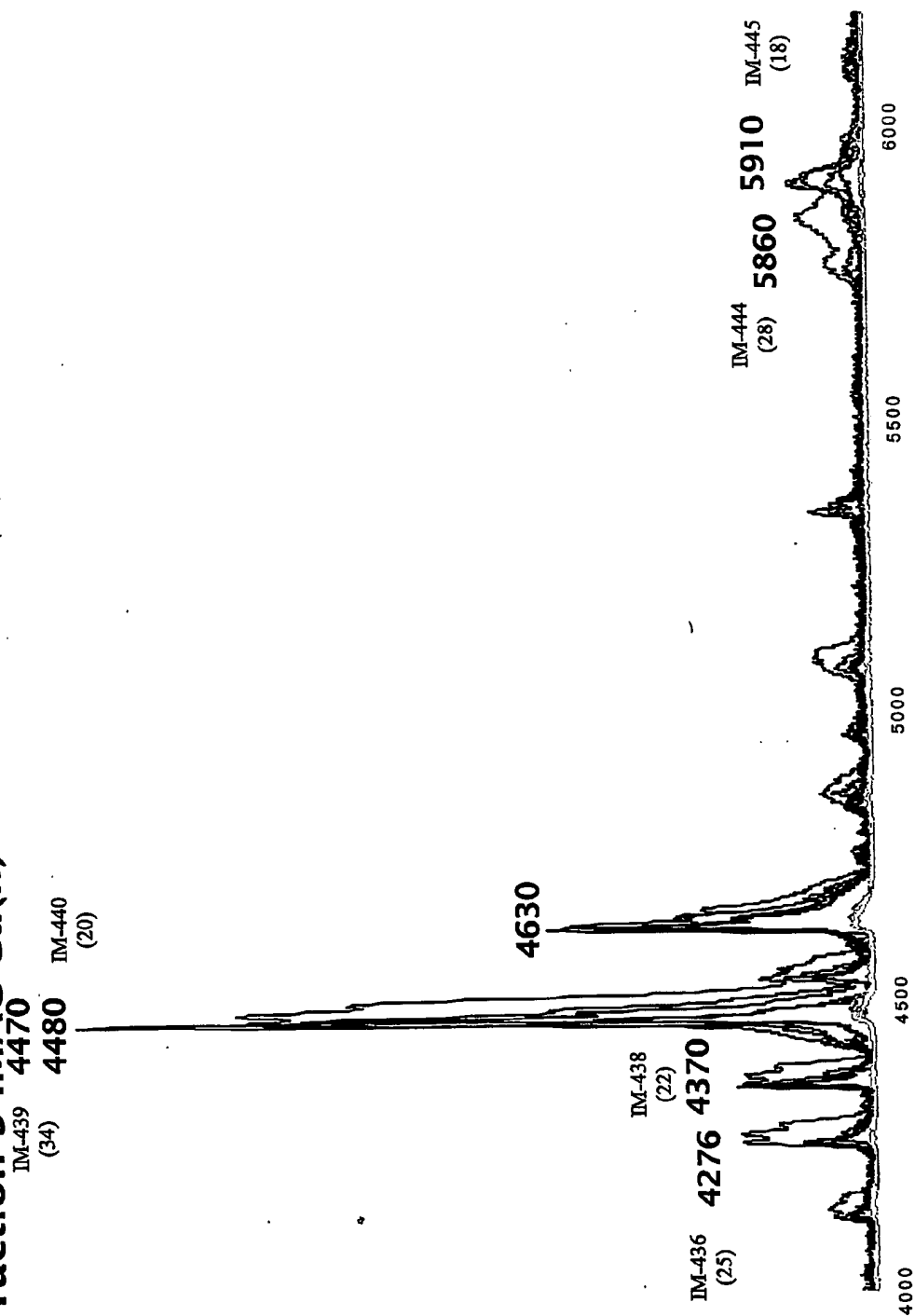


**Figure 24**  
**Protein Profile of Selected Samples**  
**Q Fraction 5 IMAC-Cu(II)**

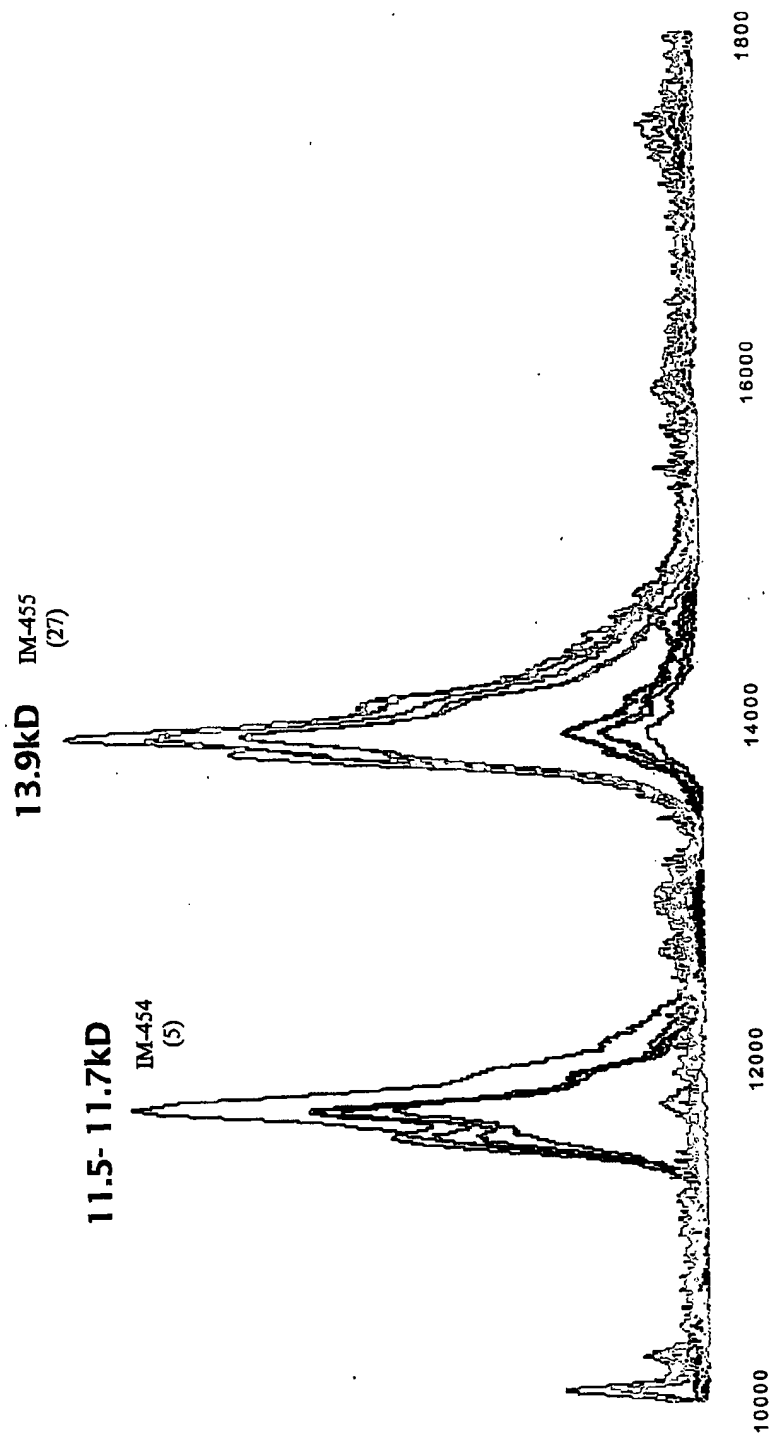


**Figure 25**  
**Protein Profile of Selected Samples**  
**Q Fraction 5 IMAC-Cu(II)**

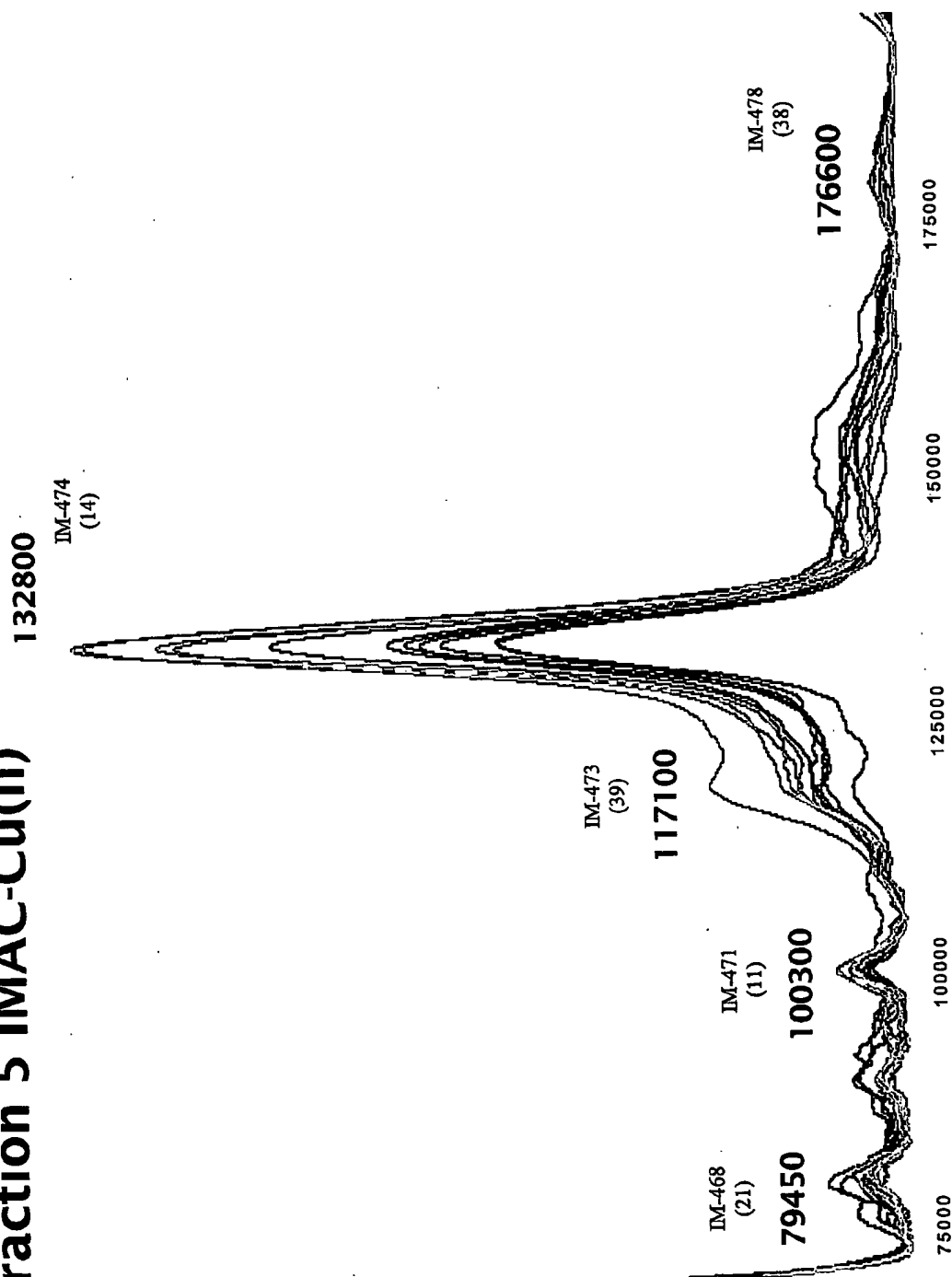


**Figure 26**  
**Protein Profile of Selected Samples**  
**Q Fraction 6 IMAC-Cu(II)**

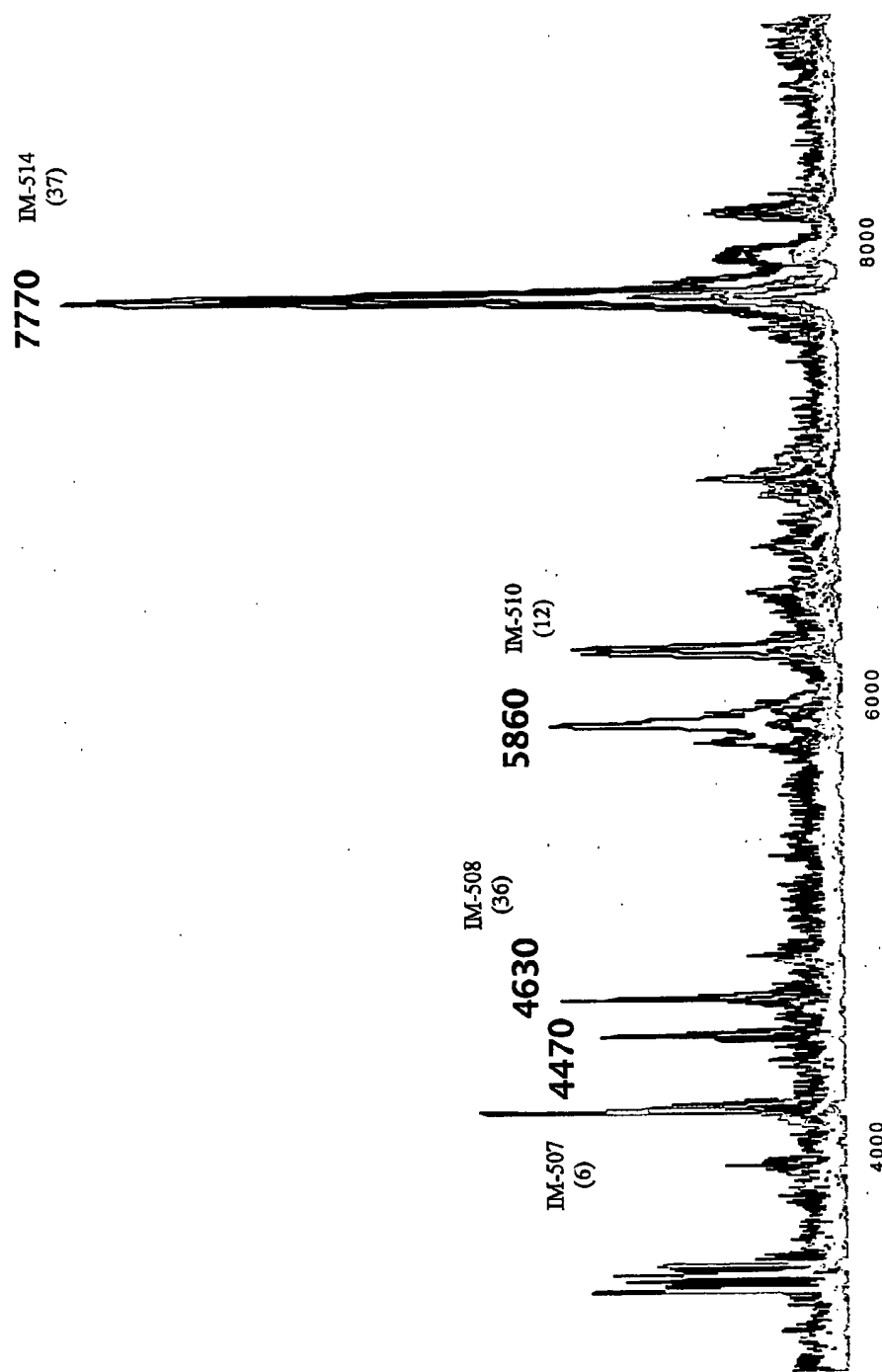

**Figure 27**  
**Protein Profile of Selected Samples**  
**Q Fraction 6 IMAC-Cu(II)**

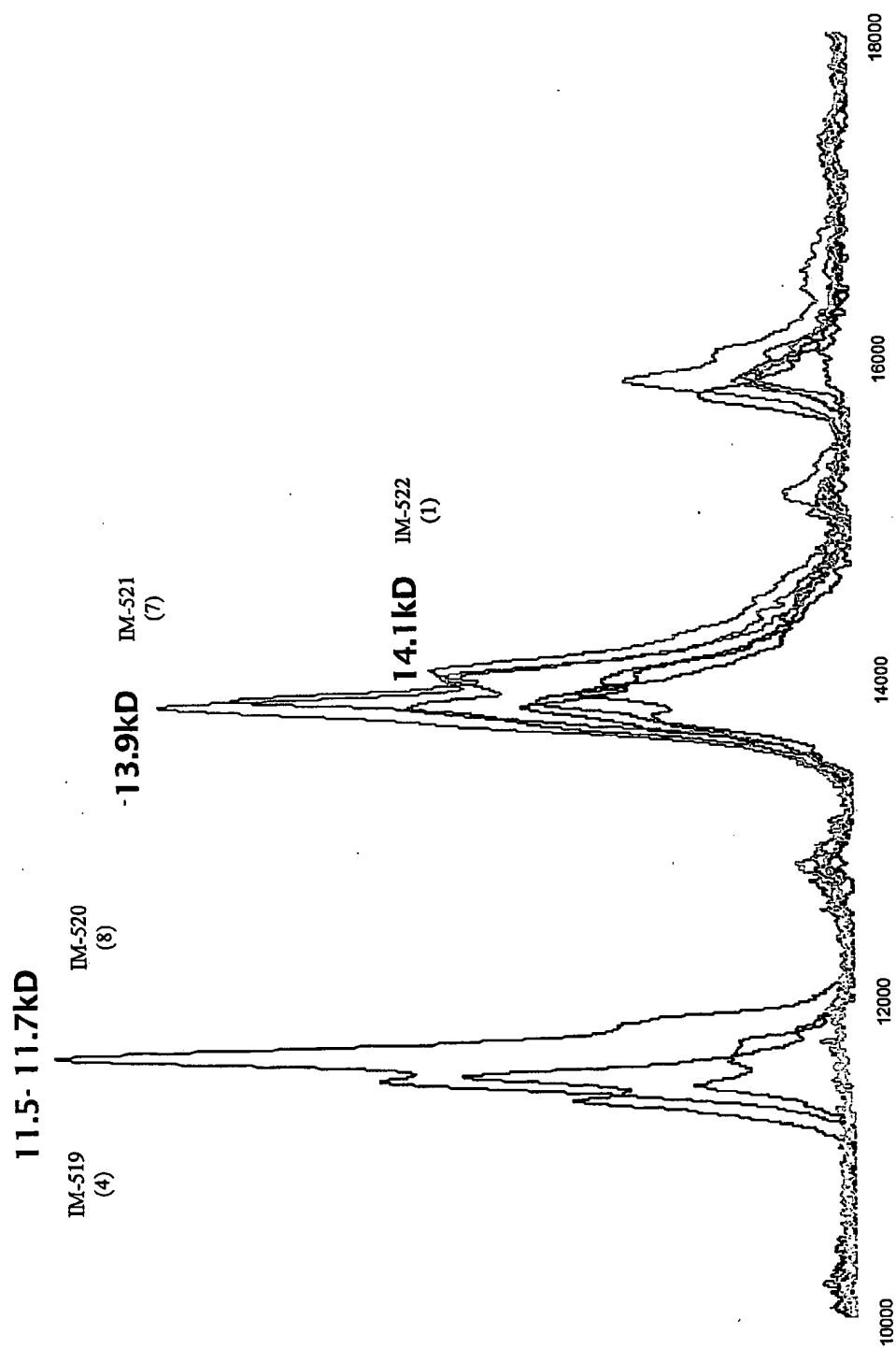

**Figure 28**  
**Protein Profile of Selected Samples**  
**Q Fraction 6 IMAC-Cu(II)**

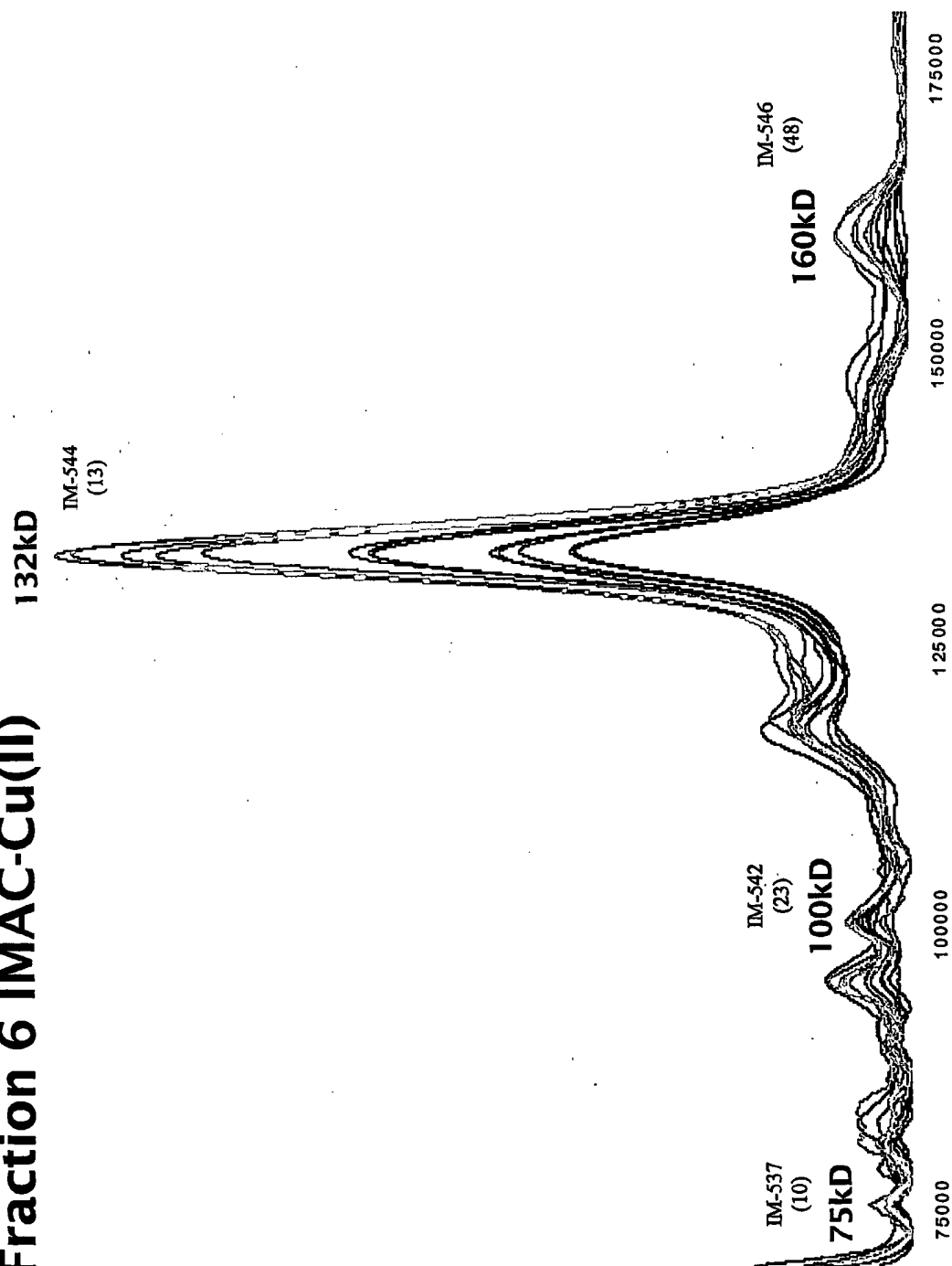



(19) World Intellectual Property  
Organization  
International Bureau



(43) International Publication Date  
22 July 2004 (22.07.2004)

PCT

(10) International Publication Number  
**WO 2004/061410 A3**

(51) International Patent Classification<sup>7</sup>: **C12Q 1/00**,  
G01N 33/53, A61K 49/00

(74) Agents: **BENT, Stephen, A.** et al.; Foley & Lardner,  
Washington Harbour, 3000 K Street, N.W. Suite 500,  
Washington, DC 20007-5101 (US).

(21) International Application Number:  
**PCT/US2003/037090**

(22) International Filing Date:  
16 December 2003 (16.12.2003)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:  
60/434,075 18 December 2002 (18.12.2002) US

(81) Designated States (*national*): AE, AG, AL, AM, AT, AU,  
AZ, BA, BB, BG, BR, BW, BY, BZ, CA, CH, CN, CO, CR,  
CU, CZ, DE, DK, DM, DZ, EC, EE, EG, ES, FI, GB, GD,  
GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR,  
KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN,  
MW, MX, MZ, NI, NO, NZ, OM, PG, PH, PL, PT, RO, RU,  
SC, SD, SE, SG, SK, SL, SY, TJ, TM, TN, TR, TT, TZ, UA,  
UG, US, UZ, VC, VN, YU, ZA, ZM, ZW.

(84) Designated States (*regional*): ARIPO patent (BW, GH,  
GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW),  
Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM),  
European patent (AT, BE, BG, CH, CY, CZ, DE, DK, EE,  
ES, FI, FR, GB, GR, HU, IE, IT, LU, MC, NL, PT, RO, SE,  
SI, SK, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA,  
GN, GQ, GW, ML, MR, NE, SN, TD, TG).

(71) Applicants (*for all designated States except US*): **CI-  
PHERGEN BIOSYSTEMS, INC.** [US/US]; 6611  
Dumbarton Circle, Freemont, CA 94555 (US). **QUEEN  
ELIZABETH HOSPITAL** [CN/CN]; 30 Gascoigne  
Road, Kowloon, Hong Kong SAR (CN).

(72) Inventors; and

(75) Inventors/Applicants (*for US only*): **YIP, Timothy, Tak,  
Chun** [CN/CN]; 1, Kapok Path, Westwood, Palm Springs,  
Yuen Long N.T., Hong Kong SAR (CN). **CHO, Chi, Shing**  
[—/CN]; Flat 9E, Block 4, The Tolo Place, Sunshine City,  
Ma On Shan, N.T., Hong Kong SAR (CN). **AU, Siu, Kie**  
[—/CN]; c/o Department of Clinical Oncology, Queen  
Elizabeth Hospital, 30 Gascoigne Road, Kowloon, Hong  
Kong SAR (CN). **YIP, Tai-Tung** [US/US]; 1532 Aster  
Court, Cupertino, CA 95014 (US). **YIP, Christine, L.**  
[US/US]; 1532 Aster Court, Cupertino, CA 95014 (US).  
**YIP, Victor, F.** [US/US]; 33008 Compton CT, Union City,  
CA 94587 (US).

Published:

— *with international search report*  
— *before the expiration of the time limit for amending the  
claims and to be republished in the event of receipt of  
amendments*

(88) Date of publication of the international search report:  
29 December 2004

*For two-letter codes and other abbreviations, refer to the "Guid-  
ance Notes on Codes and Abbreviations" appearing at the begin-  
ning of each regular issue of the PCT Gazette.*

(54) Title: SERUM BIOMARKERS IN LUNG CANCER

(57) Abstract: Certain biomarkers and biomarker combinations are useful in a qualifying lung cancer status in a subject. A diag-  
nostic methodology employing these biomarkers and combinations can detect whether a subject has lung cancer.



WO 2004/061410 A3

# INTERNATIONAL SEARCH REPORT

International application No.

PCT/US03/37090

## A. CLASSIFICATION OF SUBJECT MATTER

IPC(7) : C12Q 1/00; G01N 33/53; A61K 49/00

US CL : 435/4; 435/7.1; 424/9.1

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 435/4; 435/7.1; 424/9.1

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)  
MEDLINE, STN, WEST

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	US 2003/0091976 (BOSCHETTI et al.) 15 May 2003 (15.05.2003), entire article.	1-102
A	CHAPMAN, K. The ProteinChip Biomarker System from Ciphergen Biosystems: a novel proteomics platform for rapid biomarker discovery and validation. Biochemical Society Transactions. April 2002, Vol. 30, part 2, pages 82-87.	1-102
A	POON et al. Comprehensive Proteomic Profiling Identifies Serum Proteomic Signatures for Detection of Hepatocellular Carcinoma and Its Subtypes. Clinical Chemistry. May 2003. Vol. 49, No.5, pages 752-760.	1-102

☐ Further documents are listed in the continuation of Box C.

☐ See patent family annex.

\* Special categories of cited documents:

- "A" document defining the general state of the art which is not considered to be of particular relevance
- "E" earlier application or patent published on or after the international filing date
- "L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- "O" document referring to an oral disclosure, use, exhibition or other means
- "P" document published prior to the international filing date but later than the priority date claimed

- "T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- "X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- "Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art
- "&" document member of the same patent family

Date of the actual completion of the international search

01 November 2004 (01.11.2004)

Date of mailing of the international search report

18 NOV 2004

Name and mailing address of the ISA/US

Mail Stop PCT, Attn: ISA/US  
Commissioner for Patents  
P.O. Box 1450  
Alexandria, Virginia 22313-1450

Facsimile No. (703) 305-3230

Authorized officer

Gary B. Nickol Ph.D.

Telephone No. 571-272-1600

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☒ FADED TEXT OR DRAWING
- ☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☒ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**